
BESTRON®

DSA3500C

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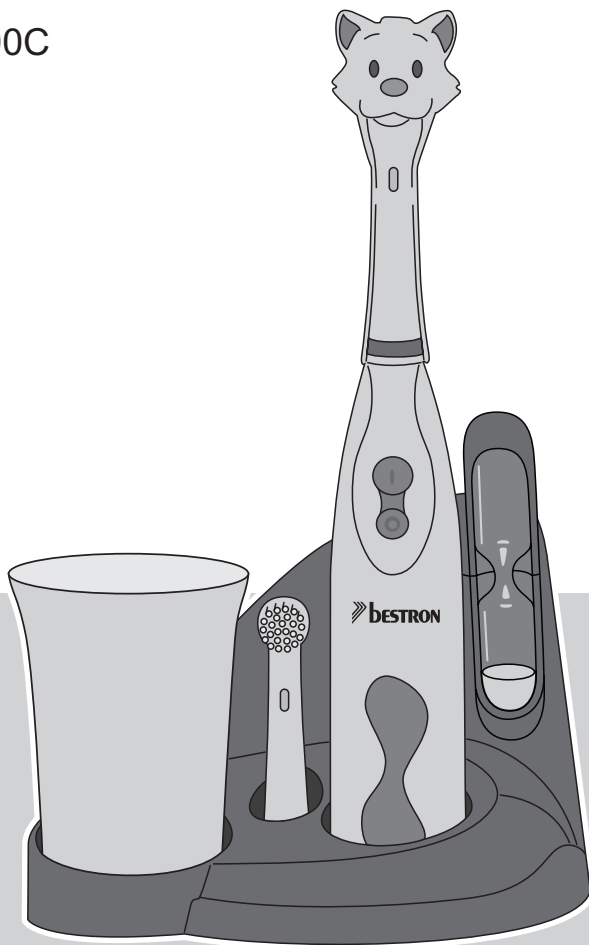Congratulations with the purchase of your toothbrush set.

SAFETY INSTRUCTIONS - General

- Please read these instructions carefully and retain them for reference.
- Use this appliance solely in accordance with these instructions.
- Without supervision this appliance may not be used by children or persons who have a physical or mental restriction and therefore are not able to use the appliance safely.
- Have any repairs carried out solely by a qualified electrician. Never try to repair the appliance yourself.
- Children younger than 3 years should not use this toothbrush without adult supervision.
- Always switch the appliance off if:
 - It malfunctions while you are using it.
 - Before you start maintenance.
 - You fit or remove accessories.
 - You have finished using the appliance.

TECHNICAL DETAILS

Type: DSA3500C
 Battery: 2x AA 1.5 V


THE ENVIRONMENT

-  • Dispose of packaging material, such as plastic and boxes, in the appropriate waste containers.
- When this product reaches the end of its useful life do not dispose of it by putting it in a dustbin; hand it in at a collection point for the reuse of electrical and electronic equipment. Please refer to the symbols on the product, the user's instructions or the packaging.
- The materials can be re-used as indicated. Your help in the re-use, recycling or other means of making use of old electrical equipment will make an important contribution to the protection of the environment.
- Contact your municipality for the address of the appropriate collection point in your neighbourhood.

OPERATION - General

Your toothbrush set is intended only for domestic use.

English

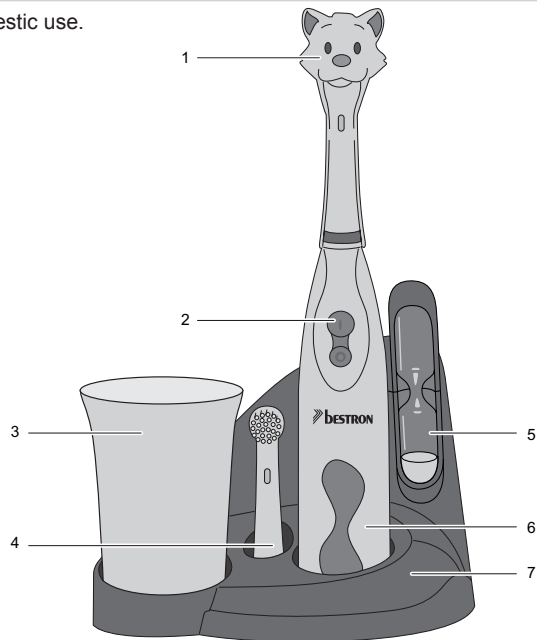1. Protective cover
2. On/Off switch
3. Cup
4. Brush
5. Sandglass
6. Housing
7. Holder

Figure 1

OPERATION - Use

- Fill the beaker (3) with water and wash your mouth with it.
- Place a brush (4) on the appliance till it fastens.
- Put toothpaste on the brush.
- Put the brush in your mouth before you switch on the appliance to prevent the toothpaste splatters around.
- Turn the sand timer (5) on switch on the appliance (2).
- Move the brush along your teeth and gum. There is no need for too much pressure.
- Also brush the other side of your teeth.
- The sand timer is at the end after two minutes. This is the minimum time to brush your teeth right.
- After brushing switch off the appliance and rinse it under the tap. After that put the cover (1) on the brush and place the appliance in the base (7).

For optimal brushing, we advise you to change the brush head every three months or sooner if the brush hairs begins to widen. The brush head is easy to change by pulling it off the toothbrush. The new brush head can be attached very easily by pushing it on the toothbrush. New brush heads can be purchased at the service center.

 The first few times you brush with an electrical toothbrush, you can suffer some bleeding gums. After your gun is used to the toothbrush, this will stop. When the bleeding keeps coming back, please contact your dentist.

OPERATION - Changing batteries

For optimal brushing, we advise you to change the batteries when the motor rotates slower than it normally does. Remove the empty batteries by open the battery compartment at the bottom of the brush. Place the full batteries in the battery compartment and close the lid.

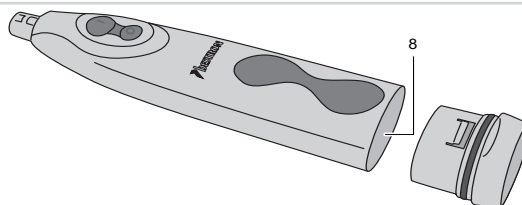
Figure 1

8. Battery compartment

CLEANING AND MAINTENANCE

- For optimal brushing, we advise you to change the batteries when the motor rotates slower than it normally does. Remove the empty batteries by open the battery compartment at the bottom of the brush. Place the full batteries in the battery compartment and close the lid.
- When not using the tooth brush for a longer period, you should remove the batteries.
- From a hygienic point of view from time to time you can clean the brush with tooth cleaner fluid.
- Don't use any aggressive cleansers.

GUARANTEE TERMS

The importer guarantees the appliance against defects caused by manufacturing and/or material faults for 60 months from the date of purchase, subject to the following conditions.

1. No charges will be made for labour costs or materials during the warranty period. Any repairs carried out during the warranty period do not extend the warranty period.
2. The importer automatically becomes the owner of any faulty parts that he replaces or a faulty appliance (in the event of replacement of the appliance).
3. All claims submitted under warranty must be accompanied by the original bill of sale.
4. Claims under warranty must be submitted to the dealer where the appliance was purchased or to the importer.
5. The warranty is issued solely to the first purchaser and cannot be transferred.
6. The warranty does not cover damage caused by:
 - a. Accidents, incorrect use, wear and/or neglect.
 - b. Faulty installation and/or use in a manner contravening the prevailing legal, technical or safety regulations.
 - c. Connection to a mains voltage other than the voltage specified on the type plate.
 - d. Unauthorized modifications.
 - e. Repairs carried out by third parties.
 - f. Careless transport, i.e. without suitable packaging materials or protection.
7. No claims may be made under this warranty for:
 - a. Losses incurred during transport.

- b. The removal or changing of the appliance's serial number.
8. The guarantee does not cover power cords, lamps or glass parts.
9. No claims can be submitted under this warranty for damage other than the repair or replacement of faulty parts. The importer can never be held responsible for any consequential loss or damage or any other consequences, resulting either directly or indirectly from the appliance supplied by the importer.
10. In case of claims under guarantee you can contact your dealer where the appliance is purchased. Bestron offers you also the possibility to send the appliance directly to our Service Department. Do not send your appliance without consulting us. The package may be refused and any any costs will be for your account. Please contact the Service Department and they will tell you how to pack and send the appliance.
11. This appliance is not suitable for professional use.

SERVICE

If a fault should occur please contact the BESTRON service department:

THE NETHERLANDS:

BESTRON NEDERLAND BV
Bestron Service
Moeskampweg 20
5222 AW 's-Hertogenbosch

Tel: +31 (0) 73 - 623 11 21
Fax: +31 (0) 73 - 621 23 96
Email: info@bestron.com
Internet: www.bestron.com

CE DECLARATION OF CONFORMITY


This product conforms to the essential requirements of the following EU safety directives:

- EMC Directive 2004/108/EC
- Low Voltage Directive 2006/95/EC

N. van de Vossenberg
Quality control

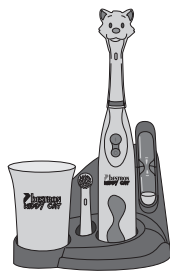BESTRON®

DSA3500C
v 030711-08