dexford

SC 200

Bluetooth 4.0 Health Scale

English

Compatible with :
iPhone 4s, iPhone 5, iPhone 5s,
iPhone 5c or later
iPod touch (5th generation)
iPad3, iPad 4, iPad air and iPad mini

Required OS: iOS6 or above

Compatible with :
Samsung Galaxy S4 (V4.2.2 or above)
Samsung Note 3 (V4.3)
Samsung Note 2 (V4.3)
HTC One (V4.2.2 or above)
HTC Butterfly (V4.2.2 or above)
Other device (V4.3)

Specification

- ◆ Bluetooth 4.0 compatible
- ◆ Measure body fat, hydration, muscle and bone percentage and record into memory and display on LCD Screen.
- ◆ Display BMI parameter
- ◆ High precision strain gauge sensors system
- ◆ Tempered glass platform of high quality and modern design.
- ◆ Smart memory of up to 10 persons' information
- ◆ Capacity: 2kg-150kg
- ◆ Division: 0.1kg
- ◆ Auto zero & Auto off
- ◆ Low battery and over load indication
- ◆ Free switch within three different units: kg/lb/st
- ◆ Powered by 4X1.5V

Keys


SET :Power on,menu select and confirm

▲ :Up

▼ :Down

LCD Display

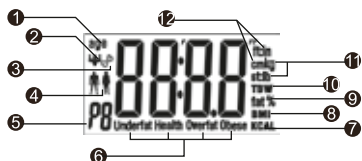- | | |
|-------------------------|-----------------|
| ① Age | ⑦ Calorie |
| ② Muscle | ⑧ BMI parameter |
| ③ Bone | ⑨ Fat% |
| ④ Gender | ⑩ Hydration% |
| ⑤ User/Memory/Code | ⑪ Weight Unit |
| ⑥ Fat status indication | ⑫ Height Unit |

Installation

Prior to first use, Download and install "My Way Fit" application to your IOS device. (Bluetooth 4.0 capabilities, e.g. iPhone 4S, iPhone 5, Newipad and iPadMini)

QUICK START GUIDE

Bluetooth connection:

- 1.Connect iOS device via Bluetooth (Turn Bluetooth "On" under the "Settings" menu on the iOS device)
- 2.Open "MyWay Fit" application, and then press the Weight icon on Measure screen.
- 3.Tag to go to Setting page.For Connection,Tag the "Health Scale" to Connection page for device pairing.
- 4.Turn on the scale,the APPs will show searched device on the display.
- 5.Tag the device for pairing.If pairing is successful, it will display "Paired and Connected".
- 6.Tag Back to return Scale Menu.
- 7.When the measurement is complete, the weight, body fat,hydration,muscle,and bone percentage will display on Scale & save on iOS application screen.

After that check the measurement history on History screen or chart screen by iOS application.

For Weight only


- Step on the scale gently and it will get ready for use when the screen shows "0.0kg" after full display,and you can press the scale back key to switch within different units(kg/lb/st)

- Step on the scale gently, position your feet evenly near the center and keep your body still and erect until the scale stabilize your weight.The scale will be auto off no further using.

For weight, body fat, hydration, muscle and bone percentage:

Press the "SET" key for seconds to set your personal I data.Press "▲" or "▼"key to select a number from P0 to P9 as your personal code when the code number flashes on the screen and once you reach your target number,please press "SET" to confirm.

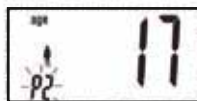
- Press "▲" or "▼" to choose "♂" Male or "♀" Female, when the gender symbol flashes on the screen and then press "SET" to confirm.






- Press "▲" or "▼" to reach your target age (10-99), when "age" and the number flash on the screen and the press "SET" to confirm.

- Press "▲" or "▼" to reach your target height (60-240cm) when the number and "cm" flash on the screen and then press "SET" to confirm.

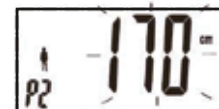
- Step on the scale gently and keep your bare feet in touch with the electrodes. Stand still and erect to get your weight measured and your body fat, hydration, muscle and bone percentage will be displayed respectively on the screen within seconds. These results will be repeated for two times and the scale will be auto off if no further using.

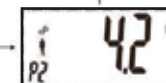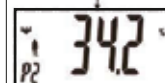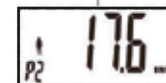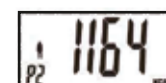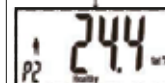
- The personal data will be stored once you set. Just confirm your code when you use the scale again and your other information will be showed to there is no need to program them again.

Caution

- 1.Personals with pacemakers should not use this product.
- 2.Personals with metal implants will be unable to obtain accurate readings.
- 3.To ensure accuracy, we suggest using scale at the same time every night (19:00-21:00)
- 4.The "Err" on scale display indicates overload. "Err2" display indicates fat percentage measurement error.
- 5.If body fat percentage is less than 5%, indicator displays "FATH". If more than 50% indicator displays "FATH". If impedance measuring makes an error, indicator will display "Err2".

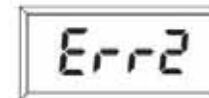- 6.The "ERR" display on iOS application screen indicates the measurement error.

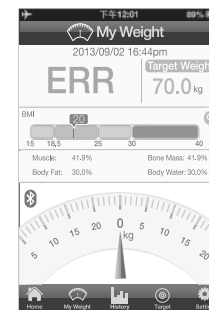Note

- 1.You should place scale on clear, flat surface before use.
- 2.Do not place scale on wet, severely hot or extremely cold surfaces.
- 3.Do not jump or stop on scale. Do not attempt to disassemble the scale.
- 4.To avoid falling off scale, do not stand on the scale to avoid glass and electrode wear out.
- 5.Do not wear shoes when standing on the scale to avoid glass and electrode wear out.
- 6.Never submerge the unit in water. Use alcohol to clean the electrodes and glass cleaner to keep them shiny, applied to a cloth first; avoid soaps.
- 7.For home use only.Not for professional use .Body fat percentage and hydration percentage is reference only. You should consult a doctor when you undertake any diet or exercise program.
- 8.Check battery if scale malfunction. If necessary, install a new battery.For other problems, please contact your local franchiser or contact our company for repair.

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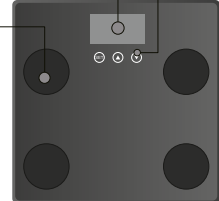
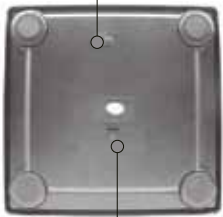
Service hotline

In the case of technical problems, contact our Service hotline. Switzerland: Tel. 0900 00 1675 (national charges, Swisscom at time of going to print: CHF 2.60/min). In the case of claims under the terms of guarantee, contact your sales outlet.

Guarantee

Dexford equipment is produced and tested according to the latest production methods. The implementation of carefully chosen materials and highly developed technologies ensure trouble-free functioning and a long service life. The terms of guarantee do not apply to the batteries or power packs used in the products. The period of guarantee is 24 months from the date of purchase. All deficiencies related to material or manufacturing errors within the period of guarantee will be redressed free of charge. Rights to claims under the terms of guarantee are annulled following tampering by the purchaser or third parties. Damage caused as the result of improper handling or operation, normal wear and tear, incorrect positioning or storage, improper connection or installation or Acts of God and other external influences are excluded from the terms of guarantee. In the case of complaints, we reserve the right to repair defective parts, replace them or replace the entire device. Replaced parts or devices become our property. Rights to compensation in the case of damage are excluded where there is no evidence of intent or gross negligence by the manufacturer. If your device does show signs of a defect within the period of guarantee, please contact the sales outlet where you purchased the Dexford device, producing the purchase receipt as evidence. All claims under the terms of guarantee in accordance with this agreement can only be asserted at the sales outlet. No claims under the terms of guarantee can be asserted after a period of two years from the date of purchase and hand-over of the product.

Disposal

In order to dispose of your device, take it to a collection point provided by your local public waste authorities (e.g. recycling centre). According to the laws on the disposal of electronic and electrical devices, owners are obliged to dispose of old electronic and electrical devices in a separate waste container. The adjacent symbol indicates that the device must not be disposed of in normal domestic waste! You are legally obliged to dispose of power packs and batteries at the point of sale or in the corresponding containers provided at collection points by local public waste authorities. Disposal is free of charge. The symbols indicate that the batteries must not be disposed of in normal domestic waste and that they must be brought to collection points provided by local public waste authorities. Packaging materials must be disposed of according to local regulations.

Declaration of Conformity

 This device fulfils the requirements stipulated in the Directive 1999/5/EC on radio equipment and telecommunications terminal equipment and the mutual recognition of their conformity. Conformity with the above mentioned directive is confirmed by the CE mark on the device. To view the complete Declaration of Conformity, please refer to the free download available on our web site www.dexford.com.

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