# INVENTUM

SINCE 1908

## MG 160

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# ENGLISH

Legislation requires that all electrical and electronic equipment must be collected for reuse and recycling. Electrical and electronic equipment marked with the symbol indicating separate collection of such equipment must be returned to a municipal waste collection point.


The device has double protective insulation and therefore complies with protection class 2.

## DESCRIPTION OF THE APPLIANCE

1. ABOUT MASSAGE
2. SAFETY INSTRUCTIONS
3. BEFORE FIRST USE
4. MASSAGE
5. USE IN CAR
6. PROBLEMS & SOLUTIONS
7. CLEANING & MAINTENANCE
8. TECHNICAL DATA  
GUARANTEE & AFTER SALES SERVICE

PAGE 16  
PAGE 16  
PAGE 16  
PAGE 17  
PAGE 17  
PAGE 18  
PAGE 18  
PAGE 19  
PAGE 19  
PAGE 22

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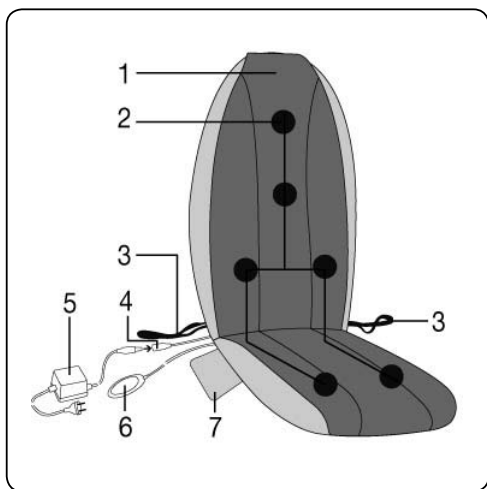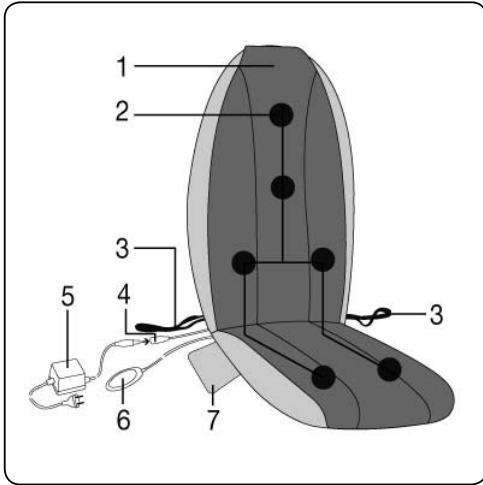
1. Protective cover
2. 6 massage motors
3. Fastening straps
4. Connecting coupling
5. Power supply unit with plug connector
6. Hand switch
7. Storage bag
8. Car adapter (not shown)

## 1 | ABOUT MASSAGE

You can have your back and thighs massaged while sitting with this electrical massage device. Massages can be both relaxing and stimulating, and are often used in the treatment of stiff or tense muscles, pains and fatigue. This device massages with the aid of 6 vibration motors and offers a wide range of setting options. The adapter supplied enables the seat cover to also be used in the car. The massage seat cover has to be used when placed on a suitable seat.

The device is only for private use and is not intended for medical or commercial purposes.

## 2 | SAFETY INSTRUCTIONS

- Please read these instructions before operating the appliance and retain for future use.
- The appliance is intended solely for domestic use within the home.
- This device is solely intended for massaging the backs and legs of human beings. It cannot replace medical treatment. Do not use the massage device if any of the following warnings apply to you. Please consult your doctor if you are unsure whether or not the massage device is suitable for you.
- Always connect the device as indicated with the voltage printed on the bottom of the appliance.
- Just like any other electrical appliance, this massage device must be handled with due care and attention in order to prevent the risk of electric shock.
- Unplug the mains cable immediately after use.
- Prevent any contact with water (except when cleaning with a slightly damp cloth)! Water must never get inside the device.
- Never immerse any part of the main body in water or any other liquid to protect against electrical hazards.

- Never use the device in the bath, in the shower, in a swimming pool or over a basin filled with water. However, if water has got inside the casing, unplug the device immediately and contact your electrical dealer or customer service.
- The device and in particular the mains cable should be checked for any damage before each use.
- Never plug the device into the mains power supply if the cable is damaged. Contact your electrical dealer or customer service.
- If a fault should occur during use, switch off the device immediately and unplug it from the socket.
- The device is only intended for individual use, and not for medical or commercial use.
- Do not use the device continuously for longer than 15 minutes to avoid overheating. Let the device cool down before using it again.
- Never cover the device when it is switched on (with blanket, pillow etc.).
- Never operate the device unsupervised, especially if there are children in the vicinity.
- Children (under the age of 15), helpless or disabled persons must use the device only under adult supervision or must be familiarised with the device beforehand.
- Do not use with animals or pets.
- The device has a warm surface. Persons who are sensitive to heat must be careful when using the device.
- Persons with cardiac pacemakers are requested to consult their doctor before using the massage device and never to massage the region of the heart.
- Do not use the device while operating a car or machine.
- Never use the device on swollen, burned or injured areas of skin. The following should also not be massaged: bones (e.g. joints, spine), the head or other sensitive areas of the body.
- Do not carry out a massage if you have a cold with a high temperature, varicose veins, thrombosis, vein inflammation, jaundice, diabetes, pregnancy, nerve

disorders (e.g. sciatica) or acute inflammation. In such cases, consult your doctor.

- Do not use the massage device in bed. You must never fall asleep during use.
- Always return the appliance after a malfunction, or if it has been damaged in any manner to the Customer Service Department for examination, repair or adjustment as special purpose tools are required.
- An electrical appliance is not a toy, and should therefore always be placed outside the reach of children. Children are not always able to understand potential risks. Teach children a responsible handling of electrical appliances.
- If you decide never to use the appliance again, we advise you to dispose of it in special containers or at a municipal waste collection point.
- Never use the appliance if it is not working properly, or if it has been damaged. Any intervention other than cleaning and normal maintenance by the customer must be carried out by an approved service centre.
- The zip fastener on the Shiatsu massage device must not be opened. It is only fitted for reasons relating to the manufacturing of the device.
- Never insert any objects into openings on the device, and never insert anything into the rotating parts. Make sure that all moving parts can move freely at all times.
- Make sure that no hair can get caught in the device.
- Avoid contact with sharp or pointed objects.
- Never sit or stand on the device with all your weight and do not place objects on the device.
- If the appliance is used for other purposes than intended, or is it not handled in accordance with the instruction manual, the full responsibility for any consequences will rest with the user. Any damages to the product or other things are not covered by the warranty.

### 3 | BEFORE FIRST USE

Before you use the appliance for the first time, please do as follows:

- Carefully unpack the appliance and remove all the packaging material.
- Keep the material (plastic bags and cardboard) out of reach of children.
- Check the appliance after unpacking for any damage, possible from transportation.
- Secure the protective cover on the surface of the massage seat cover using the Velcro strips where the seat joins with the backrest, on the reverse side at the front under the seat surface and at the top in the shoulder area. It is also possible to massage without the protective cover; this is generally found to be "harder".
- Position the seat cover upright on a suitable seat with seating surface and backrest (chair, sofa, or similar). Ensure that the seat cover is supported over its entire area and has a sufficiently high backrest.
- Use the flexible fastening straps to secure the massage seat cover.
- Connect the power supply unit and the seat cover with the plug connector.
- Connect your device to the mains by plugging the power supply unit into the mains outlet.

### 4 | MASSAGE

The massage must be pleasant and relaxing at all times. Stop the massage or change the position of the massage device if you find it is painful or unpleasant.

1. ON/OFF button with control LED  
*Switch on and off*
2. INTENSITY button with control LEDs  
*Select massage intensity*
3. HEAT button with control LED  
*Switch the heat function on and off*
4. MASSAGE ZONE button with control LEDs  
*Switch the various massage zones on and off*
5. PROGRAM button with control LEDs  
*Select massage program*
6. TIMER button with control LEDs  
*Select switch-off times*
7. Connection cable to the device

The control LEDs indicate which functions are currently selected and which massage zones are activated.

- Sit comfortably on the seat cover. Take care that you sit centered and your body is in good even contact with the seat cover.
- Switch on the massage device with the "ON/OFF" button. The massage intensity is preset to level 1.
- Press the "PROGRAM" button or at least one of the "MASSAGE ZONES" buttons to start the massage.
- If you wish to have a long massage with individual massage zones or all of them, select the massage zones with their respective buttons. The massage zones activated are displayed with the LEDs.
- The automatic switch-off function is activated by starting the massage and is preset to 20 minutes. You can set the automatic switch-off function to 10, 20 or 30 minutes by pressing the "TIMER" button. Please note that changing the timer setting makes the timer start to count again from the beginning. To avoid overheating the device, select the maximum massage duration at the beginning of the massage and do not change it during the massage.
- If you would like a massage with one of the various programs, select the required program with the "PROGRAM" button.  
The 5 massage programs are characterized by different sequences of massage zones and massage intensities.
- Press the "INTENSITY" button to select of the 3 possible massage intensities.
- Press the "HEAT" button to switch the heat function in the lumbar region on and off.
- Do not use the massage device for longer than 30 minutes. The massage can be finished at any time by pressing the "ON/OFF" button.

- After the massage is finished, connect the power supply unit plug from the mains.

**Note**

- You can place the hand switch in the storage case during and after use.
- Always wear light clothing when receiving a massage from the device.
- To receive a more gentle massage, insert a folded towel between the massage device and your body.
- Do not use the massage device for longer than 30 minutes. Otherwise, using the device for longer can over-stimulate the muscles and cause tension instead of relaxation.
- Do not use the massage device just before going to bed. The massage can also have a stimulating effect and may cause difficulties in getting off to sleep.

**5 | USE IN CAR**

The massage seat cover can also be used in the car with the car adapter supplied. Plug the car adapter connector plug into the connection jack on the device. Plug the car adapter into the cigarette lighter jack or an alternative 12 V jack (e.g. trunk). Operation proceeds as described above.

**WARNING! Never have a massage with the massage seat cover while you are driving a car yourself.**

The massage, as well as the device setting, distracts you from the road traffic. Your attention to the road traffic must never be impaired by the device. Therefore only use this device as a passenger or if you take a break and the vehicle is stationary. Please note that if the engine is off the battery can discharge under certain circumstances. For some vehicles the 12V jack only works when the engine is on. A massage with this seat cover can be very relaxing and regenerating while taking breaks.

**6 | PROBLEMS & SOLUTIONS**

Problem	Cause	Solution
The device is very loud during massage.	The seat cover is not fully in contact with the base.	Ensure that the seat cover is supported over its entire area.
Heat function cannot be activated.	Overheating protection has been tripped.	Contact your authorized dealer or the service address.
Massage stopped.	The preset maximum period of operation has elapsed.	Let the device cool down for at least 15 minutes before using again. Switch the device on again.
	The connection to the mains has been interrupted.	Check the plug connection and the contact of the mains plug with the mains supply.
No massage.	Appliance not connected to the mains power supply.	Make all connections and switch on the device.
	Device not switched on.	Switch the device on with "ON/OFF". Press "PROGRAM" or "MASSAGE ZONES" and select the setting.
	Overheating protection tripped, device is damaged.	Contact your authorized dealer or the service address.

**7 | CLEANING & MAINTENANCE**

Always disconnect the device before cleaning it. Only clean the appliance in the manner specified. Fluid must never be allowed to enter into the device or accessories. Do not use the device again until it is completely dry.

First unplug the device and allow it to cool off. Then you can clean it with a cloth which is no more than lightly moistened. You can also use a mild detergent.

**Never immerse any part of the device in water or any other liquid to protect against electrical hazards.**

The removable cover is machinewashable at 40 °C. Please note the cleaning and care symbols on the label which is sewn into the cover.

Handle the device with care, protect the device from impact with hard objects and do not drop it.

## 8 | TECHNICAL DATA

<i>Model</i>	MG160
<i>Dimensions</i>	47,5 x 107 x 7 cm
<i>Weight</i>	ca. 2,1 kg
<i>Rated output</i>	18 Watt
<i>Power supply unit</i>	Input: 230 V/ AC/ 50 Hz Output: 12 V DC/ 1,5 A
<i>Accessories</i>	Massage seat cover, protective cover, power supply unit, car adapter and instruction manual.