SHIATSU MASSAGE - MASSAGE SHIATSU

**INVENTUM**

SINCE 1908

**MG 190**

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# ENGLISH

Legislation requires that all electrical and electronic equipment must be collected for reuse and recycling. Electrical and electronic equipment marked with the symbol indicating separate collection of such equipment must be returned to a municipal waste collection point.



The device has double protective insulation and therefore complies with protection class 2.

## DESCRIPTION OF THE APPLIANCE

1. ABOUT SHIATSU MASSAGE
  2. SAFETY INSTRUCTIONS
  3. BEFORE FIRST USE
  4. SHIATSU MASSAGE
  5. PROBLEMS & SOLUTIONS
  6. CLEANING & MAINTENANCE
  7. TECHNICAL DATA
- GUARANTEE & AFTER SALES SERVICE

PAGE 16  
PAGE 16  
PAGE 16  
PAGE 17  
PAGE 18  
PAGE 18  
PAGE 18  
PAGE 19  
PAGE 22

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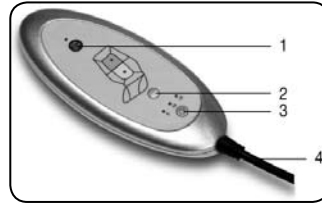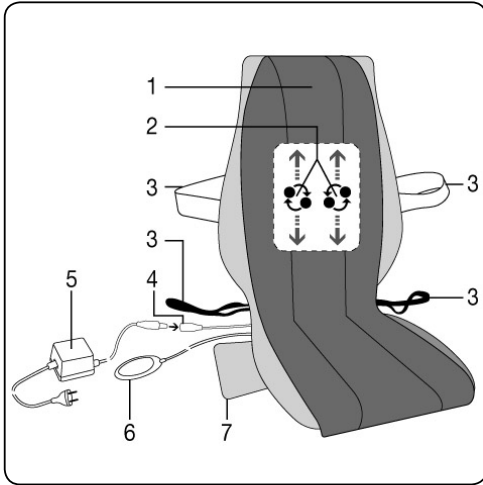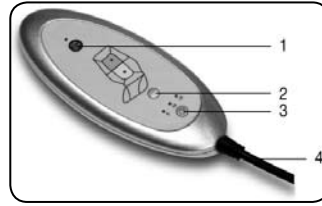
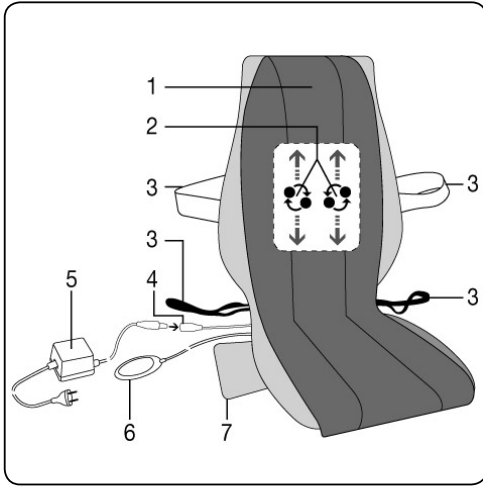
1. ON/OFF button with control LED
2. MESSAGE button with control LEDs
3. TIMER button with control LEDs
4. Connection cable to the device

1. Protective cover
2. 4 massage heads, rotating in pairs with vertical movement
3. Fastening straps
4. Connecting coupling
5. Power supply unit with plug connector and power cable
6. Hand switch
7. Storage bag

## 1 | ABOUT SHIATSU-MASSAGE

Shiatsu is a form of body massage which was originally developed in Japan around 100 years ago. It is based on ideas anchored in traditional Chinese medicine, in particular the meridian system, i.e. the paths in the human body along which energy is conveyed.

The aim of any Shiatsu massage is to promote the physical, emotional and mental wellbeing of the person being treated. To do this, energy blocks and bottlenecks in the meridians need to be dissolved and the self-regulating ability of the body needs to be stimulated.

Shiatsu therapists achieve this by exercising pressure in flowing movements along the meridians. Although Shiatsu can be literally translated as “finger pressure”, practitioners also use the heels of their hands, their elbows and knees as well as their fingers. In the process, the therapist will use his or her full body weight in order to generate the required pressure. Your Shiatsu massage device has rotating massage heads which imitate the pressure and kneading motions of a Shiatsu massage.

You can have your back massaged while sitting with this electrical Shiatsu massage device. Massages can be both relaxing and stimulating, and are often used in the treatment of stiff or tense muscles, pains and fatigue. This device, with its vertically moving and rotating massage heads, offers a powerful, intensive deep massage over the entire back.

You have to place the Shiatsu seat cover on a suitable seat pad for this purpose. The device is only for private use and is not intended for medical or commercial purposes.

In addition, the Shiatsu massage device offers the following advantages:

- Adjustable massage zones (upper, lower or the entire back),
- Timer function with selectable switch-off after 5, 10 or 15 minutes,
- Removable and washable cover,
- Disconnectable power supply unit.

## 2 | SAFETY INSTRUCTIONS

- **Please read these instructions before operating the appliance and retain for future use.**
- The appliance is intended solely for domestic use within the home.
- This device is solely intended for massaging the backs of human beings. It cannot replace medical treatment. Do not use the massage device if any of the following warnings apply to you. Please consult your doctor if you are unsure whether or not the massage device is suitable for you.
- Always connect the device as indicated with the voltage printed on the bottom of the appliance.
- Just like any other electrical appliance, this massage device must be handled with due care and attention in order to prevent the risk of electric shock.
- Unplug the mains cable immediately after use.
- Prevent any contact with water (except when cleaning with a slightly damp cloth)! Water must never get inside the device.
- **Never immerse any part of the main body in water or any other liquid to protect against electrical hazards.**
- Never use the device in the bath, in the shower, in a swimming pool or over a basin filled with water.

However, if water has got inside the casing, unplug the device immediately and contact your electrical dealer or customer service.

- The device and in particular the mains cable should be checked for any damage before each use.
- Never plug the device into the mains power supply if the cable is damaged. Contact your electrical dealer or customer service.
- If a fault should occur during use, switch off the device immediately and unplug it from the socket.
- The device is only intended for individual use, and not for medical or commercial use.
- Do not use the device continuously for longer than 15 minutes to avoid overheating. Let the device cool down before using it again.
- Never cover the device when it is switched on (with blanket, pillow etc.).
- Never operate the device unsupervised, especially if there are children in the vicinity.
- Children (under the age of 15), helpless or disabled persons must use the device only under adult supervision or must be familiarised with the device beforehand.
- Do not use with animals or pets.
- Persons with cardiac pacemakers are requested to consult their doctor before using the massage device and never to massage the region of the heart.
- Never use the device on swollen, burned or injured areas of skin. The following should also not be massaged: bones (e.g. joints, spine), the head or other sensitive areas of the body.
- Do not carry out a massage if you have a cold with a high temperature, varicose veins, thrombosis, vein inflammation, jaundice, diabetes, pregnancy, nerve disorders (e.g. sciatica) or acute inflammation. In such cases, consult your doctor.
- Do not use the massage device in bed. You must never fall asleep during use.
- Always return the appliance after a malfunction, or if it has been damaged in any manner to the Customer Service Department for examination, repair or adjustment as special purpose tools are required.
- An electrical appliance is not a toy, and should therefore always be placed outside the reach of children. Children are not always able to understand potential risks. Teach children a responsible handling of electrical appliances.
- If you decide never to use the appliance again, we advise you to dispose of it in special containers or at a municipal waste collection point.
- Never use the appliance if it is not working properly, or if it has been damaged. Any intervention other than cleaning and normal maintenance by the customer must be carried out by an approved service centre.
- The zip fastener on the Shiatsu massage device must not be opened. It is only fitted for reasons relating to the manufacturing of the device.
- Never insert any objects into openings on the device, and never insert anything into the rotating parts. Make sure that all moving parts can move freely at all times.
- Make sure that no hair can get caught in the device.
- Avoid contact with sharp or pointed objects.
- Never sit or stand on the device with all your weight and do not place objects on the device.
- If the appliance is used for other purposes than intended, or is it not handled in accordance with the instruction manual, the full responsibility for any

consequences will rest with the user. Any damages to the product or other things are not covered by the warranty.

### 3 | BEFORE FIRST USE

Before you use the appliance for the first time, please do as follows:

- Carefully unpack the appliance and remove all the packaging material.
- **IMPORTANT: Deactivate the transport lock on the massagemechanism before starting-up for the first time. Remove the screw on the reverse side of the seat cover using the tool supplied. The Shiatsu seat cover will be irreparably damaged and unusable if the transport lock screw is not removed and the Shiatsu seat cover is operated immediately. No guarantee is accepted in this case.**
- Insert the transport fixation screw each time the device is transported again. The device may be irreparably damaged each time it is transported without the transport lock.
- Keep the material (plastic bags and cardboard) out of reach of children.
- Check the appliance after unpacking for any damage, possible from transportation.
- Secure the protective cover on the surface of the Shiatsu seat cover using the Velcro strips where the seat joins with the backrest, on the reverse side at the front under the seat surface and at the top in the shoulder area. It is also possible to massage without the protective cover; this is generally found to be "harder".
- Position the seat cover upright on a suitable seating arrangement with seating surface and backrest (chair, sofa, or similar). Ensure that the seat cover is supported over its entire area and has a sufficiently high backrest. As a result of the solid back cover, the Shiatsu seat cover can also be operated with small backrests.
- Use the flexible fastening straps and Velcro strips to secure the Shiatsu seat cover.
- Connect the power supply unit and the seat cover with the plug connector.
- Connect your device to the mains by plugging the power supply unit into the mains outlet.

## 4 | SHIATSU MASSAGE

The massage must be pleasant and relaxing at all times. Stop the massage or change the position of the massage device if you find it is painful or unpleasant.

- Sit comfortably and upright on the seat cover. Ensure you sit centered and such that the massage heads can move to the left and right of your spinal column.
- Switch on the massage device with the ON/OFF button. Select the region you wish to massage with the MESSAGE button.

1. Entire back
2. Lower back
3. Upper back
4. Stop

The LED lamp indicates the massage region set.

- The Shiatsu massage unit with rotating massage heads slowly moves upwards from the starting position in the lowest region of the back and moves up and down in the massage region set. Application only at a single point is not possible, as the rotational movement and the displacement up and down are coupled.

- The automatic switch-off function is activated with selection of the massage region and is preset to 15 minutes. You can set the automatic switch-off function to 5, 10 or 15 minutes by pressing the TIMER button.

Please note that changing the timer setting makes the timer start to count again from the beginning.

To avoid overheating the device, select the maximum massage duration at the beginning of the massage and do not change it again during the massage.

- Do not use the massage device for longer than 15 minutes. The massage can be finished at any time by pressing the ON/OFF button. Please note that on switching off, the Shiatsu massage unit firstly returns to the starting position at the lowest region of the back. This is indicated by the ON/OFF LED blinking. Never interrupt the return to the starting position. Once this position is attained, the device switches off.
- If the device is stopped using the MESSAGE button, it also returns to the starting position after 15 minutes and then switches off.
- DO NOT disconnect the device from the mains supply during operation. Always switch it off only once it has reached the starting position at the bottom.
- Disconnect the power supply unit plug from the mains.

### Note

- You can place the hand switch in the storage case during and after use.
- Always wear light clothing when receiving a massage from the device.

- To receive a more gentle massage, insert a folded towel between the massage device and your body.
- Do not use the massage device for longer than 15 minutes. Otherwise, using the device for longer can over-stimulate the muscles and cause tension instead of relaxation.
- Do not use the massage device just before going to bed. The massage can also have a stimulating effect and may cause difficulties in getting off to sleep.

## 5 | PROBLEMS & SOLUTIONS

Problem	Cause	Solution
Massage heads rotating at reduced speed.	The load on the massage heads is too great.	Reduce the contact pressure on the head.
Massage heads only rotate in the upper or lower region.	The massage has been activated for the limited region.	Press the MESSAGE button to change/expand the region.
Massage heads are not moving.	Appliance is not connected to the power supply.	Connect the plug to the mains supply and switch on the appliance.
	Device not switched on.	Switch the device on with ON/OFF. Select the massage region.
	Transport lock has not been removed, the massage mechanism is damaged.	Contact your authorized dealer or the service address.
	Overheating protection tripped, device is damaged.	Contact your authorized dealer or the service address.

## 6 | CLEANING & MAINTENANCE

Always disconnect the device before cleaning it. Only clean the appliance in the manner specified. Fluid must never be allowed to enter into the device or accessories. Do not use the device again until it is completely dry.

First unplug the device and allow it to cool off. Then you can clean it with a cloth which is no more than lightly moistened. You can also use a mild detergent.

**Never immerse any part of the device in water or any other liquid to protect against electrical hazards.**

The removable cover is machinewashable at 40°C. Please note the cleaning and care symbols on the label which is sewn into the cover.

Handle the device with care, protect the device from impact with hard objects and do not drop it.

## 7 | TECHNICAL DATA

<i>Model</i>	MG190
<i>Dimensions</i>	50 x 107 x 7,5 cm
<i>Weight</i>	ca. 3,8 kg
<i>Rated output</i>	24 Watt
<i>Power supply unit</i>	Input: 230 V/ AC/ 50 Hz Output: 12 V DC/ 2,0 A