PERSONENWEEGSCHAAL - PERSONENWAAGE - PÈSE PERSONNE - SCALE

INVENTUM

SINCE 1908

PW 535 G

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

ENGLISH

Legislation requires that all electrical and electronic equipment must be collected for reuse and recycling. Electrical and electronic equipment marked with the symbol indicating separate collection of such equipment must be returned to a municipal waste collection point.

DESCRIPTION OF THE APPLIANCE	PAGE 10
1. SAFETY INSTRUCTIONS	PAGE 10
2. WEIGHING CORRECTLY	PAGE 10
3. BEFORE FIRST USE	PAGE 11
4. BATTERIES	PAGE 11
5. WEIGHING	PAGE 11
6. CONVERT KG/LB	PAGE 11
7. CLEANING & MAINTENANCE	PAGE 11
GUARANTEE & AFTER SALES SERVICE	PAGE 14

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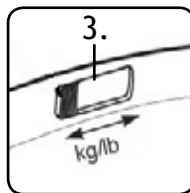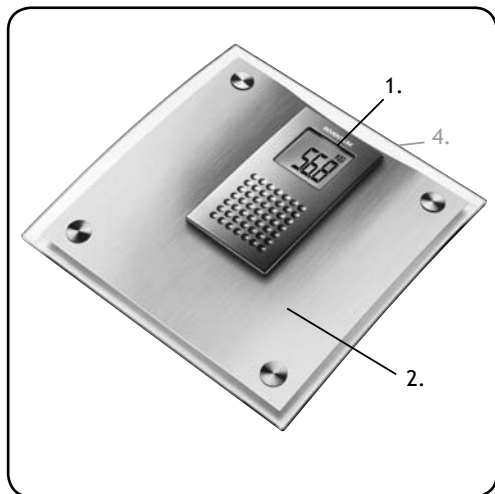
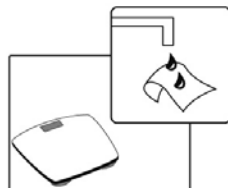
fig. 4

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

7 | NETTOYAGE & ENTRETIEN

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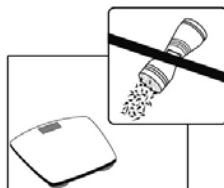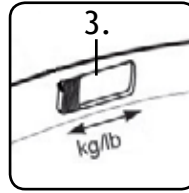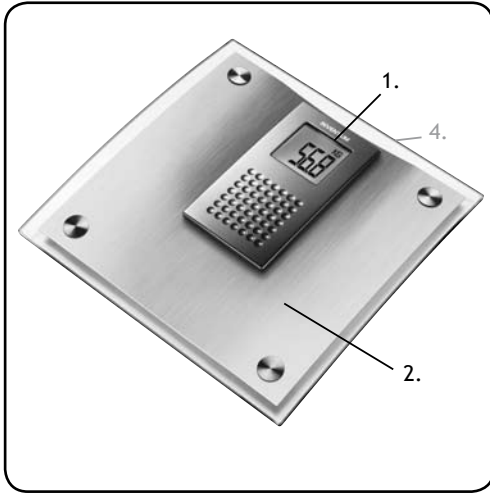1. Display
2. Safety glass surface
3. Convert button (kg/lb)
4. Battery compartment (underneath)

Maximum capacity: 150 kilograms

Minimum capacity: 10 kilograms

Indication per 100 grams

1 | SAFETY INSTRUCTIONS

- Please read these instructions before operating the appliance and retain for future use.
- The appliance is intended solely for domestic use within the home.
- Always use batteries which corresponds with the voltage printed on the bottom of the appliance.
- If possible, always weigh yourself at the same time of day (ideally in the morning), after going to the toilet, on an empty stomach and without clothes to achieve comparable results.
- Minor deviation are normal. For example when you drink 1 glass of water a deviation of approx. 300 grams can occur.
- Do not stand on the scale with wet feet and do not step on the scale when its surface is wet - danger of slipping!
- Protect the scale against hard knocks, moisture, dust, chemicals, major temperature fluctuations and heat sources which are too close (stove, heating radiators).
- If you want to store the appliance for a long period, please keep it in the original box.
- Always use the appliance on a secure, dry, flat surface.
- Never use the appliance if it is not working properly, or if it has been damaged. Any intervention other than cleaning and normal maintenance by the customer must be carried out by an approved service centre.
- If you decide never to use the appliance again, we advise you to dispose of it in special containers or at a municipal waste collection point.
- Never immerse the appliance water or any other liquids.
- If the appliance is used for other purposes than intended, or is it not handled in accordance with the instruction manual, the full responsibility for any

consequences will rest with the user. Any damages to the product or other things are not covered by the warranty.

- Spent batteries and rechargeable batteries do not constitute normal household waste! They are considered to be toxic waste and, as such, should be disposed of in special containers, toxic waste collection points or brought to electrical goods dealers. You are legally obliged to dispose of the batteries. **Note:** Batteries containing pollutant substances are marked as follows: Pb = Battery contains lead, Cd = Battery contains cadmium, Hg = Battery contains mercury.

2 | WEIGHING CORRECTLY

Weight problems, particularly preponderance, can cause or accelerate diseases. Therefore regular checks are advisable. The body weight depends both on fixed and influenceable factors. To the first belong age, sex as well as hereditary-conditioned influences.

To the influenceable factors belong the kind and quantity of the nourishment, physical and physiological activity, as well as the process and the treatment of illness-related processes. Of special importance is the temporal change of your body weight. In order to seize this, you should weigh regularly. The maximum capacity of the scale is max. 150 kg (330 lb). During the weight measurement the results are indicated in 100 g steps (0.2 lb).

3 | BEFORE FIRST USE

Before you use the appliance for the first time, please do as follows: Carefully unpack the appliance and remove all the packaging material. Keep the material (plastic bags and cardboard) out of reach of children. Check the appliance after unpacking for any visible damage, possible from transportation.


4 | BATTERIES

If present, pull the battery insulating strip off the battery compartment cover or remove the batteries protective film and insert the batteries according to the polarity.

If the scale fails to operate, remove the batteries completely and insert them again. Your scale is equipped with a replace battery indicator.

When the scale is operated with batteries which are too weak, Lo appears in the display field and the scale automatically switches off. The battery must be replaced in this case (3 x 1,5 V, AAA, Mikro). The used, completely drained standard and rechargeable batteries must be disposed of in specially marked collection containers, at toxic waste collection points or electrical product retailers.

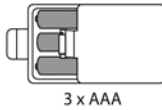5 | WEIGHING

Place the scale on a firm level floor (no carpet); a firm floor covering is required for correct measurement.

Press quickly and forcefully with your foot on the platform of your scale!

The entire display (Fig. 1) appears up to 00.0 (Fig. 2) as a self-test (tap-on technology).

Now the scale is ready to measure your weight.

Stand on the scale without moving and distribute your weight evenly on both legs.

The scale immediately begins to measure your weight (Fig. 3).

The scale switches off 10 seconds after you step down from it.

The measuring result remains visible during this time.

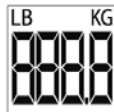
fig. 1

fig. 2

fig. 3


If the scale detects an error during weighing, Err appears in the display (Fig. 4). Make sure that you always switch on the scale first and wait for the display 00.0 (Fig. 2) before stepping on. Switch the scale on by tapping the scale surface once with your foot.



fig. 4

6 | CONVERT KG/LB

Measuring of the weight is possible in kg (kilograms) or in lb (pounds).

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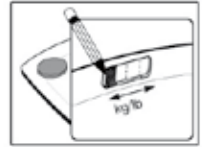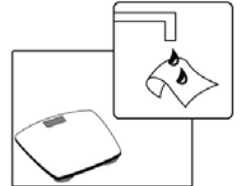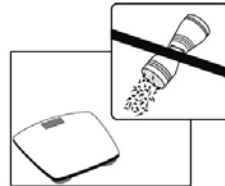GUARANTEE & AFTER SALES SERVICE

- This product is guaranteed for a period of 2 years from the date of purchase against mechanical and electrical defect.
- Guarantee can only be given if proof of purchase (e.g. till receipt) can be shown. Therefore keep it in a safe place. The guarantee period starts at date of purchase.
- Guarantee is not valid in case of:
 - normal wear;
 - damage through misuse use or when the instructions were not followed correctly;
 - lack of maintenance;
 - not following the operation- and safety instructions;
 - incompetent montage or repairs made by yourself or others;
 - using accessories and spare parts which are not suitable for the appliance; for assessment by our technical department.
- The guarantee is not valid for:
 - normal wear;
 - external cables;
 - batteries;
 - (interchangeable) signal- and controllamps for assessment by our technical department.

Return to your supplier or send the appliance, including a copy proof of purchase and description of the complaint, stamped to:

Inventum Holland BV
Generatorstraat 17
3903 LH VEENENDAAL
The Netherlands

or to:
new E.T.P. bvba
Passtraat 121
9100 Sint Niklaas
BELGIUM

- Pack it carefully (preferably in the original carton). Ensure unit is clean.
- Enclose your name and address and quote model number on all correspondence.

When the guarantee is no longer valid or when the period of guarantee has passed, we will charge any costs that flow out of an assessment and/or repairing.

On request we will supply a statement before repairing. Without requesting a statement, we assume that you, by sending the appliance, automatically order to repair the appliance.

A guarantee period of two months is valid for repairs done by our technical department and the used spare parts, counting from the repair date. When a malfunction occurs within this period which directly flows out of the repairs and the used spare parts, the consumer has the right of repairs without any costs.

Inventum Holland BV
Generatorstraat 17, 3903 LH , Veenendaal, Nederland
Tel.: 0318-542222, Fax.: 0318-542022, Mail: helpdesk@inventum.eu

Inventum Belgium BVBA
Lindenstraat 32, 3140, Keerbergen, België
Tel: 015-512653, Fax.: 015-529718, Mail: info@inventumbelgium.be

Website: www.inventum.eu