
LANAFORM®

Foot Massager

Massager

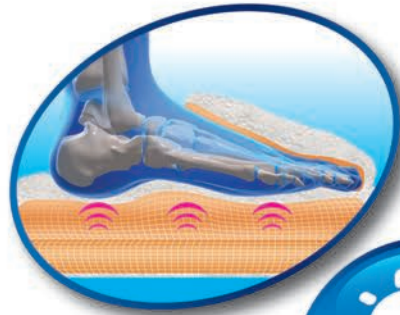Eases feet

Soulage les pieds

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INSTRUCTION MANUAL

Thank you for purchasing a LANAFORM® FOOT MASSAGER. The FOOT MASSAGER can be used to refresh tired feet through its gentle vibrations.

PLEASE READ ALL THE INSTRUCTIONS BEFORE USING YOUR FOOT MASSAGER, IN PARTICULAR THESE FEW FUNDAMENTAL SAFETY RULES:

- Only use this appliance for the purposes described in this manual.
- This appliance is not intended for use by persons (including children) whose physical, sensory or mental capabilities are impaired, or by persons lacking experience or knowledge, unless they have had the benefit of prior supervision and instruction in using the device from a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- Please do not use accessories which are not recommended by LANAFORM or which are not supplied with this appliance.
- Never use this appliance in a room where aerosol products (sprays) are used or in a room where oxygen is administered.
- Never obstruct the appliance's air vents and never place the device on an unstable surface such as a carpet or mattress. This could prevent air circulation. Make sure that the air vents are not clogged with hair, lint, dust, etc. Do not insert any objects into the air vents.
- Do not use this product before going to bed. Massage has a stimulating effect which may inhibit sleep.
- Do not use this massage appliance while you are asleep.
- Do not place the Foot Massager near a radiator or any other source of heat.
- Never use this appliance on a part of the body which is swollen or inflamed or if you have any form of skin rash.
- Never allow anything to fall into or insert any object (drawing pin or other metal clip) into any of the openings on this device.
- If the Foot Massager stops operating, have it inspected and repaired by the supplier or its after-sales department.
- Do not use this appliance under a cover or under a cushion. Excessive heat may cause a fire, electrocution of the individual or other injuries.
- This product is a non-professional massage appliance designed to relax tired muscles. Do not use this appliance as a replacement for medical care, it is intended exclusively for family use.

COMPONENTS

1. Battery compartment

INSTRUCTIONS FOR USE

- Insert the batteries via the zip fastener located at the back of your Foot Massager.
- Place 2 “LR20” batteries (not supplied) in the compartment provided for this purpose. Replace the vibrating unit into the slipper in the correct location, and close it.
- Insert your feet into the upper part of the slipper.
- To start the vibrations, press down gently towards the bottom of the slipper.

MAINTENANCE

1. Remove the vibrating unit before washing the slipper.
2. Do not wash the vibrating unit.
3. Wash the slipper at a maximum of 30°.
4. Dry the slipper at a maximum of 30°.

TECHNICAL CHARACTERISTICS

Voltage: 2x1.5 Volts – LR 20

COMPOSITION OF THE PRODUCT

35% acrylic - 65% polyester

ADVICE REGARDING THE DISPOSAL OF WASTE

- All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.
- When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law. Prior to disposal, please remove the battery and dispose of it at a collection point so that it can be recycled.

Please note: under no circumstances should spent batteries be disposed of with household waste!!!

LIMITED WARRANTY

LANAFORM® guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM® guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM® product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM®'s control.

LANAFORM® may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM® will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM®'s Service Centre®. Any attempt to maintain this product by a person other than LANAFORM®'s Service Centre will render this guarantee void.

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