LANAFORM®

Full Mass

Massage Device

Vibration
technology

Technologie
vibratoire

Triltechnologie

FULL MASS

INTRODUCTION

We would like to thank you for purchasing a «FULL MASS» by LANAFORM®. The «FULL MASS» is a quality massage appliance that uses vibration technology. Easy to handle and practical, the «FULL MASS» will allow you to take advantage of a massage at any time of the day regardless of where you are. Comfortable and sophisticated, the «FULL MASS» offers 6 different operating speeds and a pre-programmed mode that allows you to adapt the massage as you wish. This appliance will improve your well-being: reduced fatigue, muscular tension and stiffness. In addition, thanks to its specially researched design, the «FULL MASS» not only adapts to your belly, but also to your back, buttocks, thighs and calves.

PLEASE READ ALL THE INSTRUCTION PRIOR TO USING YOUR MASSAGER, IN PARTICULAR THESE FEW BASIC SAFETY INSTRUCTIONS

- Only use this appliance for the purposes described in this manual.
- This appliance is not intended for use by persons (including children) whose physical, sensorial or mental capabilities are impaired, or by persons lacking experience or knowledge, unless they are supervised or instructed in the use of the appliance first by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- Please do not use accessories which are not recommended by LANAFORM® or which are not supplied with this appliance.
- If the power lead is damaged, it must be replaced by a special or similar lead available from the supplier or after-sales department.
- Do not use this appliance if the plug is damaged, if it is not working correctly, if it has been dropped or if it is damaged or if it has fallen into water. If this is the case, have the appliance inspected and repaired by the supplier or its after-sales department.
- Do not transport this appliance by carrying it by its power lead and do not use this lead as a handle.
- Always disconnect the appliance from the mains after use or before cleaning it.
- Do not try to retrieve an appliance which has fallen into water; simply disconnect it from the mains immediately.
- An electrical appliance should never be left plugged in and unsupervised. Disconnect it from the mains if you are not using it.
- Keep the power lead away from hot surfaces.
- Never let anything fall into and never insert anything into one of the air vents.
- Do not use this appliance with the air vents obstructed.
- Never pour any liquid into one of the openings on this product.
- Do not use this appliance less than 2 hours before and after meals or before going to bed. This is because the stimulating effect of the massage can disturb the intake and digestion of food or delay falling asleep.
- Do not use this massage appliance while you are asleep.
- Never use this appliance in contact with water or in a humid environment. This could result in the risk of electric shock.
- This appliance is intended solely for indoor use.
- Never use this appliance in a room where aerosol products (sprays) are used or in a room where oxygen is administered.
- Do not use this appliance under a cover or under a cushion. Excessive heat may cause a fire, electrocution of the individual or other injuries.
- Never use this appliance in a room where the temperature exceeds 40°C.
- If the «FULL MASS» is stored in a room where the temperature is low, it is recommended that you wait several minutes before operating the appliance in a warmer environment. Indeed, water droplets could form caused by water condensation, which could alter the operation of your appliance.

- To disconnect the appliance, set the device to "OFF" mode \odot , and then unplug it.
- If you are worried about your health or if you are under medical treatment, consult a doctor before using this appliance.
- If you feel any pain while using this appliance, stop using it immediately and consult your doctor.
- Do not use this appliance for more than 10 minutes in the same position and on the same area of the body. Do not use for more than a total of 20 minutes. You must allow for a rest period of 30 minutes between each use.
- Never use this appliance on a part of the body that is swollen or inflamed or if you have any form of skin rash.
- Do not use this appliance if you are pregnant or during menstruation.
- This product is a non-professional massage appliance intended exclusively for family, designed to relax tired muscles. Do not use this appliance as a substitute for proper medical care.
- In case of malfunction, do not attempt to repair the "FULL MASS" yourself. It must be taken back to the supplier or its after-sales department.

INSTRUCTIONS FOR USE

1) Characteristics :

- Massage produced by helical vibratory movements.
- 6 massage intensities (= Manual mode).
- 5 massage areas (belly/back/buttocks/thighs/calves).
- 1 pre-programmed mode combining the different massage intensities (= Automatic mode).

2) Description of the "FULL MASS":

Belts:

The belt is fastened using attachment belts located on each side of the motor. Easy to use, you can enjoy a massage while reading, while resting or even while working.

Manual control:

Simply press the manual control button on the "FULL MASS" to define the massage of your choice (according to the various cycle speeds) and to target specific areas: belly/back/buttocks/thighs/calves. An adhesive strip has been placed on the latter so that it can be stored on the side of the belt.

3) Use of the "FULL MASS":

- Connect the adapter to a plug socket and connect the other end to the female input connector located on the side of the belt.
- Turn the appliance on or off by pressing the ON/OFF button (1).


Choice of two modes (Manual or Automatic):

Manual Mode:

- With the FULL MASS, you can select the desired direction of rotation (clockwise or anti-clockwise) by pressing the button (3).
- 6 different intensities can also be selected. To do this, press the corresponding buttons (4 & 5). Each press corresponds to one level of intensity.

Automatic Mode:

- When the appliance is turned on, it automatically reverts to the pre-programmed mode. If you want to change the mode, just press the button (2).

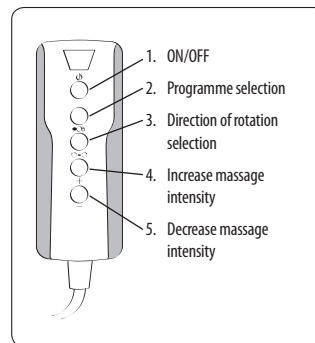MAINTENANCE

Storage:

- Place the massage belt in its box or in a dry place.
- Prevent it from coming into contact with sharp edges and other pointed objects which might cut or pierce its surface.
- In order to avoid any breakage of the power lead, do not wind it around the belt (the same applies to the control cable).
- Do not hang the belt up by the power lead or by the manual control lead.

Cleaning:

- Disconnect the belt from the power socket and wait for a few minutes before cleaning it.
- Clean it with a soft, slightly damp sponge.
- Never let water or any other liquid come into contact with the belt.
- Do not immerse it in liquid to clean it.
- Never use abrasive detergents, brushes or a thinning agent to clean it.

LIMITED WARRANTY

LANAFORM® guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM® guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM® product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM®'s control.

LANAFORM® may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM® will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM®'s Service Centre®. Any attempt to maintain this product by a person other than LANAFORM®'s Service Centre will render this guarantee void.

TECHNICAL SPECIFICATIONS

- Voltage: Input: AC220~240V
Output: DC24V
- Frequency: 50-60Hz
- Power: 40W

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