LANAFORM®

Hot & Cold Mass

Handheld massager

6
MASSAGE ACCESSORIES

Operates
hot and cold

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Thank you for choosing Hot and Cold MASS by LANAFORM®. The technology used in the device combines relief and massage through the application of heat and cold.

Why hot and cold treatment? This concept is particularly indicated to quickly relieve and ease muscle fatigue and pain. The application of heat and cold stimulates the circulation of blood much quicker than a simple massage. This device includes six massage accessories that allow you to choose the type of massage you want.

PLEASE READ ALL INSTRUCTIONS BEFORE USING YOUR "HOT & COLD MASS", PARTICULARLY THOSE THAT CONCERN FUNDAMENTAL SAFETY INSTRUCTIONS:

- Verify that the voltage of your outlet is the one used in the device.
- This device is designed to be used for short periods of time. It can be operated continuously for 30 minutes, but it is essential to switch it off to cool down and recharge the motor before using it again.
- Do not leave the device unattended while it is on or while it is plugged in.
- Without medical advice, the use of the "HOT & COLD MASS" is strongly advised against for women who are pregnant, people with thrombosis or undefined body pains, as well as those with a pacemaker.
- Never use this device if the plug or power cord are damaged, if they do not operate properly, if the device has fallen, has broken, or has been placed in water or another liquid.
- Do not try to recover a device that has fallen in water. Unplug it immediately.
- If the power cord is damaged, it must be replaced by a special cord or similar available from a post-sales service provider.
- Do not wrap the cord around the device.
- Do not carry this device by its power cord.
- Do not insert any object into one of its openings.
- Never obstruct the air ventilation grilles.
- As with any electrical device, it is not advised to use the "HOT & COLD MASS" during thunderstorms.
- Never use this device near humid areas, such as bathtubs, sinks, etc.
- Never place the device in water, in your bath, or in the shower.
- Keep the cable away from all hot or sharp objects.
- Remove the divider before cleaning the device and after recharging it.
- Do not use the "HOT & COLD MASS" near your eyes or any other sensitive, infected, or injured part of the body.
- Do not use this device in a room where aerosol sprays are used or in a room where oxygen is administered.
- Do not use the device under a cover or cushion. This can cause overheating that can cause a fire, electric shock, or injury.
- This device is not intended for use by people, including children, whose physical, sensory, or mental capacities are reduced, or by people without experience of knowledge, unless they can benefit, through the assistance of a person who is responsible for their safety, from monitoring or instructions on the use of the device given beforehand. Children should be monitored to ensure that they do not play with the device.
- Store this device out of the reach of children.
- Use the "HOT & COLD MASS" only for its intended uses and what is described in this manual.

FEATURES

- Hot and cold massage head.
- 180 degree rotational massage head.
- 6 custom massage attachments.
- Light weight. Long handle with hand loop design for easy reach of the body back.
- For easier use, the device is rechargeable and has a life of around 30 minutes.

INSTRUCTIONS FOR USE

Charging the massage unit

- Before the first use, you have to recharge the unit for 24 hours to ensure the built-in rechargeable battery is fully activated.
- After the first use, you only need to recharge for 7 hours.
- For each full recharge, the unit can operate for 30 minutes.
- Turn the unit switch to off position. Plug charger unit to the electric outlet and connect the receptacle into the end of the massage unit.
- A red warning light appears during.
- After full recharge, unplug the charger.
- Recharge the unit after using approximately 20 minutes.

HOW TO USE

- For hot and cold massage, simply select the Hi/Lo massage by pushing the massage switch to high or low.
- Push the switch to cold position.
- Press gently the hot and cold massage head to the desired area of muscle for 5 seconds for cold massage treatment. It will take about 30 seconds for the metal head to turn cold. (approximately 13 °C)
- For the "hot" massage, move the button to the "Heat" position. After about 1 minute, the massage head will be warm. (around 55 °C)
- Then press gently the hot and cold massage head to the same spot area for 30 seconds hot massage.
- We recommend performing 5 care cycles per day. A care cycle consists of 5 seconds of cold massage and 30 seconds of hot massage. Do not perform more than 5 massage cycles on the same area of muscle.
- Select the desired massage attachment heads. There are 6 custom massage heads for your selection.
- It is recommended you not to use the unit for more than 20 minutes for each massage unit.
- Having finished the massage, slide the switch to «Off» to turn off the unit.

FUNCTIONS OF THE ACCESSORIES

A. Accupoint massage attachment

This accessory allows you to perform a deep, vigorous massage of tired muscles in order to relieve them.

B. Double rolling massage attachment

A massage performed with two wheels leads to an increase in blood circulation. During the use, move the unit up and down on the body directly.

C. Gentle massage attachment

This attachment is for the whole body. It relaxes sore and stiff muscles effectively.

D. Shiatsu massage attachment


This attachment is designed to simulate the human fingers so that it provides the deep penetrating massage feeling of human hands to get rid of muscles tiredness effectively.

E. Thump massage attachment

Simulates hand thump. It massages muscles deeply and calms tired muscles.

F. Multi-nodes massage attachment

With micro-finger points, this attachment provides deep and strong penetrating massage to large muscle areas. You are relaxed and your muscles' pain is reduced.

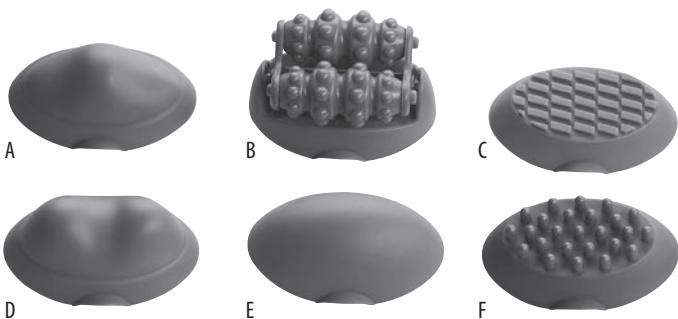MAINTENANCE

To clean: be sure to unplug the unit and allow it to cool down before cleaning. Use a soft, slightly damp cloth to wipe. Never allow water or any other liquids to come into contact with the handle. Attachments can be cleaned in soap water.

Storage: Unplug the appliance from the outlet and allow it to cool down before storing it in its box or in a clean and dry place. Do not hang the unit by the power cord.

LIMITED GUARANTEE

LANAFORM® guarantees that the product is free from any material or manufacturing defect for a period of two years as from the date of purchase, except as stipulated herebelow.

The guarantee on this LANAFORM® product does not cover damage caused by inappropriate or misuse, accident, attachment of any unauthorized accessory, any modification made to the product or any circumstance of any kind beyond the control of LANAFORM®.

LANAFORM® cannot be held liable for indirect, consequential or special damage of any kind.

All implicit guarantees of the aptitude of the product are limited to two years as from the initial date of purchase.

On reception, LANAFORM® will repair or replace your appliance as appropriate and will return it to you. This guarantee must be carried out by LANAFORM® Service Centre. Any maintenance of this product entrusted to anyone other than the LANAFORM® Service Center voids this guarantee.

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