**LANAFORM®**

# Multi Mass

Handheld massager

Removable  
ergonomic handle

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# INTRODUCTION


Thank you for choosing the MULTI MASS from LANAFORM®, a high-quality massage device that guarantees complete satisfaction.

Thanks to its handheld design, detachable handle and ultra-light weight, you will not tire of using it and will enjoy the extremely comfortable massage for the various parts of the body.

The two massage accessories allow you to choose the type of massage that best suits your needs: either the gentle massage fitting or the acupressure fitting for a vigorous massage.

## DESCRIPTION OF THE APPLIANCE

- 1) Stop-start switch / massage / massage + infrared heat
- 2) Massage intensity adjustment button
- 3) Adhesive adjustable strap for the back of the hand
- 4) Infrared heat massage head
- 5) Detachable handle
- 6) Handle release buttons
- 7) Cover for the device
- 8) Gentle massage accessory
- 9) Acupressure massage accessory
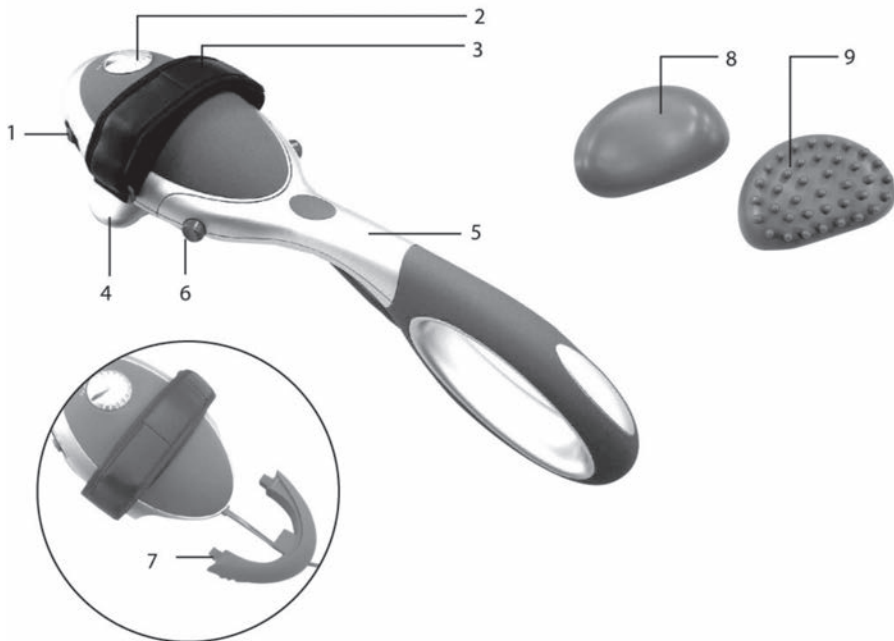
## FEATURES

- Light
- Extra long detachable handle, to easily reach all parts of the body
- Large heating surface to relax muscles
- Adjustable massage intensity
- Adhesive strap for the back of the hand

### **PLEASE READ ALL THE INSTRUCTIONS PRIOR TO USING YOUR «MULTI MASS», IN PARTICULAR THESE FEW BASIC SAFETY INSTRUCTIONS:**

- Only use this appliance for the purposes described in this manual.
- Check that your mains voltage is compatible with that of the appliance.
- Do not use this device for more than 10 minutes in this same position or on the same area of the body. Usage should not exceed 20 minutes in total. A rest period of 30 minutes should be respected between each use.
- Do not leave the device unattended when switched on or plugged in at the socket.
- Unless advised by a doctor, the Multi Mass is not recommended for pregnant women, people prone to thrombosis or suffering from body pain, the cause of which is not yet known as well as people who wear a pacemaker.
- Never use this device on a part of the body that is swollen or inflamed or where there is a rash.
- Store this device out of reach of children.
- If the power lead is damaged, it must be replaced by a special or similar lead available from the supplier or after-sales department.
- Do not use this appliance if the plug is damaged, if it is not working correctly, if it has been dropped or damaged or if it has fallen into water. If this is the case, have the appliance inspected and repaired by the supplier or its after-sales department.
- Always disconnect the appliance from the mains after use and before cleaning it.
- Should the appliance fall into water, unplug it immediately before retrieving it.
- Keep the cable away from any hot or sharp objects.
- Check the temperature of the device before using on any sensitive area.
- Do not use the Multi Mass around the eyes or near any other sensitive area of the body or around any area of the body that is infected or injured.
- Do not use the device while sleeping or without constant attention.
- Do not use before going to sleep. This device has a stimulating effect and may delay the onset of sleep.
- Never use this appliance in a room where aerosol products (sprays) are used or in a room where oxygen is administered.
- This product is a non-professional massage appliance designed to soothe tired muscles. Do not use this appliance as a substitute for proper medical care.
- Do not use the device under a cover or cushion. This may cause overheating that could lead to a fire, electric shock or personal injury.
- Do not carry this appliance by its power supply lead or use this lead as a handle. The device should not be carried by its handle.
- Never obstruct the appliance's air vents and never place the device on an unstable surface such as a carpet or mattress. This could prevent air circulation. Make sure that the air vents are not clogged with hair, lint, dust, etc. Do not insert any objects into the air vents.
- To unplug the device, place all controls in the "Off" position.


- This appliance is not intended for use by persons (including children) whose physical, sensorial or mental capabilities are impaired, or by persons lacking experience or knowledge, unless they are supervised or instructed in the use of the appliance first by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

## INSTRUCTIONS FOR USE

- When you plug in the device, make sure that the device is in “Off” position.
- To switch it on, simply slide the switch into the “M & M+H” position (see device description, point 1).
  1. M = Massage only
  2. M+H = Massage with infrared heat.
- You can use the two accessories for different types of massage. For a strong vigorous massage, use the fixed round massage head.
  1. Gentle massage fitting (see device description, point 8).
  2. Acupressure fitting for a vigorous massage (see device description, point 9).
- It is not recommended to massage the same area for more than 5 minutes. However, the duration of the massage depends on how you feel.
- It is not recommended to exceed a total massage time of 20 minutes, to allow the device to cool down. A rest period of 30 minutes should be respected between each use.
- When the massage is finished, switch the device off by sliding the switch to the «0» position, unplug it, and leave it to cool down then store it away in a dry place.
- To select the massage intensity, slowly turn the intensity adjustment button, from “Lo” to “Hi”, until the desired intensity is obtained (see device description, point 2).

This infrared and impact manual massage device has a detachable handle and handheld design.

For manual operation, press the two handle release buttons (see device description, point (6)) and remove. Then, click the “cover” onto the end of the device, as shown in the photo, and adjust the adhesive strap to the size of your hand (see device description, point (3)). You can then use the device to comfortably massage any part of the body, except the back.

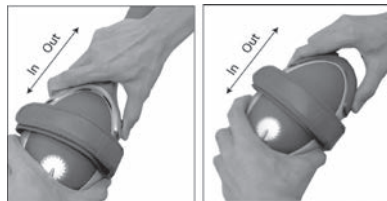To massage the back, you should replace the handle, by carrying out the reverse of the instructions above.

## MAINTENANCE

### CLEANING:

1. Disconnect the Multi Mass from the power socket and wait for a few minutes before cleaning it.
2. Clean it with a soft, slightly damp sponge.
3. Never let water or any other liquid come into contact with the Multi Mass.
4. Do not immerse it in liquid to clean it.
5. Never use abrasive detergents, brushes or a thinning agent to clean it.

**STORAGE:**

Unplug the device at the socket and leave it to cool down before storing it away in its box in a clean dry place, out of reach of children. Do not wrap the power lead around the appliance. They may use and break it. Do not hang the device by its cable.

**ADVICE REGARDING THE DISPOSAL OF WASTE**

- All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.
- When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

**LIMITED WARRANTY**

LANAFORM® guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM® guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM® product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM®'s control.

LANAFORM® may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM® will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM®'s Service Centre®. Any attempt to maintain this product by a person other than LANAFORM®'s Service Centre will render this guarantee void

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Grazie per aver scelto MULTI MASS di LANAFORM©, un apparecchio per massaggi di qualità che vi garantirà la massima soddisfazione.