
MEDISANA®

home of wellness

- DE** Handmassagegerät HM 850
GB Hand Massage Device HM 850
FR Appareil de massage manuel HM 850
IT Massaggiatore manuale HM 850
ES Equipo para masaje manual HM 850
PT Massajador manual HM 850
NL Handmassage-toestel HM 850
FI Käsikäyttöinen hierontalaite HM 850
SE Handmassageapparat HM 850
GR Χειροκίνητη συσκευή μασάζ HM 850

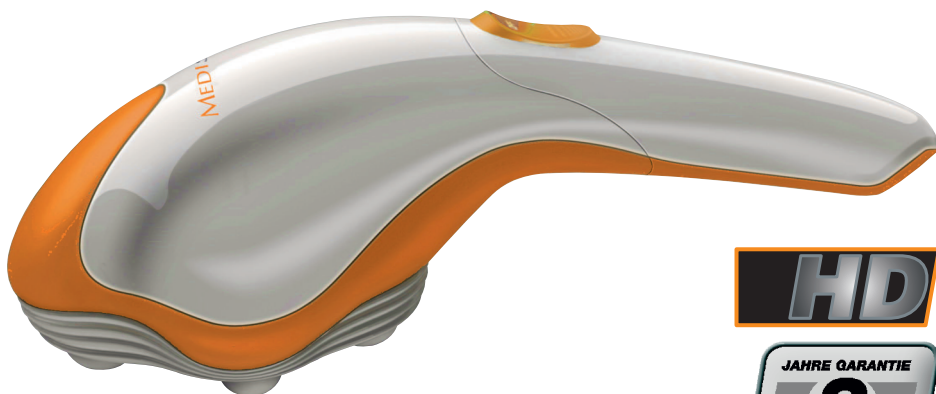Art. 88295

Gebrauchsanweisung
Instruction Manual
Mode d'emploi
Istruzioni per l'uso
Instrucciones de manejo

Bitte sorgfältig lesen!
Please read carefully!
A lire attentivement s.v.p.!
Da leggere con attenzione!
¡Por favor lea con cuidado!

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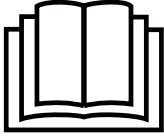
IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols

This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

LOT

Lot number

Manufacturer

Safety Information

- This unit is intended for private use only. If you should have concerns about your health, consult a doctor before using.
- The device must only be used for its intended purpose as described in the instruction manual. Using the unit for any other purpose invalidates the warranty.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Do not use the hand massage device **HM 850**, or only after consulting your doctor, if:
 - you are pregnant,
 - you have a cardiac pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- Special attention is necessary, when the unit is used on, from or in the reach of children, ill or helpless people.
- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- Children must be supervised to ensure that they do not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- Do not use the appliance near your eyes or near any other sensitive areas.
- Very long periods of use can lead to overheating. So after continuously operating for 15 minutes, user should stop it for **at least** 15 minutes before using it again.
- Do not use the appliance if it shows signs of damage, if it is not in perfect working order, or if it has been dropped or dipped into water. In order to prevent hazards, always send the unit in to the service centre.
- Only use suitable and in the scope of delivery originally included accessories together with this item.
- Keep the cable and the appliance away from heat, hot surfaces, moisture and liquids. Do not touch the mains plug when you are standing in water and always take hold of the plug with dry hands.
- Do not allow the appliance to come into contact with pointed or sharp objects.
- Do not cover up the appliance when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.
- If the unit is not going to be used for a long period, please remove the batteries.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using the device.

Only use the massager indoors.

Do not use the massager in damp rooms (e.g. when bathing or showering).

2 Useful Information

Thank you very much

Thank you for your confidence in us and congratulations! With the **Hand Massage Device HM 850** you have acquired a quality product from **MEDISANA**. To ensure that your **Hand Massage Device HM 850** provides successful results for a long time, we recommend that you read carefully through the following notes for the use and care.

2.1 Items supplied and packaging

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your dealer or your service centre.

The following parts are included:

- 1 **MEDISANA** Hand Massage Device **HM 850**
- 2 Batteries 1,5V, C (Baby / UM2)
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required.

If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

2.2 Scope of application

The tapping massage with hand massage device **HM 850** is intended for targeted, deep-efficiently tapping massage treatment of the whole body. The tapping massage is a relaxing massage applied to the nerve tracts to loose blockings and improve the blood circulation. The well-being and the relaxation of muscles are supported. The device has 2 levels of intensity.

CAUTION

Make sure that the maximum operating time of 15 minutes is not exceeded.

3 Operating

3.1

Insert / replace batteries

To insert / replace batteries, press the latching nose on the underside of the unit to the inside and take the lid off. Insert 2 new batteries 1.5V, C (Baby / UM2). Pay attention to the correct polarity (as marked in the battery compartment). Put back the lid into its original position (latching nose need to "click" in place).

Replace the batteries when the massage intensity of the massage heads diminishes.

WARNING BATTERY SAFETY INFORMATION

- Do not disassemble batteries!
- Clean the battery / item contact points before inserting the batteries if necessary!
- Remove discharged batteries from the device immediately!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes! If battery acid comes in contact with any of this parts, rinse the affected area with copious amounts of fresh water and seek medical attention immediately!
- If a battery has been swallowed seek medical attention immediately!
- Replace all of the batteries simultaneously!
- Only replace with batteries of the same type, never use different types of batteries together or used batteries with new ones!
- Insert the batteries correctly, observing the polarity!
- Remove the batteries from the device if it is not going to be used for an extended period!
- Keep batteries out of children's reach!
- Do not attempt to recharge these batteries! There is a danger of explosion!
- Do not short circuit! There is a danger of explosion!
- Do not throw into a fire! There is a danger of explosion!
- Only store unused batteries in the original packing and not in reach of metallic items in order to avoid short circuit!
- Do not throw used batteries into the household refuse; put them in a hazardous waste container or take them to a battery collection point, at the shop where they were purchased!

**3.2
Use**

With the function switch ❶ you can adjust the intensity of the massage (position 1 = low; position 2 = high). Apply the massage heads ❷ lightly on the part of the body to be massaged. Do not concentrate on one part of the body for too long; rather change the area of massage frequently. After treatment, switch the device off by moving the function switch ❶ to position 0.

 **WARNING**

**Only apply gentle pressure during massage!
Too strong pressure may lead to injuries and/or damages
on the device!**

**3.3
Duration of
massage
treatment**

In general, you can allow approximately 10 to 15 minutes operation for each massage session, one to two times a day. Make sure the appliance is cooled down completely before using it again!

4 Miscellaneous

**4.1
Cleaning and
maintenance**

- Remove the batteries before cleaning the unit.
- Clean the unit using a soft cloth lightly moistened with a mild soapy solution. Rub the unit dry with a clean dry cloth.
- Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol.
- Never immerse the appliance in water to clean it and ensure water cannot get into the appliance.
- Wait until the device is completely dry before using it again.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

**4.2
Disposal**

This product should not be disposed of together with domestic waste. Please hand in all electrical or electronic devices at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority for information about disposal. Please remove the batteries before disposing of the equipment. Do not dispose of old batteries with your household waste, but dispose of them at a battery collection station at a recycling site or in a shop. Consult your municipal authority or your dealer for information about disposal.

| | | |
|---------------------------------|----------------------|--|
| 4.3 | Name and model | : MEDISANA Hand Massage Device HM 850 |
| Technical specifications | Voltage supply | : 3 V $\overline{\text{---}}$, 2 x 1,5 V batteries, type C (Baby / UM2) |
| | Massage heads | : 4 |
| | Operating conditions | : Temperature 10°C - 40°C (50°F - 104°F), rel. humidity up to 90% |
| | Storage conditions | : Temperature 2°C - 30°C (35.6°F - 86°F), dry and dust-free |
| | Dimensions | : Approx. L 29 x W 9.6 x H 9.3 cm |
| | Weight | : Aprox. 296 g (without batteries) |
| | Article No. | : 88295 |
| | EAN Code | : 40 15588 88295 1 |

In accordance with our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

Warranty and repair terms Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

MEDISANA AG

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
41468 NEUSS

Germany

eMail: info@medisana.de

Internet: www.medisana.com

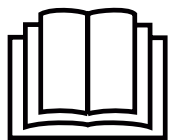The service centre address is shown on the attached leaflet.

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