## DE Büro-Massagestuhl MSO

GB Massage office chair MSO

FR Chaise de bureau massante MSO

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## GB

- 1 4 three-dimensional, movable massage heads
- 2 2 vibration massage motors in the thigh area
- 3 Seat inclination
- 4 Bogie
- 5 Weight adjustment
- 6 Height adjustment
- 7 Button for automatic cable rewind
- 8 On/Off button for massage function
- 9 Auto-button
- 10 Complete back massage button
- 11 Shoulder massage button
- 12 Hip massage button
- 13 Shiatsu position buttons (height adjustment by shiatsu massage)
- 14 Massage adjustment button (width adjustment by roll massage)
- 15 Vibration massage button
- 16 Heat function button
- 17 Roll massage button
- 18 Shiatsu massage button

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## IMPORTANT INFORMATION! DO NOT THROW AWAY!

Read the instruction manual, especially the safety instructions, thoroughly before using the product and keep the instruction manual for later reference. If you give the unit to a third party be sure they also receive the instruction manual.

**Key:**

This instruction manual belongs to this unit. It contains important information for initial operation and use. Read the entire instruction manual. Failure to observe these instructions can lead to severe injury or damage to the unit.

**WARNING!**

Follow warning notices to prevent possible injury to the user.

**CAUTION!**

Follow these instructions to prevent possible damage to the unit.

**HINT**

These helpful hints provide you with additional information about installation or use.

Protection class II

LOT number

Manufacturer

**Safety  
instructions*****Power supply***

- Before connecting the unit to a power source make sure that the information on the rating label matches the mains voltage.
- Do not operate in the vicinity of high-frequency electromagnetic transmitters.
- Keep the power cable and the unit away from heat, hot surfaces, moisture and liquids. Never touch the power plug or On/Off button with wet or damp hands or while standing in water.
- Unit components carrying voltage may not come into contact with liquid.
- Do not reach for the unit if it falls into water. Unplug the power cable immediately.
- Connect the unit so that the power plug is freely accessible.
- Always unplug the power plug from the power outlet immediately after use.
- Disconnect the unit from the power grid by pulling on the plug. Never pull on the cable!
- Do not carry, pull or turn the unit by the power cable.

***for special persons***

- This unit is not intended for use by persons (including children) with restricted physical, sensory or mental abilities, lack of experience and/or lack of knowledge, unless they are supervised by someone responsible for their safety and have received instructions on how to use the unit.
- Children must be supervised to ensure that they do not play with the unit.
- Do not use this unit as an aid to or a replacement for medical applications. This may exacerbate chronic illness and symptoms.
- You should refrain from using the Office Massage Chair MSO or consult your physician prior to use if:
  - you are pregnant,
  - you have a pacemaker, artificial joints or electronic implants,
  - you suffer from one or more of the following illnesses or medical conditions: circulatory disorder, varicose veins, open wounds, contusions, skin cracks, phlebitis.
- The unit has a hot surface. Persons sensitive to heat must be careful when using the unit.
- Discontinue use and consult your physician if you feel pain or if the massage feels uncomfortable.

***before using the unit***

- Check the cable and the unit carefully before every use for damage. Do not put a defective unit into operation.
- Do not use the unit if damage to the unit or the cable is visible, if it does not function properly, if the cable or the unit have become wet. Send the unit to the service point for repairs to avert hazards.

***for using the unit***

- Use the unit only for its intended purpose as according to the instruction manual.
- In case of misuse the warranty is no longer valid.
- Do not leave the unit unattended while it is connected to the power grid.
- Avoid letting the unit come into contact with pointed or sharp objects.
- Never place the unit directly next to an electric oven or other heat source.
- Do not stand on the unit.
- Discontinue use immediately if using the unit feels uncomfortable or painful.
- Do not use the massage function for longer than 15 minutes. A too lengthy massage duration can lead to overstimulation of muscles and to nerve irritations.
- Verify that seat supports and screw connections are in their correct place before every use.

**Do not rotate on the chair when the power cable is plugged in. This can break off the cable.**

**The unit is not intended for household or commercial purposes or for the medical sector. Consult your physician before using the chair if you have health concerns.**

**Use the Office Massage Chair only in closed spaces.**

**Do not use the Office Massage Chair in damp spaces (e.g. while bathing or showering).**

***for care and maintenance***

- The chair is maintenance-free.
- In case of malfunction do not repair the unit yourself. This does not only annul the warranty, but can also cause serious danger (fire, electric shock, injury). Only let authorized service points carry out repairs.

## 2 Useful Information

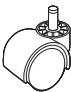
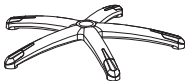
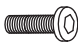
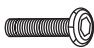
### Thank you for your confidence and congratulations!

In purchasing the **Office Massage Chair MSO** you have acquired a quality product of **MEDISANA**. In order to achieve desired success and long-lasting pleasure with your **MEDISANA** Office Massage Chair we recommend that you carefully read the following tips for use and maintenance.

### 2.1 Shipment and packaging

First verify that the shipment is complete and that the article is not damaged in any way. If in doubt, do not use the article and contact a service point. The shipment includes:

- 1 **MEDISANA** Office Massage Chair **MSO**, consisting of the following components:
  - 1 backrest cushion with automatic cable rewind
  - 1 seat cushion
  - 1 headrest
  - 1 left armrest
  - 1 right armrest
  - 1 adjustable seat
  - 1 gas spring with plastic coating
  - 1 stand (five-armed)
  - 5 wheels
- 1 bag for small parts with 4 caps, 8 screws M8x25mm, 4 screws M8x35mm, and 1 Allen key M8
- 1 instruction manual

	4 caps	  		
	8 screws M8x25mm			
	4 screws M8x35mm			
	1 Allen key M8			
				
1 adjustable seat	1 seat cushion	1 headrest	1 backrest cushion with automatic cable rewind	2 armrests (1 left & 1 right)

Recycle packaging materials or return them to the raw material cycle. Please dispose properly of packaging materials you no longer need. Contact your distributor immediately in the event of damage caused during transit.


## WARNING

Do not allow children to play with packaging films. Danger of suffocation!

### 2.2 Assembly


#### Step 1:

Place the 5 wheels into the ends of the stand.

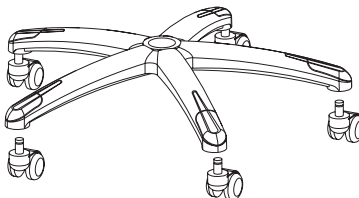
#### Step 2:

Place the gas spring in the designated hole in the stand.

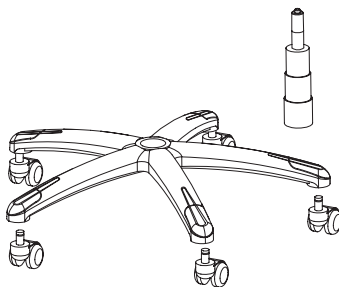
#### Step 3:

Mount the backrest on the seat with 4 screws M8x25mm as illustrated. Initially, tighten the screws by hand.

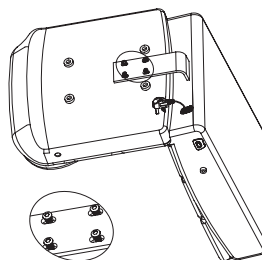
#### Step 4:

Attach the connector to the right armrest and then secure both the right and the left armrests with 2 screws M8x25mm respectively, as illustrated. Tighten the screws. Place the 4 provided caps onto the screw heads.

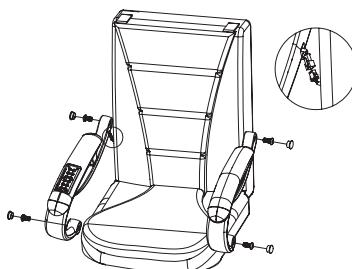**Step 5:**


Now tighten the 4 screws from step 3 with the Allen key. Secure the adjustable seat on the lower surface of the seat cushion with 4 screws M8x35mm and tighten them.

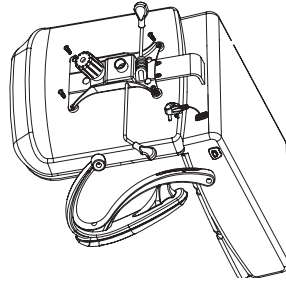**Step 6:**

Now connect the two completely assembled chair parts (upper and lower part). After connecting both chair parts shake them lightly so that all the parts adjust to each other. Now place the headrest on the backrest cushion using the Velcro strip.

**Step 7:**

Connect all cables as shown in the illustration. The Office Massage Chair MSO is now completely assembled. Before use check to make sure that all screw and plug connections between the seat and the stand are tight and correctly placed.
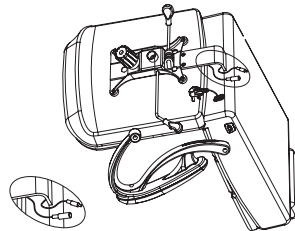
## 3 Operating

### 3.1 Operation

The Office Massage Chair offers ergonomic seating comfort. Application of the Office Massage Chair's massage function heightens your sense of well-being after a strenuous day. The massage effect proves itself to be especially helpful in dealing with cramped muscles and tired tissue. The unit assists in overcoming everyday stress by offering comfortable relaxation. Correctly setting up the chair's adjustable components to fit your body proportions is a prerequisite for its positive use. The Office Massage Chair offers ergonomic seating comfort.

### 3.2 The right seat height

The seat height is ideal when your underarms can be laid on the work space (e.g. desktop) approximately in a right angle to your upper arms. You can adjust the seat height as follows:

**Higher:** Relieve the seat and pull the height setting lever ⑥ upwards.

**Lower:** While seated, pull the height setting lever ⑥ upwards.

### 3.3 Changing the seat inclination

Change the seat inclination to conform optimally to your personal needs. To do so, push the seat inclination lever ③ forwards while seated. By shifting your body weight the seat rocks forward or backward—in this way, you can change the seat angle. Return the seat inclination lever ③ to its original position when the desired angle is reached.

### 3.4 Weight adjustment

Adjust the chair's rocking to your body weight by turning the weight adjustment ⑤. Turning clockwise loosens the chair's rocking mechanics, which is ideal for persons with lower body weight. Turning counter-clockwise tightens the rocking mechanics. Please observe that the chair's max. weight is 120kg!

### 3.5 Massage functions

The Office Massage Chair MSO has 4 functions for targeted massage application: Finger pressure massage (shiatsu massage), vibration massage and roll massage. Additionally, when using the finger pressure (shiatsu) or roll massage, you can activate a heat function or you can also active the vibration massage. Never use the unit for longer than 15 minutes. The unit shuts off automatically after ca. 15 minutes. Always let the unit cool completely down before using it again. Operate the various functions using the buttons on the control unit on the right armrest:

- Pull the power cable out of the automatic cable rewind in the backrest. Connect the plug to a power outlet. Turn on the unit by pushing the On/Off button ⑧. The LED starts to glow. If you wish to turn off the unit again, simply push this button again. While the massage heads are being retracted the LED blinks and goes out when the unit is turned completely off. Unplug the power cable from the outlet and press the automatic cable rewind button ⑦. The cable winds back.
- Pushing the Auto-button ⑨ begins a 15-minute application of all available functions (shiatsu, roll and vibration massage incl. activated heat function). Afterwards the unit turns off automatically.

## WARNING

**Do not rotate on the chair when the power cable is plugged in. This can break off the cable.**

- Pushing the “Complete back massage” button **10** causes the LED to light up and the massage heads begin a massage that expands over the entire back. Please choose the desired massage type as well (shiatsu or roll massage).
- Pushing the “Shoulder massage” button **11** causes the LED to light up and the massage heads begin a massage concentrating on the upper back area. Please choose the desired massage type as well (shiatsu or roll massage).
- Pushing the “Hip massage” button **12** causes the LED to light up and the massage heads begin a massage concentrating on the lower back area. Please choose the desired massage type as well (shiatsu or roll massage).
- Use the Shiatsu position buttons **13** during a shiatsu massage to concentrate the massage on a specific height. Push and hold the desired button to move the massage heads during a shiatsu massage up or down. After releasing the button the heads stay in this position and massage the according spot permanently. Push the On/Off button **8** to end the massage. After 15 minutes the unit shuts off automatically.
- The massage adjustment button **14** (only usable during a roll massage) changes the range of the massage position. Push and hold the button to massage a broader or smaller area. When the desired range is reached, release the button to set this massage setting.
- The Vibration massage button **15** turns on the vibration massage motor in the seat. This massages the legs in the thigh area by means of vibration. The massage begins gently. Pushing the button **15** a second time increases the vibration (medium massage intensity). Pushing the button **15** a third time turns on the highest vibration level (highest massage intensity). Pushing the button **15** a fourth time turns off the vibration massage and the LED goes out.
- The heat function can be added to a shiatsu or roll massage, or you can add a vibration massage. Pushing the button **16** turns on the heat function, the LED starts to glow. The heat radiates off of the massage heads **1**. Pushing the button again turns off the heat function, the red LED goes off.
- Pushing the Roll massage button begins **17** the roll massage function. The roll massage is characterized by more or less heavy pressure to the massaged tissue, which is comparable to a massage by two thumbs. Thereby the massage heads **1** assume alternating different positions. This facilitates a holistic massage. The roll massage is especially usable for stabilizing connecting tissue, promoting blood circulation, detoxifying and stimulating metabolism. Pushing the Roll massage button **17** again deactivates the roll massage.
- Pushing the “Shiatsu massage” button **18** starts the shiatsu massage. Shiatsu, a form of finger pressure massage, belongs to the important massage techniques like manual lymph drainage and reflexology massage. This form of massage was developed in Japan and originates from traditional Chinese massage. Soft touch and soothing pressure awaken the vital energy and cause it to flow. The goal of the massage is not only to release tension but also to harmonize body and mind through optimal flow of the vital energy. Pushing the Shiatsu massage button **18** again deactivates the shiatsu massage.

## 4 Miscellaneous

### 4.1 Malfunction and Repair

Should the unit fail to react after being turned on or after you have chosen a massage function, please check that the electrical power connection was made correctly and that all cables were properly connected. Unplug the unit, wait a few seconds and reconnect the unit to the mains electricity supply. Push the On/Off button  and choose a massage function. Contact a service point if no further reaction is visible.

### 4.2 Cleaning and Maintenance

- Check to make sure that the unit is turned off and unplugged before cleaning. Only clean the unit with a damp and soft rag.
- Never use aggressive cleaning agents, hard brushes, solvents or alcohol. That could damage the surface.
- Do not submerge the unit in water to clean it and take care that no water enters the unit.
- Wait until the unit is completely dry before further use.
- Untangle the cable in case it is twisted.

### 4.3 Information for disposal

Do not dispose of this unit with household garbage. Every consumer is obligated to return all electric or electronic appliances, whether or not they contain pollutants, to a collection point in his or her city or to a local retailer to ensure an environmentally friendly disposal. Consult your local authority or retailer for more information regarding disposal.

### 4.4 Technical data

Name and model	: <b>MEDISANA Office Massage Chair MSO</b>
Power supply	: 100 - 240V~ / 50 - 60Hz
Performance	: 60Watt
Automatic shutoff	: after ca. 15 minutes
Dimensions ca.	: 61 x 76 x 113cm
Weight ca.	: 23kg
Cable length ca.	: 2.50m
Operating conditions	: Only in dry spaces, use as per instruction manual
Storage conditions	: Clean and dry, storage temperature 0°C - max. 40°C
Item No.	: 88920
EAN No.	: 40 15588 88920 2

**In the course of constant product improvement we reserve the right to make technical and design modifications.**

## Warranty and repair terms and conditions

In case of a warranty claim, contact your specialty store or a service point.

Should you have to send the unit in, note the malfunction and include a copy of the receipt.

### The following warranty conditions apply:

1. A warranty of 3 years starting on the date of sale is given on MEDISANA products. In case of a warranty claim, the date of sale must be verified by a proof of purchase.
2. Flaws due to material or manufacturing defects will be resolved without cost during the warranty period.
3. The warranty period is not extended in case of a warranty claim, neither for the unit nor for replaced parts.
4. The warranty does not cover:
  - a. any damages caused by incorrect use, e.g. by failure to observe the instruction manual;
  - b. damages caused by maintenance or repair carried out by the consumer or unauthorized third parties;
  - c. damages caused during transit from the manufacturer to the consumer or during transmittal to the service point;
  - d. accessories subject to normal wear and tear.
5. There is no liability for indirect or direct consequential damages caused by the unit if the damage to the unit is accepted as a warranty claim.

### MEDISANA AG

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41468 NEUSS

Germany

eMail: [info@medisana.de](mailto:info@medisana.de)

Internet: [www.medisana.de](http://www.medisana.de)

The service address is found on a separately enclosed sheet.

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