# **SALTER**

# **MiBABY**

**MiBaby scale with PC tracking software**

---

## **Instructions and Guarantee**

## MiBaby

The MiBaby scales, height chart and PC software are intended to provide a precious record of a child's growth and development as well as provide convenience and reassurance to parents wishing to check and monitor their baby's growth and development in the comfort of their own home. MiBaby should be used in addition with and not as a substitute to professional medical guidance and postnatal care. If you are concerned about any aspect of your child's health please seek professional medical guidance.

### SOFTWARE INSTALLATION

Please go to [www.saltermibaby.com](http://www.saltermibaby.com) and download the latest version of the MiBaby software. The full operating instructions will be included with this download. The MiBaby computer software is compatible with Windows 2000, XP, Vista, Windows 7 and Windows 8. It is not compatible with apple computers. To install the software follow this sequence:


1. Close all open programs.
2. Double-click on the application file that you have just downloaded.
3. Follow the on screen installation instructions.


4. Click Next.


5. Select a folder or click next.

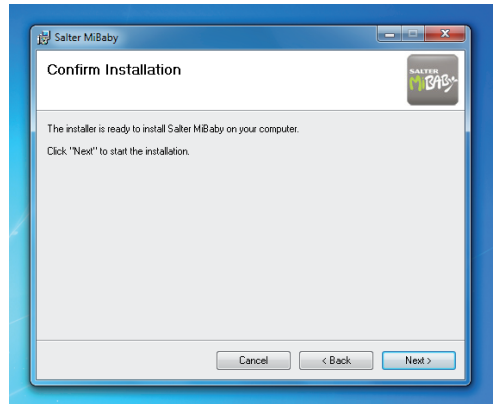6. Click Next.

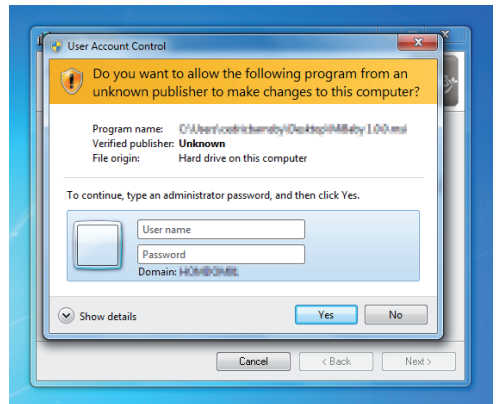7. If requested enter your computer administrator User name and Password.

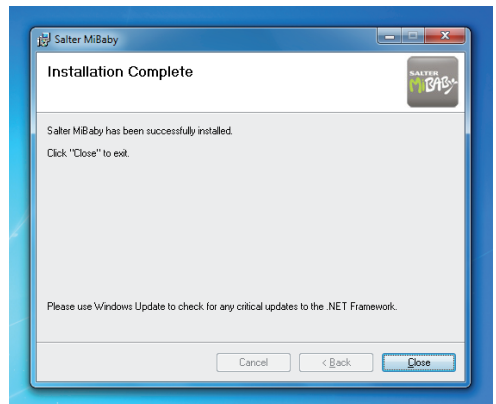8. Click Close.

9. The Salter MiBaby icon should appear on the desktop. The software installation is now complete.

## MiBaby SOFTWARE

Double click on the MiBaby icon to open the software.


Click on a user icon for an existing user, new user to create a new user or import to import existing user data.

### New User

Enter the baby's details and select your preferred units. Click ok to confirm.

### Summary

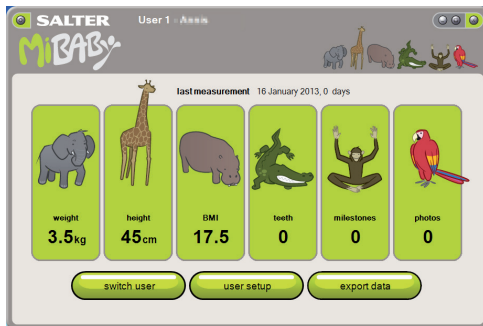The summary page shows the baby's last measurements. Click on any of the animal icons to enter the measurement data.


Click on switch user to select a different child's profile. The MiBaby software supports up to 4 users.

Click user setup to return to the user setup page.

Click export data to save the child's data to another location on the PC. This can then be transferred to another PC and shared with friends and family. Exported data can be imported using the import button on the user selection page.

**Note:** Only the user number selected data is exported. For multiple users data must be exported individually.

### Weight Height and BMI

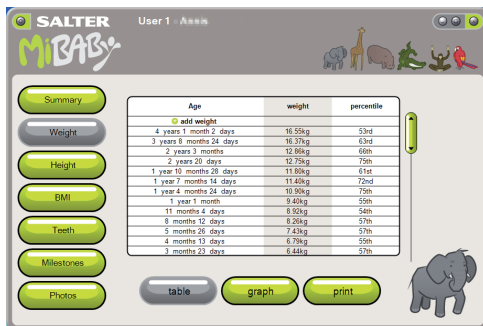You can navigate to Weight, Height, BMI, Teeth, Milestones, Photos or back to the Summary page using the buttons on the left of the page.


Weight Height and BMI measurements are displayed in a table. You can manually add a new weight or height by clicking on the add weight/add height button at the top of the table.

Enter the new weight/height and date and click ok.

Double click on an existing weight/height to edit it.


**Note:** BMI measurements are calculated from the weight and height measurements so there is no need to enter these. Regularly updating both weight and height measurements will give the most accurate BMI measurements.

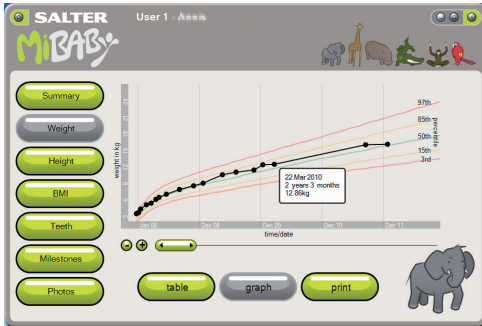
Click on graph to display the measurements as a graph. The results are shown against percentile curves. You can hover over a data point to see details. You can zoom in or out and travel along the timeline using the +, - buttons and slide bar.

Percentile curves are derived from: World Health Organisation (2006). WHO child growth standards: length/height-for-age, weight-for-age, weight-for-length, weight-for-height and body mass index-for-age: methods and development.


### Teeth

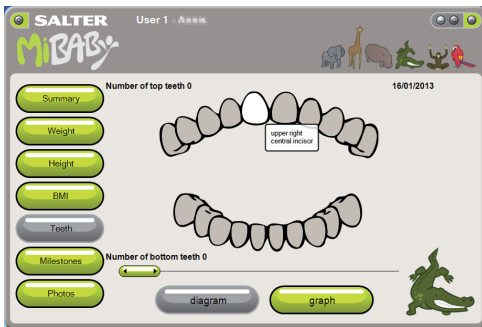Teeth are shown as a diagram.

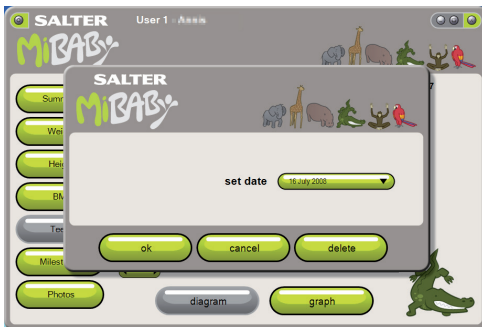
Click on a tooth to enter/edit the date it appeared. Use the slide bar to see the order the teeth appeared.

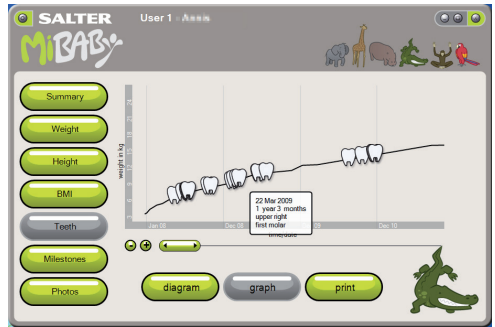
Click on graph to see the order the teeth appeared on a timeline. You can hover over a data point to see details. You can zoom in or out and travel along the timeline using the +, - buttons and slide bar.


### Milestones

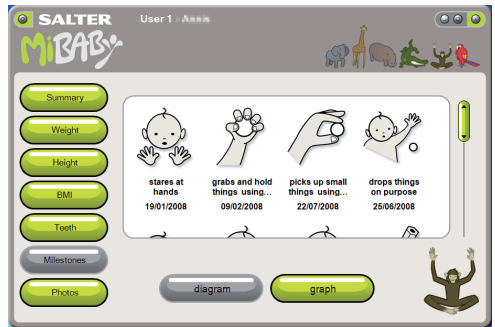Milestones are shown as diagrams.

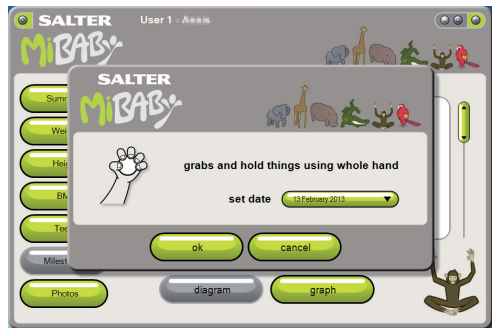Click on the + button to add dates the milestones were achieved. There is also a custom milestone icon with free text. You may add multiple custom milestones to record individual key events in your child's development. Click on an existing milestone to edit.

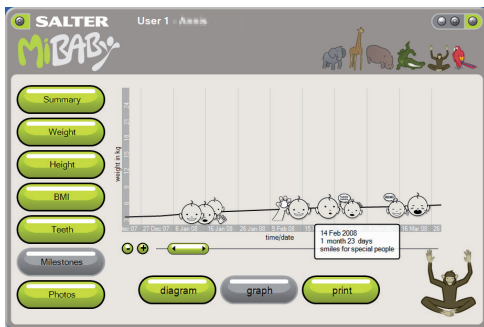
Click on graph to display the measurements as a graph. You can hover over a data point to see details. You can zoom in or out and travel along the timeline using the +, - buttons and slide bar.

### Photos

Photos are shown as a gallery. Click on a photo to see a large version.

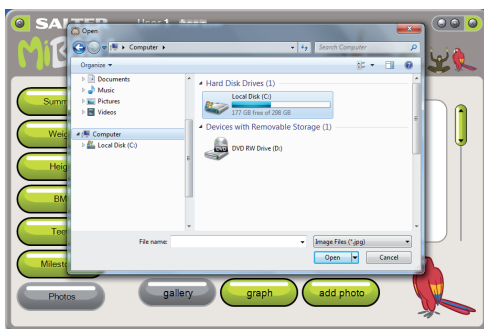Click on add photo and locate the photo file location on your computer to add a photo. Click on add date and add comment to add/edit details.

Click on graph to see the photos on a timeline. You can hover over a data point to see details. You can zoom in or out and travel along the timeline using the +, - buttons and slide bar. Click on the camera icon to see a large version of the photo.

### MiBaby HEIGHT CHART

The printed MiBaby height chart provides your child with a visual record of their height engaging the child in understanding how they grow. The stickers are intended to encourage regular height checks but if they run out the height chart can be written on directly. To fit the height chart first select a suitable location on an unobstructed section of wall. Measure 60cm or 2ft from the floor and mark this position on the wall. Ensure the 60cm or 2ft marker of the height chart aligns with the marked position. Affix the height chart to the wall with a suitable adhesive such as blue tack or double sided tape. Take care not to damage the decoration. When measuring your child ensure they stand up straight, feet together and still and mark their height on the chart. The child may enjoy sticking the sticker on themselves. This height can easily be entered into the MiBaby PC software to provide a digital record of the child's height.

### NEW FEATURE!

This scale features our convenient step-on operation. Once initialised the scale can be operated by simply stepping straight on the platform – no more waiting!

### PREPARING YOUR SCALE

1. Open the battery compartment on the scale underside.
2. Remove isolating tab from beneath the battery (if fitted) or insert batteries observing the polarity signs (+ and -) inside the battery compartment.
3. Close the battery compartment.
4. Select kg, st or lb weight mode by the switch on the underside of the scale.
5. For use on carpet attach enclosed carpet feet.
6. Position scale on a firm flat surface.

### INITIALISING YOUR SCALE

1. Press the platform centre and remove your foot.
  2. '0.0' will be displayed.
  3. The scale will switch off and is now ready for use.
- This initialisation process must be repeated if the scale is moved.**  
**At all other times step straight on the scale.**

### OPERATING YOUR SCALE – ADULT ONLY

1. Select adult mode by the switch inside the battery compartment.
2. Step on and stand very still while the scale computes your weight.
3. Your weight is displayed.
4. Step off. Your weight will be displayed for a few seconds.
5. The scale will switch off.

### OPERATING YOUR SCALE – BABY MODE

1. Select baby mode by the switch inside the battery compartment.
2. Step on and stand very still while the scale computes your weight.
3. Your weight and the icon will be displayed. Then '0.0' will be displayed.
4. Step off the scale and carefully pick up the baby.
5. Step back on the scale holding the baby.
6. The baby weight only and icon will be displayed.
7. Step off. The baby weight will be displayed for a few seconds.
8. The scale will switch off.

## WARNING INDICATORS

- O-Ld** Weight exceeds maximum capacity.
- Lo** Replace battery.
- Err** Unstable weight. Stand still.

## ADVICE FOR USE AND CARE

- Always weigh yourself on the same scale placed on the same floor surface. Do not compare weight readings from one scale to another as some differences will exist due to manufacturing tolerances.
- Placing your scale on a hard, even floor will ensure the greatest accuracy and repeatability.
- Weigh yourself at the same time each day, before meals and without footwear. First thing in the morning is a good time.
- Your scale rounds up or down to the nearest increment. If you weigh yourself twice and get two different readings, your weight lies between the two.
- Clean your scale with a damp cloth. Do not use chemical cleaning agents.
- Do not allow your scale to become saturated with water as this can damage the electronics.
- Treat your scale with care - it is a precision instrument. Do not drop it or jump on it.
- Caution: the platform may become slippery when wet.

## TROUBLESHOOTING GUIDE

If you experience any difficulty in the use of your scale:

- Check the battery is correctly fitted.
- Check that you have selected your preference of stones/pounds, kilograms or pounds.
- Check that the scale is on a flat, level floor and not touching against a wall.
- Repeat the 'Initialisation' procedure every time the scale is moved.
- If when you use the scale there is no display at all or if **Lo** is displayed - try a new battery.
- If **O-Ld** is displayed the scale has been overloaded.

## WEEE EXPLANATION

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

## BATTERY DIRECTIVE

This symbol indicates that batteries must not be disposed of in the domestic waste as they contain substances which can be damaging to the environment and health. Please dispose of batteries in designated collection points.

## GUARANTEE

This product is intended for domestic use only. Salter will repair or replace the product, or any part of this product, free of charge if within 15 years of the date of purchase, it can be shown to have failed through defective workmanship or materials. This guarantee covers working parts that affect the function of the scale. It does not cover cosmetic deterioration caused by fair wear and tear or damage caused by accident or misuse. Opening or taking apart the scale or its components will void the guarantee. Claims under guarantee must be supported by proof of purchase and be returned carriage paid to Salter (or local Salter appointed agent if outside the UK). Care should be taken in packing the scale so that it is not damaged while in transit. This undertaking is in addition to a consumer's statutory rights and does not affect those rights in any way. For UK Sales and Service contact HoMedics Group Ltd, PO Box 460, Tonbridge, Kent, TN9 9EW, UK. Helpline Tel No: (01732) 360783. For Ireland, please contact Petra Brand Masters, Unit J4 Maynooth Business Campus, Maynooth, Co. Kildare, Ireland. Tel +00 353 (0) 1 6510660. e-mail sales@petrabrandmasters.ie.

[www.salterhousewares.com/servicecentres](http://www.salterhousewares.com/servicecentres)

**Register your product today at:  
[www.homedicsgroup.com/register](http://www.homedicsgroup.com/register)**

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