# SANITAS

## PERSONENWAAGE SGS 07, SGS 25

- |                        |                     |                    |                     |
|------------------------|---------------------|--------------------|---------------------|
| GB Bathroom scale      | P Balança de banho  | N Personvekt       | SK Osobná váha      |
| F Pèse-personne        | TR Terazi           | FIN Henkilövaaka   | SLO Osebna tehtnica |
| E Báscula de baño      | GR Ζυγαριά ανθρώπων | RUS Напольные весы | H Személymérleg     |
| I Bilancia pesapersona | DK Personvægt       | PL Waga osobowa    | ميزان أشخاص         |
| NL Personenweegschaal  | S Personvåg         | CZ Osobní váha     |                     |

### Garantie

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
### Important notes

- On delivery the scale is set to the unit “kg”. A switch located on the back of the scale enables you to change to “pound” and “stones” (lb, St).
- Place the scale on a firm level floor; a firm floor covering is required for correct measurement.
- Audibly, briefly and firmly tap the scale with your foot to activate the vibration sensor!
- Stand on the scales without moving and with equal weight distribution on both feet!
- All our personal scales comply with EU guideline 2004/108 + relevant supplements.
- Should you have any questions regarding operation and usage of our appliances, please contact your local retailers or Customer Service.
- Spent batteries and rechargeable batteries do not constitute normal household waste!
- They are considered to be toxic waste and, as such, should be disposed of in special containers, toxic waste collection points or brought to electrical goods dealers.  
You are legally obliged to dispose of the batteries.

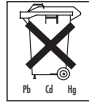- Note: Batteries containing pollutant substances are marked as follows: Pb = Battery contains lead, Cd = Battery contains cadmium, Hg = Battery contains mercury.
- Please dispose of the device in accordance with the directive 2002/96/EG – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please refer to the local authorities responsible for waste disposal.
- Repairs may only be performed by Customer Service or by accredited retailers.

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