SANITAS

SZA 50

- Ⓓ **Gebrauchsanleitung**
Elektrische Zahnbürste
- ⒼⒷ **Instruction for Use**
Electric toothbrush
- Ⓕ **Mode d'emploi**
Brosse à dents électrique
- Ⓔ **Instrucciones para el uso**
Cepillo de dientes eléctrico
- Ⓐ **Инструкция по применению**
Электрическая зубная щетка

Service-Adresse:
Hans Dinslage GmbH
Riedlinger Straße 28
88524 Uttenweiler
GERMANY
Tel.-Nr.: +49 (0)7374-915766
Fax-Nr.: +49 (0)7374-920723
E-Mail: kd@dinslage.net

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Items included in the package:

- Easy grip electrical toothbrush with built-in rechargeable battery
- 4 brush heads
- Induction charging station with mains cord and mains plug, as well as a mounting set for optional wall mounting (2 screws and 2 dowels)
- These instructions for use

Explanation of symbols

The following symbols are used in the operating instructions and on the type plate.

WARNING Warning of risks of injury or health hazards

CAUTION Safety information about possible damage to device/accessories

Note Important information

The device has double protective insulation

1. Getting to know your device

The electric toothbrush has been developed for daily dental cleaning. The bristles are rounded and are of 2 lengths. They help keep the dental enamel and the gums in good condition by allowing deep cleaning in inaccessible places. The result is a long-lasting fresh feeling and clean teeth.

The toothbrush also offers you the following advantages:

- Wall and flat surface mount
- 2 speeds and massage effect
- Timer function with vibration signal
- Charge indicator lamp

2. Proper use

WARNING

- Children under 3 must not use the electric toothbrush.
- This device is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety or they are given instructions on how to use the device. Supervise children to ensure that they do not use the device as a toy.

Only use the device to clean teeth and for gum massage. The device is only for household/private use, not for commercial use. In the case of improper use or incorrect handling, no liability can be assumed for damages resulting.

3. Safety information

Please read these operating instructions completely and carefully. Failure to observe the information provided below can lead to personal injury or material damage. Keep these instructions in a safe place where they can also be read by other persons using the device. Always pass on these instructions together with the device if you sell it or pass it on.

WARNING Electric Shock

The charging station is live as long as it is connected to the mains supply. Electrocutation represents a life-threatening hazard! Therefore observe without fail:

- The charging station must be located at a safe distance from bathtubs, showers, swimming pools, wash basins and water-filled vessels. According to VDE standard 0100 [German Electrical Engineers Association], the safety distance must be at least 60 cm. We recommend a distance of at least 1 m. If the charging station should, nevertheless, fall into water, immediately disconnect the mains plug and have the device inspected by a specialist before using again. This also applies if the mains cord or the charging station are damaged.
- Never touch the powered charging station with moist or wet hands. Never touch the charging station while connected if the charging station itself is moist or wet or if it is standing on a wet floor. Never submerge the device in water and protect it against all forms of moisture.
- The mains voltage must correspond with the voltage specified on the type plate of the device. The device must only be operated from a power outlet with fault-current protection.
- Regularly check the charging station and mains cord, as well as the easy grip electrical toothbrush, for damage. A damaged device must not be used. Never open the device yourself. In this case, please contact our customer service or an authorized dealer.
- The device and the lead must not come into contact with hot surfaces or sharp objects.
- Never pull the plug out of the power outlet with the cord.
- Disconnect the mains plug if the device is not working properly, for cleaning and if the charging station is not in use.
- Do not kink the mains cord or wind it around the charging station.
- Do not operate the device outdoors.
- Never attempt to repair the device yourself. Electrocutation represents a life-threatening hazard!

WARNING Health Hazards

- Slight bleeding of the gums may occur when used for the first time. Should bleeding of the gums persist for longer than two weeks, consult your dentist.
- In the case of surgical procedures on the teeth or gums performed within approx. 2 months, ask your dentist before using this toothbrush.
- Replace the brush head if the bristles are snapped off. Also replace the brush head for hygienic reasons if you drop the device.
- Consult your doctor before using this device if you are fitted with a cardiac pacemaker or other electrical implant.

4. Device description

1 **ON/OFF** switch, speed / massage effect setting (non-slip pressure switch)

2 **OFF** button

3 Charge indicator

4 Grip (waterproof)


5 Charging station

6 Brush head (replaceable)

7 Lid for brush holder

8 Brush holder with removable brush insert

9 Wall mount

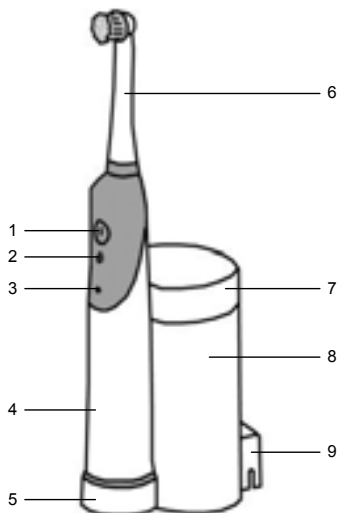5. Start-up

Unpacking

WARNING

- Ensure that the packaging foils do not fall into children's hands. Risk of suffocation!
- Unwind the mains cord on the device and pull it to its full length.

Remove all packaging and keep it for subsequent storage. If you notice transportation damage while unpacking, contact your dealer.

Wall mounting


You can wall mount the charging station (6) with the aid of the mounting set.


WARNING

When mounting, ensure that you do not drill into electrical lines, cables or water pipes.

1. Choose a suitable place in the bathroom.
2. Mark 2 vertically aligned drilling holes each with a diameter of 6 millimeters. The distance between the fixture holes is given by the wall mount.
3. Now drill the holes in the wall and insert the dowels (A2).
4. Remove the wall mount from the charging station by sliding it downward (**Figure 1**).
5. Firmly secure the wall mount to the wall with the two screws (A1) (**Figure 2**).

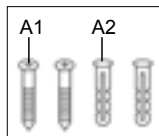6. Now slide the charging station onto the secure wall mount from above (**Figure 3**).

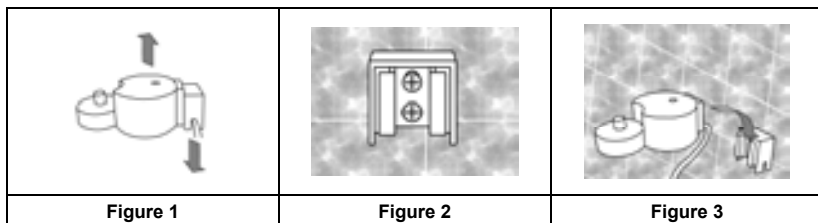Charging the rechargeable toothbrush

Before using for the first time the built-in rechargeable battery in the easy grip electrical toothbrush must be charged. Charging takes at least 24 hours. The easy grip electrical toothbrush is then available for around 60 minutes of operation.

1. Connect the charging station (5) to the mains voltage specified on the type plate.
2. Place the easy grip electrical toothbrush (4) on the charging station (5) such that the control buttons face forward. The charge indicator (3) lights red as soon as the easy grip electrical toothbrush is charging. If the charge is very weak at the start of charging, the charge indicator blinks up to several hours long and then changes to light continuous red. When fully charged, the charge indicator lights green. The inbuilt overcharging protection means that the battery cannot be overcharged.
3. The easy grip electrical toothbrush always remains fully operational if you leave it in the charging station for daily use. The inbuilt overcharging protection prevents damage to the battery.

You obtain the maximum battery capacity if you discharge the battery at least once every 6 months. For this purpose, disconnect the charging station (5) from the mains and discharge the battery through normal use. Then fully recharge the battery.

6. Using the toothbrush

Brushing teeth

WARNING

Slight bleeding of the gums may occur when used for the first time. Should bleeding of the gums persist for longer than two weeks, consult your dentist.

Note

- Dentists recommend cleaning for 2 minutes.
- It is recommended to replace the brush head every 2-3 months.
- For safety reasons and to protect your teeth and gums, the device automatically switches off after 10 minutes.

CAUTION

- Always switch off the easy grip electrical toothbrush when replacing the brush head.

1. Clean the brush head with warm water before using for the first time. Push the brush head onto the easy grip electrical toothbrush. Ensure that the bristles point toward the control buttons and that the brush head is firmly attached.

2. Rinse your mouth out with water before cleaning your teeth.
3. Moisten the toothbrush and apply toothpaste.
4. To avoid splashing, hold the toothbrush at an angle to the tooth surface. Only then switch on the easy grip electrical toothbrush by pressing the **ON/OFF** switch (1). The easy grip electrical toothbrush automatically switches to the highest speed level (High Speed).
5. By repeatedly pressing the **ON/OFF** switch (1) you can select between: high speed, low speed or massage effect.
6. Slowly go from tooth to tooth starting with the outside surface. Then clean the back surfaces of the back teeth, then the inside surfaces and finally the occlusal surfaces. To protect your gums and teeth avoid brushing with excessive pressure.
7. An acoustic vibration signal sounds after 2 minutes of cleaning. Press the **OFF** button (2) to switch off the easy grip electrical toothbrush.
8. Remove the brush from your mouth. Rinse your mouth out well with water.

Replacement brush heads are available in retail or from our customer service. Use the instrument only with the accessories provided.

7. Cleaning and care

Cleaning brush heads and the easy grip electrical toothbrush


WARNING

- Never use aggressive solvents or cleaning agents.
- Never immerse the device in water. Do not put the device in the dishwasher.

Clean the brush head and easy grip electrical toothbrush each time after use.

1. Rinse the brush head under running water. Here the brush is attached to the toothbrush and the toothbrush is switched on (**Figure 4**). Be careful about any splashing water.
2. Remove the brush head from the toothbrush and wash the brush well under running water (**Figures 5 and 6**).
3. Separately rinse the tip of the toothbrush under running water while it is switched on (**Figure 7**).
4. Dry both parts well.

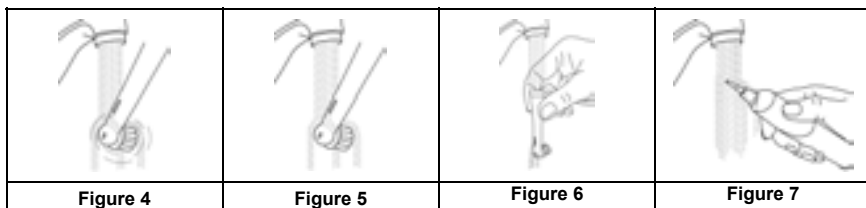Cleaning the charging station and brush holder

WARNING

The charging station with cord and plug must not be submerged in water and must be protected from all forms of moisture.

Clean the charging station and brush holder from time to time.

1. Disconnect the charging station from the power supply.

2. Remove the brush holder and clean it and the removable brush insert with a damp cloth. You can use a mild detergent.
3. Clean the charging station with a slightly damp cloth. Leave the charging station to dry completely before using the next time.

Disposal

This device contains a rechargeable battery fitted into the handle. The battery contains nickel-metal hydride (NiMH). In the interests of environmental protection the easy grip electrical toothbrush including its battery must not be disposed of with domestic waste. Disposal can take place via the relevant collection points in your country.

Follow local regulations regarding the disposal of materials. The appliance should be disposed of according to EC Directive **2002/96/EC-WEEE** (Waste Electrical and Electronic Equipment). In case of queries, contact the municipal authorities responsible for waste disposal in your area.

8. Technical data

Mains connection	100 V-240 V ~, 50/60 Hz
Power consumption	2 W
Rechargeable battery	2 x 1.2 V AA Ni-MH, 900 mAh

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