


SOEHNLE

BODY BALANCE EASY CONTROL

BEDIENUNGSANLEITUNG
OPERATING INSTRUCTIONS
MODE D'EMPLOI
ISTRUZIONI PER L'USO
GEBRUIKSAANWIJZING
INSTRUCCIONES DE MANEJO
MANUAL DE INSTRUÇÕES
BRUKSANVISNING
BRUGSANVISNING
KÄYTTÖOHJEET
KEZELÉSI ÚTMUTATÓ
INSTRUKCJA OBSŁUGI
NÁVOD K POUŽITÍ
ИНСТРУКЦИИ ПО ЭКСПЛУАТАЦИИ
KULLANMA KILAVUZU
Οδηγίες λειτουργίας
NAVODILA ZA UPORABO
UPUTA ZA UPORABU

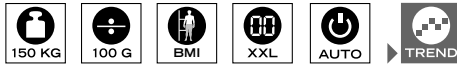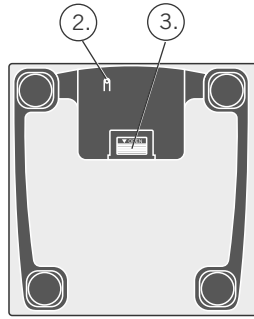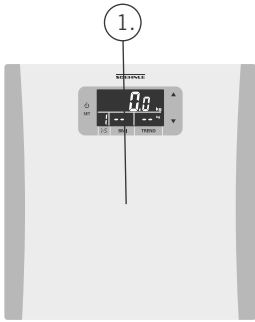DE

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






EN

Entering personal data and weighing

1. Push ON/SET key.
Use ▲/▼ keys to select memory slot and confirm (SET).
2. Then use ▲/▼ keys to set body height and confirm (SET).
3. Finally enter the target

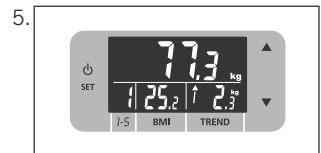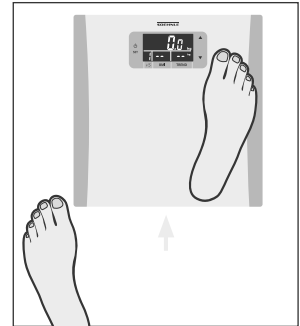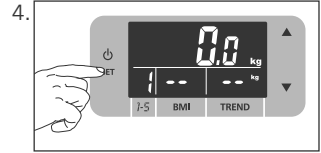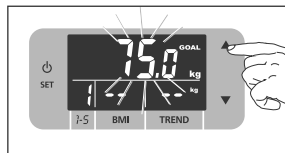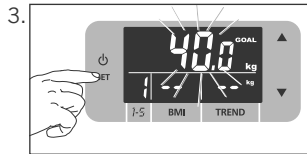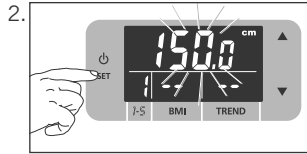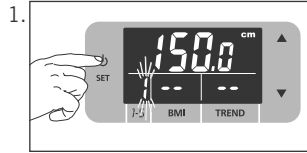- weight using the ▲/▼ keys and confirm (SET).
- The 0-display will appear. Step on the scale for initial weighing. Please stand still for as long as the display flashes.
- After determination of weight, the scale will display the body weight, the memory slot, the calculated body mass index and the

difference of your current to your target weight.

To change or delete data simply write over the existing data.

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EN

Weighing with BMI calculation and display of weight difference

1. Step on the scale.
2. The body weight will be determined and the automa-

tic recognition of persons will be performed. Please stand still for as long as the display flashes.

3. The scale will then display the body weight, the memory slot, the calculated body mass index and the difference of your current to your target weight.
4. If the body weight of other persons that are stored on the scale is approximately the same, a manual allocation will be required. This will be necessary if the memory slot number flashes. In that case simply select the memory slot using ▲/▼ keys and confirm (SET).

The scale will then display the body weight, the memory slot, the calculated body mass index and the difference of the current to the target weight of the correct person.

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EN**Guarantee**

For 3 years as of the date of purchase, Soehnle guarantees the removal of defects due to material or processing errors, free of charge, by way of repair or replacement. In case of a guarantee claim, please return your scale to your seller, together with purchase receipt and guarantee card.

EC Conformity

 This device complies with the applicable EC Directive 2004/108/EC. You will find

the EU compliance declaration on www.soehnle.com.

Battery disposal EC Directive 2008/12/EC

Batteries are not a part of your regular household waste. You must return batteries to your municipality's public collection or wherever batteries of the respective type are being sold.

Batteries containing hazardous materials are marked with these symbols:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

Disposal of electric and electronic devices EC Directive 2002/96/EC

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

Consumer service

If you have questions or comments, the following points of contact are available:

Germany Ph: (08 00) 5 34 34 34
International Ph.: +49 26 04 97 70
Monday to Friday 09:00 to 12:00 am

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Leifheit AG