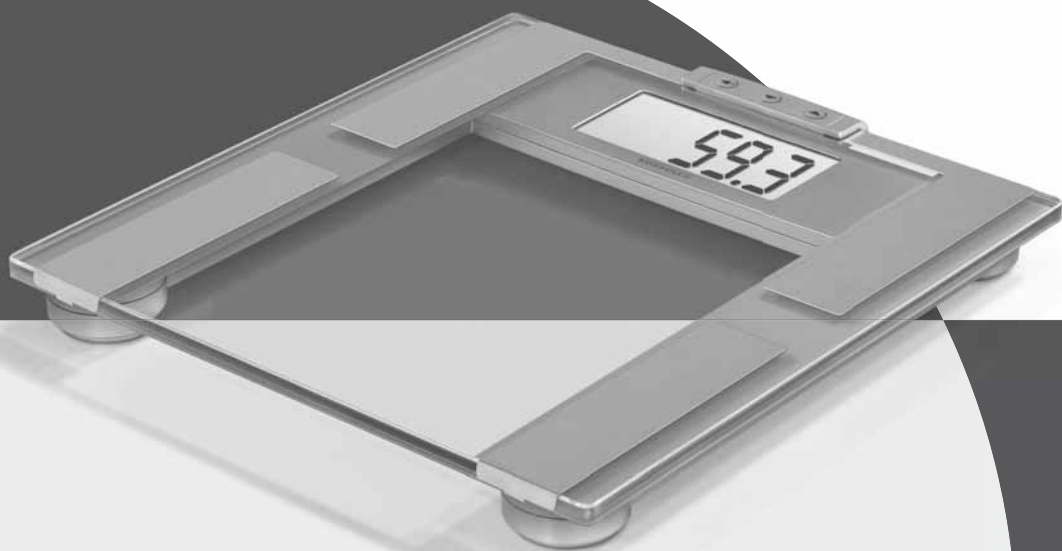
SOEHNLE

Pharo 200 Analytic

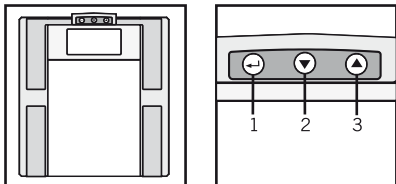
BIO IMPEDANCE
TECHNOLOGY

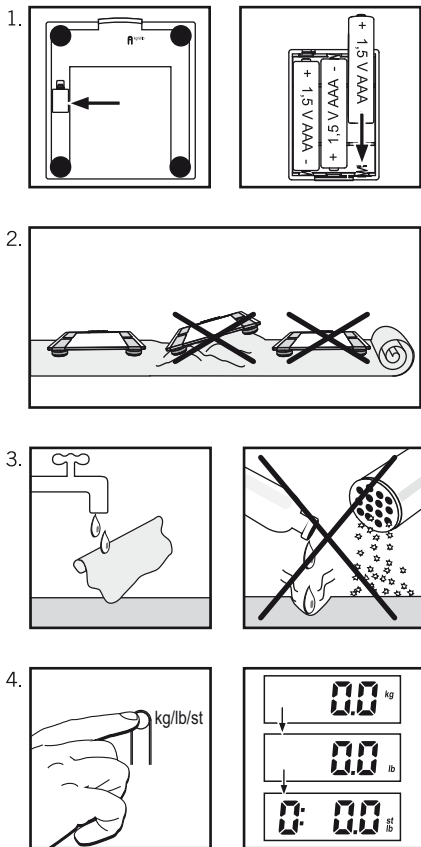
BEDIENUNGSANLEITUNG	3
USER MANUAL	5
MODE D'EMPLOI	7
ISTRUZIONI PER L'USO	9
GEBRUIKSAANWIJZING	11
INSTRUCCIONES DE MANEJO	13
MANUAL DE INSTRUÇÕES	15
BRUKSANVISNING	17
BRUGSANVISNING	19
ΚΑΥΤΤÖΗJΕΕΤ	21
KEZELÉSI ÚTMUTATÓ	23
INSTRUKCJA OBSŁUGI	25
NÁVOD NA OBSLUHU	27
NÁVOD K POUŽITÍ	29
ИНСТРУКЦИ ПО ЭКСПЛУАТАЦИ	31
KULLANMA KILAVUZU	33
ΟΔΗΓΙΕΣ ΛΕΙΤΟΥΡΓΙΑΣ	35
NAVODILA ZA UPORABO	37
UPUTA ZA UPORABU	39

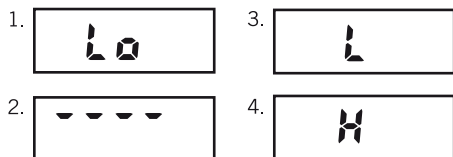
A

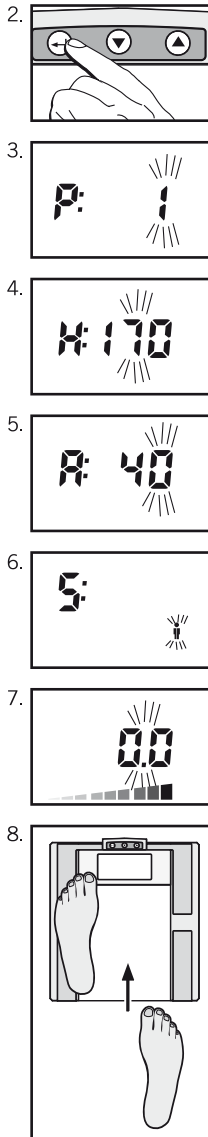
B

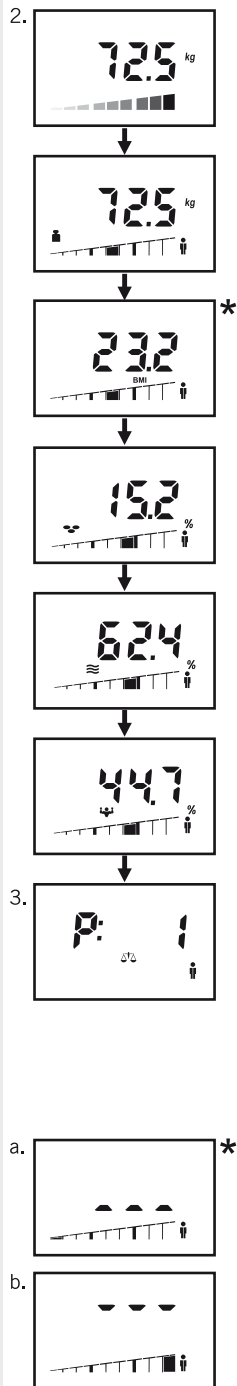
E

C

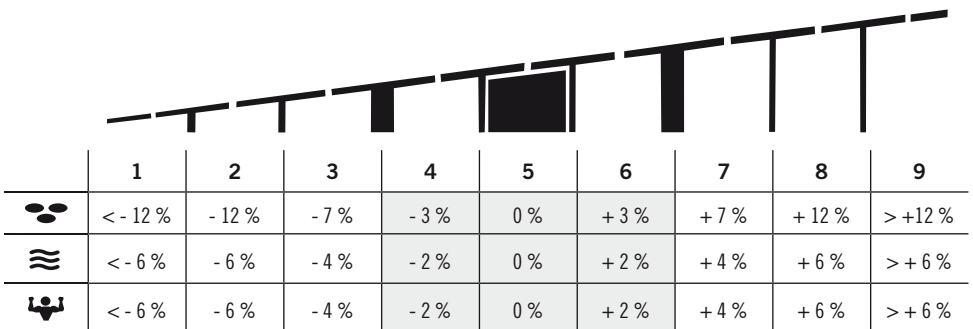
D

- Nicht geeignet für Personen mit elektronischen Implantaten (Herzschrittmacher).
- It is not suitable for persons with electronic implants (pacemaker, etc.).
- Les personnes qui portent des implants médicaux (p. ex. un stimulateur cardiaque etc.) ne doivent pas utiliser cet appareil.
- Non è adatto per persone con impianti elettronici (pacemaker, ecc.).
- Niet geschikt voor personen met elektronische implantaten (pacemaker, enz.).
- No es apropiado para personas con implantes electrónicos (marcapasos, etc.).
- Não é adequado para pessoas com implantes electrónicos (pacemakers, etc.).
- Ej lämplig för personer med elektroniska implantat (Pace-maker etc).
- Tuote ei sovellu henkilöille, joilla on elektronisia istutteita (sydämen tahdistin tms.).
- Nem alkalmazható olyan személyeknél, akik elektronikusan implantátumokkal (például szívritmusszabályozóval stb.) rendelkeznek.
- Nie nadaje się dla osób, którym wszczepiono elektroniczne implanty (np. rozruszniki serca).
- Prístroj nie je vhodný pre osoby s elektronickými implantátmi (kardiostimulátorom).
- Není vhodný pro osoby s elektronickými implantáty (kardiostimulátory, atd.).
- Никогда не используйте данный прибор вместе с электронными медицинскими приборами, например, кардиостимуляторами.
- Elektronik parçalar (kalp pili, vs.) taşıyan insanlar için uygun değildir.
- Δεν είναι κατάλληλο για άτομα με ηλεκτρονικά εμφυτεύματα (βηματοδότες, κ.λπ.).
- Ni primerno za osebe z elektronskimi implantacijami (spodbujevalec srca itd.).
- Nije pogodan za lica sa elektroničnim implantatima (bajpas, itd.).

Abweichung zum Normalwert • Deviation from the normal value • Ecart par rapport à la valeur normale • Deviazione rispetto al valore normale • Afwijking van de normale waarde • Variación con respecto al valor normal • Desvio do valor normal • Afvigelse fra normalværdien • Avvikelser från normalvärdet • Poikkeama normaaliarvosta • Eltérés a normál értékhez képest • Odchylka od normálních hodnot • Odchýlka od normálnej hodnoty • Odchylenie od wartości normalnej • Отклонение от нормального значения • Normal degerlerden sapma • Απόκλιση από την κανονική τιμή • Odstopanja od standardnih vrednosti • Odstupanje od normalne vrijednosti

A. Bedienelemente

1. Bestätigen (↵)
2. Minus (▼)
3. Plus (▲)

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1. Batterien einlegen.

Hinweis:

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B. Preparation

1. Remove battery strip.

Note: Immediately (as long as „0.0“ is still being displayed) place the body analysis scale on a level surface and wait until the scale powers down automatically.

Then begin data input. Otherwise, it might be possible that the weight of the scale would erroneously be added to the initial weighting process.

The same might occur during regular operation if the scale is held in the hand prior to the weighting process.

2. For all measurements, place scale level on a hard surface (not carpet).
3. Cleaning and Care: clean only using slightly damp cloth. Do not use solvents or abrasives. Never submerge scale in water.
4. A country-specific option to switch from kg/cm to lb/in or st/in is possible by pushing the key on the bottom of the scale to switch units.

Attention! Slipping hazard in case of wet surface.

C. Data input

For the correct functionality of the body analysis scale it is necessary to enter the personal data of the respective person and then immediately step on the scale with bare feet.

1. Place scale on a table to enter data.
2. Start data input (push ↵ key for 1 second)
3. Select memory slot (P1 ... P8) (▼ = minus, ▲ = plus) and confirm (↵).
4. Enter body height (▼ = minus, ▲ = plus) and confirm (↵).
5. Enter age (▼ = minus, ▲ = plus) and confirm (↵).
6. Select gender using ▼/▲ (♂ = male, ♀ = female) and confirm (↵).
7. The number of the memory slot is being displayed.

Push ↵ key and hold for **3 seconds** until „0.0“ appears on the screen.

8. Now immediately place the body analysis scale on the floor, wait until the screen displays „0.0“ and then step on the scale with bare feet.

The initial measurement is indispensable and serves to store the data for later automatic recognition of persons. If no person steps on the scale, data input must be repeated.

The first body analysis is performed after a person steps on the scale.

This scale is equipped with a Time-Out mode. It will automatically switch off after approx. 40 seconds if no keys are pressed.

Data must be entered again after battery change.

D. Body analysis

1. For most precise results, please step on the scale while unclothed. The weight and - for the duration of the analysis (5 to 10 seconds) - a moving bar indicator will be displayed. Please stand still until the analysis is complete.

2. After automatic recognition of persons, the scale will again display the body weight, BMI, body fat content (👤), body water content (≍), and muscle mass (👤) along with the corresponding status bar.

The average range for the respective person is displayed by the median 3 segments of the bar indicator.

If the values are outside of the body analysis scale's measurement range, the screen displays the bottom segments (image a) if the values are lower. The upper segments will be displayed (image b) if the values are higher.

The precise measurement ranges are listed under "F. Technical Data".

3. The scale will then display the memory slot and one or more of the following symbols. This means that one or more body values are not within average range.

- 👤 Body weight
- 👤 Body fat content
- ≍ Body water content
- 👤 Muscle mass

If the scale symbol (⚖️) is displayed, your analysis values are within average range!

4. Step off the scale. It switches off automatically.

In rare cases it is not possible to automatically allocate the analysis results. In this case, the scale will alternately display the memory slots of persons coming closest to the analysis result. Since weight measurement is already completed, you may now step off the scale and use the ▲ or ▼ key (▼ key for the first memory slot displayed, ▲ key for the second memory

slot displayed) to confirm your individual memory slot. Following this, body weight, BMI, body fat and body water content, as well as muscle mass will be displayed.

In case of weight differences of more than

+/- 3 kg compared to the last measurement, the user will not be recognized, which means that data input must be repeated.

Please note that when stepping on the scale while wearing footwear or socks, only the weight is being determined and a body analysis is not performed or performed incorrectly!

You will achieve continuously precise results if you always stand in the same position on the scale if possible.

E. Messages

1. Batteries spent.
Insert new batteries
2. Overload: Over 200 kg.
3. Implausible analysis result: review personal programming
4. Improper foot contact: clean scale or feet and ensure sufficient moisture on skin. If skin is too dry, slightly wet feet or perform measurement preferably after bath or shower.

F. Technical data

Load capacity x graduation =
max. 200 kg x 100 g

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International Ph.: +49 26 04 97 70

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Guarantee

Leifheit AG grants you 3 years guarantee as of the date of purchase (or for order upon receipt of the goods) for the quality product at hand. The guarantee claims must be asserted immediately upon occurrence of the defect within the guarantee period. The guarantee is limited to the quality of the products.

The following are excluded from the guarantee:

- (1) wear caused by proper use or natural causes;
- (2) damage caused by improper use or handling (e.g. impact, bump, fall);

(3) damage caused by non-compliance with stated operating instructions;

(4) the rechargeable battery

In case of a guarantee claim, Leifheit will at its own discretion provide either repair of defective parts or replace the product. If a repair is not possible and if an identical product is no longer available within the product range for replacement, you will receive a replacement product that is, insofar as possible, of equal value. In case of a guarantee claim, the purchase price cannot be reimbursed. This guarantee also does not grant compensation rights. This guarantee furthermore does not grant any claims for damages.

To assert the guarantee, please contact the seller where you purchased the product and present the defective product, the purchase receipt (copy) and the guarantee slip.

The guarantee is valid globally.

Your statutory rights, in particular warranty rights, continue to be valid and are not limited by this warranty.

EC Conformity

 This device complies with the applicable EC Directive 2004/108/EC. You will find the EU compliance declaration on www.soehnle.com.

Battery disposal EC Directive 2008/12/EC

Batteries are not a part of your regular household waste. You must return batteries to your municipality's public collection or wherever batteries of the respective type are being sold.

Batteries containing hazardous materials are marked with these symbols:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

Disposal of electric and electronic devices

EC Directive 2002/96/EC

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

Consumer service

If you have questions or comments, the following points of contact are available:

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LEIFHEIT AG