
# SOEHNLE

## SOLAR FIT

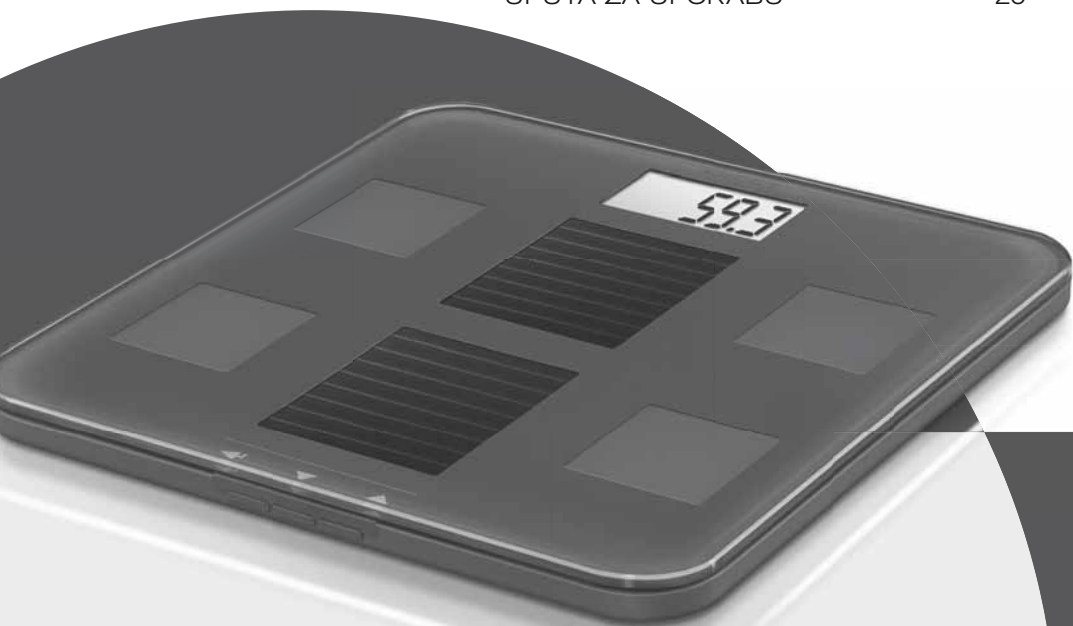BEDIENUNGSANLEITUNG	2
OPERATING INSTRUCTIONS	3
MODE D'EMPLOI	4
ISTRUZIONI PER L'USO	5
GEBRUIKSAANWIJZING	6
INSTRUCCIONES DE MANEJO	7
MANUAL DE INSTRUÇÕES	8
BRUKSANVISNING SVENSKA	9
BRUGSANVISNING DANSK	10
BRUKSANVISNING NORSK	11
KÄYTTÖOHJEET	12
KEZELÉSI ÚTMUTATÓ	13
INSTRUKCJA OBSŁUGI	14
NÁVOD K POUŽITÍ	15
ИНСТРУКЦИИ ПО ЭКСПЛУАТАЦИИ	16
KULLANMA KILAVUZU	17
ΟΔΗΓΙΕΣ ΛΕΙΤΟΥΡΓΙΑΣ	18
NAVODILA ZA UPORABO	19
UPUTA ZA UPORABU	20

Alter/Age	 					 				
	-	✓	+	++	✓	-	✓	+	++	✓
10-12	< 12	12-23	23-30	> 30	> 60	< 8	8-18	18-24	> 24	> 64
12-18	< 15	15-25	25-33	> 33	> 58,5	< 8	8-18	18-24	> 24	> 63,5
18-30	< 20	20-29	29-36	> 36	> 56	< 8	8-18	18-24	> 24	> 62,5
30-40	< 22	22-31	31-38	> 38	> 53	< 11	11-20	20-26	> 26	> 61
40-50	< 24	24-33	33-40	> 40	> 52	< 13	13-22	22-28	> 28	> 60
50-60	< 26	26-35	35-42	> 42	> 51	< 15	15-24	24-30	> 30	> 59
60 +	< 28	28-37	37-47	> 47	> 50	< 17	17-26	26-34	> 34	> 58


Körperfettanteil, body fat content, taux de graisse, percentuale di grasso, lichaamsvetaandeel, procentaje de la grasa corporal, procentagem de gordura corporal, fedtprocent, andelen kroppsfett, Fettprozent, ehon rasvapitoisuutesi, test zsírtartalmát, zawartość tłuszczu w ciele, podíl tělesného tuku, Содержание жира в организме, үцүттәки yağ oranını, το ποσοστό λίπους, delež telesne maščobe, udjel tjelesne masnoće.





Körperwasseranteil, body water content, taux d'eau dans le corps, percentuale acqua corporei, lichaamsweteraandeel, procentaje del agua corporal, procentagem de água corporal, vandprocent, andelen kroppsvatten, Body vannprosent, kehon nestepitoisuutesi, test víztartalmát, zawartość wody w ciele, podíl vody v těle, Содержание воды, үцүттәки su oranını, το ποσοστό νερού, delež telesne vode, udjel tjelesne tekućine.




# A




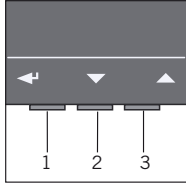
# B

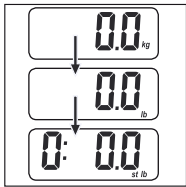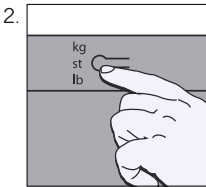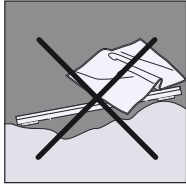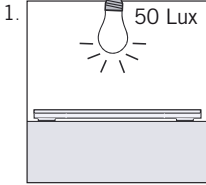
# C

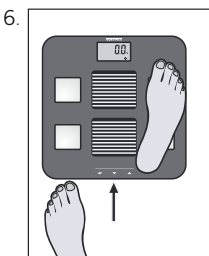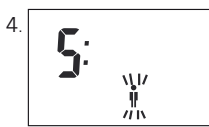
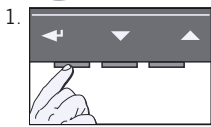
# D
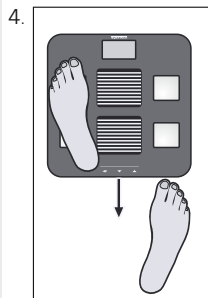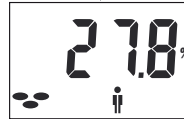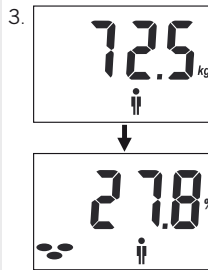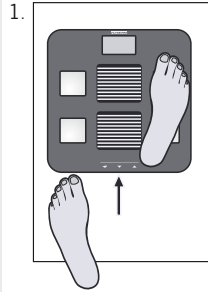
# E
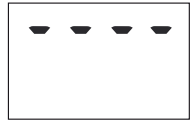
# F
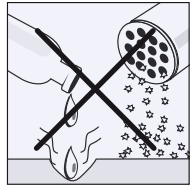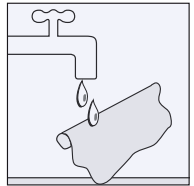
## A. Bedienelemente

1. Bestätigen (◀)
2. Minus (▼)
3. Plus (▲)

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1. Confirm (◀)
2. Minus (▼)
3. Plus (▲)

## B. Preparation

1. For all measurements place scale level on a hard surface (not carpet). Ensure sufficient illumination. Do not cover solar cells.
2. A country-specific option to switch from kg/cm to lb/in or st/in is possible by pushing the key on the bottom of the scale to switch units.

**Attention!** Slipping hazard in case of wet surface.

## C. Data input

**Note:** Place the body analysis scale on a level surface. Then begin with data entry. Otherwise, it might be possible that the weight of the scale would erroneously be added to the initial weighting process.

The same might occur if the scale is held in the hand prior to the weighting process.

1. Start data entry (push ◀ - key for 1 second)
2. Select memory slot (P1 ... P8) (▼ = minus, ▲ = plus) and confirm (◀).
3. Enter body height and age (▼ = minus, ▲ = plus) and confirm (◀).

4. Select gender using ▼/▲ (♂ = male, ♀ = female) and confirm (◀).

5. The number of the memory slot is being displayed. Push ◀ key.

0.0 appears on the screen.

6. Then immediately step on the scale with **bare feet** for the first body analysis. If you step on the scale with socks or shoes, then only the body weight will be displayed.

The initial measurement is indispensable and serves to store the data for later automatic recognition of persons. If no person steps on the scale, data entry must be repeated.

This scale is equipped with a Time-Out mode. It will automatically switch off after approx. 40 seconds if no keys are pressed.

## D. Body analysis

1. Step on the scale with bare feet; stand still for approx. 5 sec.
2. Automatic recognition of persons. Select personal memory slot (▼/▲) if person cannot be automatically recognized.
3. Display weight, body fat content, body water content, and personal memory slot.

**Note:** Body fat and water are calculated based on the determined weight and the personal data.

4. Step off the scale.

## E. Messages

1. Overload: Over 150 kg.

## F. Cleaning and Care

Clean only using slightly damp cloth. Do not use solvents or abrasives. Never submerge scale in water.

## G. Technical data

Load capacity x graduation = max. 150 kg x 100 g

Graduation of body fat content: 0.1%

Graduation of body water content: 0.1%

Body height: 100-250 cm

Age: 10-99 years

8 individual memory slots

## H. Garantie

Leifheit AG grants you 5 years guarantee as of the date of purchase (or for order upon receipt of the goods) for the quality product at hand. The guarantee claims must be asserted immediately upon occurrence of the defect within the guarantee period. The guarantee is limited to the quality of the products.

The following are excluded from the guarantee:

- (1) wear caused by proper use or natural causes;
- (2) damage caused by improper use or handling (e.g. impact, bump, fall);
- (3) damage caused by non-compliance with stated operating instructions;
- (4) the rechargeable battery

In case of a guarantee claim, Leifheit will at its own discretion provide either repair of defective parts or replace the product. If a repair is not possible and if an identical product is no longer available within the product range for replacement, you will receive a replacement product that is, insofar as possible, of equal value. In case of a guarantee claim, the purchase price cannot be reimbursed. This guarantee also does not grant compensation rights. This guarantee furthermore does not grant any claims for damages.

To assert the guarantee, please contact the seller where you purchased the product and present the defective product, the purchase receipt (copy) and the guarantee slip. The guarantee is valid globally.

Your statutory rights, in particular warranty rights, continue to be valid and are not limited by this warranty.

## I. EC Conformity

This device complies with the applicable EC Directive 2004/108/EC.

## K. Disposal of electric and electronic devices EC Directive 2002/96/EC

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

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LEIFHEIT AG