


# taurus

## Dental Care

Cepillo Eléctrico Recargable  
Raspall Elèctric Recarregable  
Rechargeable Electric Toothbrush  
Brosse à dent Électrique Rechargeable  
Elektrische Zahnbürste Aufladbar  
Spazzolino Elettrico Ricaricabile  
Escova Eléctrica Recarregável  
Oplaadbare Elektrische Tandborstel  
Elektryczna Szczoteczka do Zębów  
Επαναφορτιζόμενη ηλεκτρική  
οδοντόβουρτσα  
Электрическая зубная щетка на  
аккумуляторе  
Perie Electrică Reîncărcabilă  
Презареждаща се ел.четка за зъби  
فرشاة كهربائية قابلة للشحن

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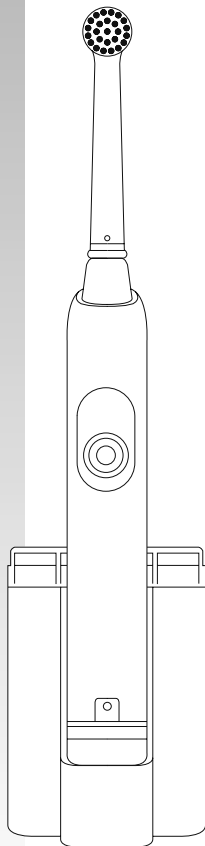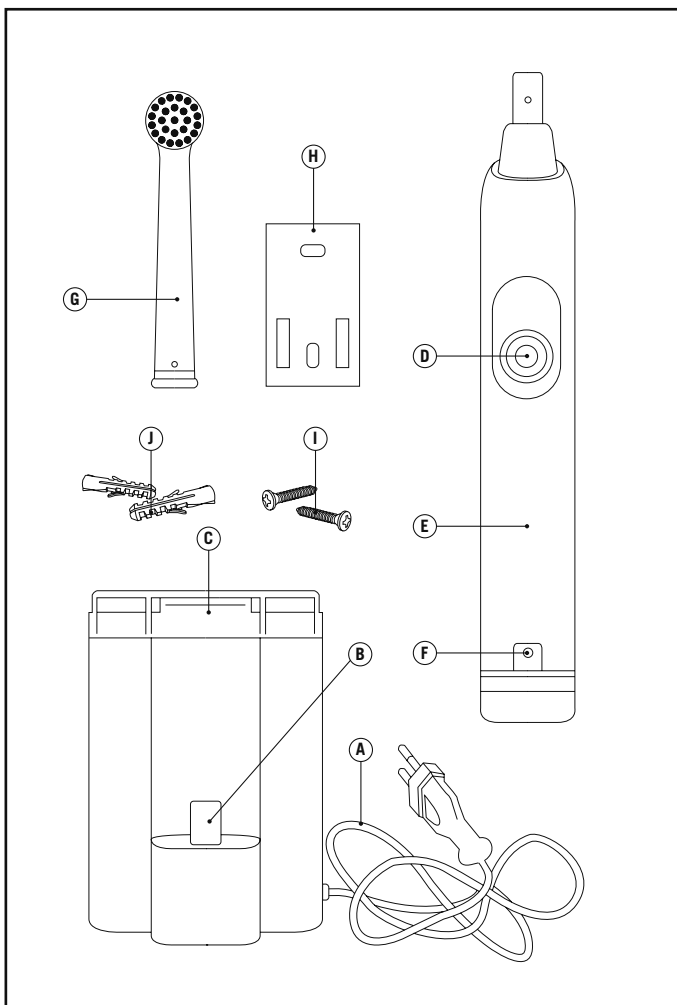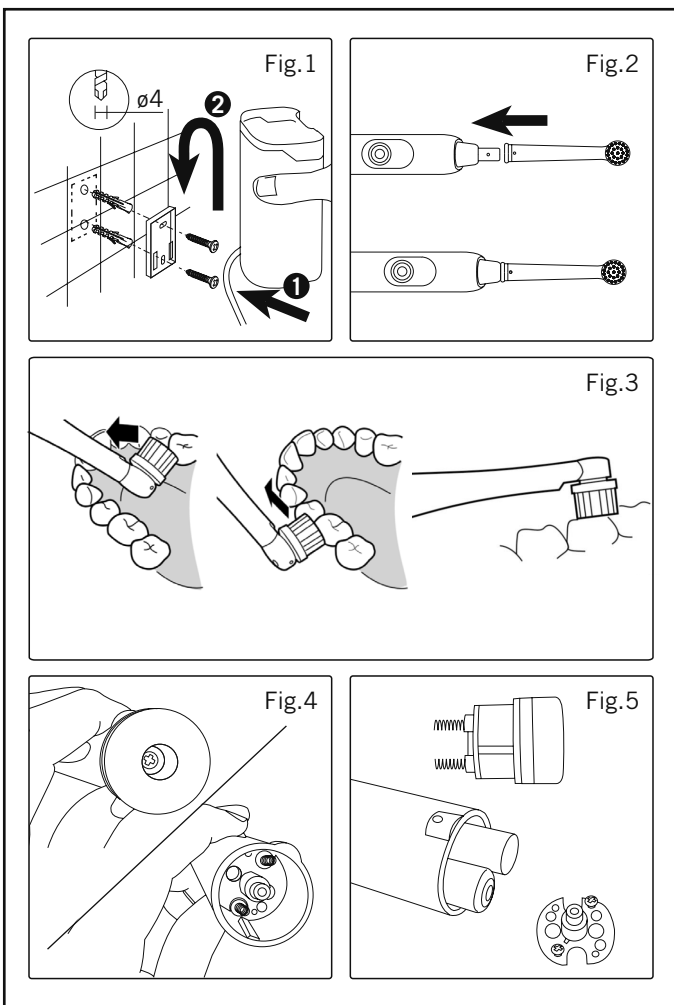## Rechargeable Electric Toothbrush Dental Care

### Dear customer,

Many thanks for choosing to purchase a TAURUS brand product.

Thanks to its technology, design and operation and the fact that it exceeds the strictest quality standards, a fully satisfactory use and long product life can be assured.

### Description

- A Electric cable
- B Charger
- C Toothbrush container
- D On/off switch
- E Top housing
- F Charge light
- G Detachable head
- H Wall mount
- I Screws
- J Raw plugs

If your toothbrush model does not have the accessories outlined above, they can be acquired separately from the Technical Assistance Service.

### Safety advice and warnings

- Read these instructions carefully before switching on the appliance and keep them for future reference. Failure to follow and observe these instructions may result in an accident.

### Use or working environment:

- The appliance must only be used with the stand that is provided.
- The appliance should be used and kept over a flat and stable surface.
- Make sure that the charger is placed on a stable surface before putting the appliance on the stand.
- The appliance must be used and placed on a flat, stable surface.

### Electric safety:

- Do not use the appliance if the cable or plug is damaged.
- Ensure that the voltage indicated on the charger rating label matches the mains voltage before plugging in the charger.
- Connect the appliance to a base with an earth socket withstanding a minimum of 10 amperes.
- Only use the appliance with the specific electric base / connector, supplied with the appliance.
- If any of the appliance casings breaks, immediately disconnect the appliance from the mains to prevent the possibility of an electric shock.
- Do not use the appliance if it has fallen on the floor, if there are visible signs of damage or if it has a leak.
- WARNING: Keep the appliance dry.
- Do not use the appliance with damp hands or feet, or when barefooted.
- Do not force the power cord. Never use the power cord to lift up, carry or unplug the battery charger.
- Do not wrap the cable around the appliance.
- Check the condition of the electrical connection cable. Damaged or tangled cables increase the risk of

electric shock.

- As an additional protection to the electrical supply for the appliance, it is advisable to have a differential current device with a maximum sensitivity of 30mA.

Ask for advice from a competent installer.

- Do not touch the plug with wet hands.

**Personal safety:**

- Avoid touching liquid that may have leaked from the battery. If this liquid comes into contact with your eyes, wash thoroughly and seek medical assistance from a doctor. Leaked battery liquid may cause irritation or burns.

Use and care:

- Fully unroll the appliance's power cable before each use.

- Do not use the appliance if the accessories or parts are not fastened together properly.

- Do not use the appliance if the on/off switch does not work.

- This appliance is for household use only, not professional or industrial use.

- This appliance can be used by people unused to its handling, disabled people or children aged from 8 years and above if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

- Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised

- Keep the appliance and its cable

out of reach of children under the age of eight.

- This appliance is not a toy.

Children should be supervised to ensure that they do not play with the appliance.

- This appliance should be stored out of reach of children and/or the disabled.

- Keep the appliance in a dry, dust-free place, out of direct sunlight.

- Keep the appliance in good condition. Check that the moving parts are not misaligned or stuck, that there are no broken parts, or other conditions that may affect the correct functioning of the appliance.

- Use this appliance, its accessories and tools in accordance with these instructions, while taking into account the conditions in which it is being used and what it is being used for. Using the appliance for anything other than the stipulated use may provoke hazardous situations.

- Never rest the appliance into a surface while in use.

- Do not use the appliance with pets or animals.

- Recharge the battery/ies only with the charger supplied by the manufacturer.

A charger that is suitable for one type of battery may cause explosion or fire hazard when used with other battery/ies.

- Only use the appliance with the battery/ies that have been designed specifically for this purpose. Using any other type of battery may cause explosions or a fire hazard.

### Service:

- It is recommended that the detachable head is replaced every 3 months in order to obtain optimum brushing results.
- If the bristles are misshapen in less than 3 months it is possible that the appliance is being pressed too hard on the teeth and gums during brushing.
- Any misuse or failure to follow the instructions for use renders the guarantee and the manufacturer's liability null and void.

### Charging the battery

The new battery does not come completely charged and must be charged up to the maximum before using the appliance for the first time.

- Charge the battery at room temperature between 4 and 40°C.
  - Plug in the charger to the mains.
  - The time needed to charge completely is from 16 to 20 hours approximately. Never try to overcharge the batteries (although it is recommended that the batteries are charged for 18 hours for the first time).
  - Place the appliance on the charger.
  - While the battery is charging, the pilot light will come on indicating that the charger is charging the battery.
- Note: Do not use the appliance while it is connected to the charger.

### Wall mount (Fig 1):

- Mount the charger on the wall using the screws and raw plugs included with the appliance.
- Mark points on the wall for drilling, taking into account the distance between

- the holes on the wall mount.
- Drill holes and insert raw plugs.
- Fit the wall mount making sure that the holes coincide with the raw plugs and then screw in place.
- Make sure that the wall mount is attached correctly so that it can support at least double the weight of the appliance.
- Place the appliance on the wall stand

### Instructions for use

#### Before use:

- Before using this product for the first time, it is important to make a full charge of batteries.

#### Insertion of detachable head:

- Place the detachable head on the top housing, making sure that the brush is on the same side as the on/off switch (Fig. 2)
- Pull detachable head upwards to remove

#### Use:

- The appliance can be used by plugging into the mains, or with batteries
- Unroll the cable completely before plugging it in.
- Place the appliance correctly over the charger.
- Connect the appliance to the mains.
- Charge up the batteries fully.
- Before using the toothbrush rinse out your mouth and wet the toothbrush
- Put toothpaste on the toothbrush.

This appliance can be used with any type of toothpaste or mouthwash

- To avoid splashing, place the toothbrush in your mouth before switching on.
- Switch the appliance on, by pressing the on/off switch, and brush teeth.

#### **Brushing techniques:**

- Place the toothbrush in your mouth in a parallel position to your teeth and move slowly, holding it for a few seconds over each tooth.
- Brush the outside of your teeth gently, followed by the inside of your teeth and finally the top chewing surfaces, as well as the gum line, without pressing too hard. Use rotating movements to brush teeth. (Fig 3)
- A brushing time of 2 minutes is recommended to completely remove plaque.
- During the first days of use it is possible that light bleeding of the gums may occur, if you suffer from this condition normally. This problem is generally resolved after a few days. If it continues for more than 2 weeks, seek assistance from a dentist.
- Once you have finished brushing your teeth, take the appliance out of your mouth before switching off. After brushing, rinse your mouth thoroughly again.

#### **Once you have finished using the appliance:**

- Turn the appliance off, using the on/off switch.
- Clean the appliance

#### **Cleaning**

- To keep the detachable head in perfect working condition, rinse under water for a few seconds after brushing with the toothbrush turned on.
- Switch off the appliance and remove the detachable head. Rinse both parts under running water and dry.
- The brush may change shape due to use as a result of the user's individual brushing technique, but this shouldn't affect the effectiveness of the toothbrush.
- It's recommended that the detachable head is changed every 3 months, depending on use, to obtain optimum brushing results.
- If the appliance is not in good condition of cleanliness, its surface may degrade and inexorably affect the duration of the appliance's useful life and could become unsafe to use.
- Dry all parts before assembling and storing

#### **Supplies**

- Supplies can be acquired from distributors and authorised establishments (detachable heads) for your appliance model.
- Always use original supplies, designed specifically for your appliance model.
- These supplies are available in specialist shops.

#### **Anomalies and repair**

- Take the appliance to an authorised technical support service if problems

arise. Do not try to dismantle or repair without assistance, as this may be dangerous.

- In case of breakdown call the customer service telephone number indicated in the guarantee sheet. Do not attempt to disassemble or repair the appliance yourself as this may be dangerous.

- If the connection to the mains has been damaged, it must be replaced and you should proceed as you would in the case of damage.

For EU product versions and/or in case that it is requested in your country:

#### **Ecology and recyclability of the product**

- The materials of which the packaging of this appliance consists are included in a collection, classification and recycling system. Should you wish to dispose of them, use the appropriate public recycling bins for each type of material.

The product does not contain concentrations of substances that could be considered harmful to the environment

- This appliance contains, or could contain, oil. Remember that one should abide by the regulations concerning the elimination of this type of contaminant.

This symbol means that in case you wish to dispose of the product once its working life has ended, take it to an authorised waste agent for the selective collection of waste electrical and electronic equipment (WEEE)

This symbol means that the product may include a battery or batteries; the user must remove them before disposing of the product. Remember that batteries must be disposed of in duly authorised containers. Do not dispose of them in fire.

#### **How to remove the batteries from the interior of the device:**

It is important that the batteries are fully discharged prior to its withdrawal from the apparatus.

#### **To remove the batteries, once the life of the appliance, proceed as follows:**

- Open the top shell of the apparatus, removing previously with the help of a screwdriver fixing screws (Fig 4)
- Remove the following part by unscrewing the screws, with the help of a screwdriver (Fig. 5)
- Safely remove the battery from its housing.

This appliance complies with Directive 2006/95/EC on Low Voltage and Directive 2004/108/EC on Electro-magnetic Compatibility.

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