Family


Coach

Terraillon®

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
www.terraillon.com


A

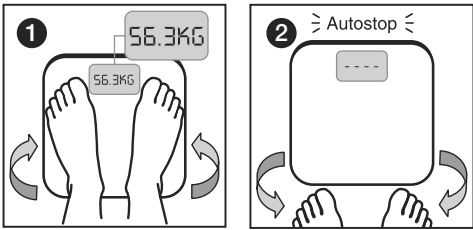
B.1


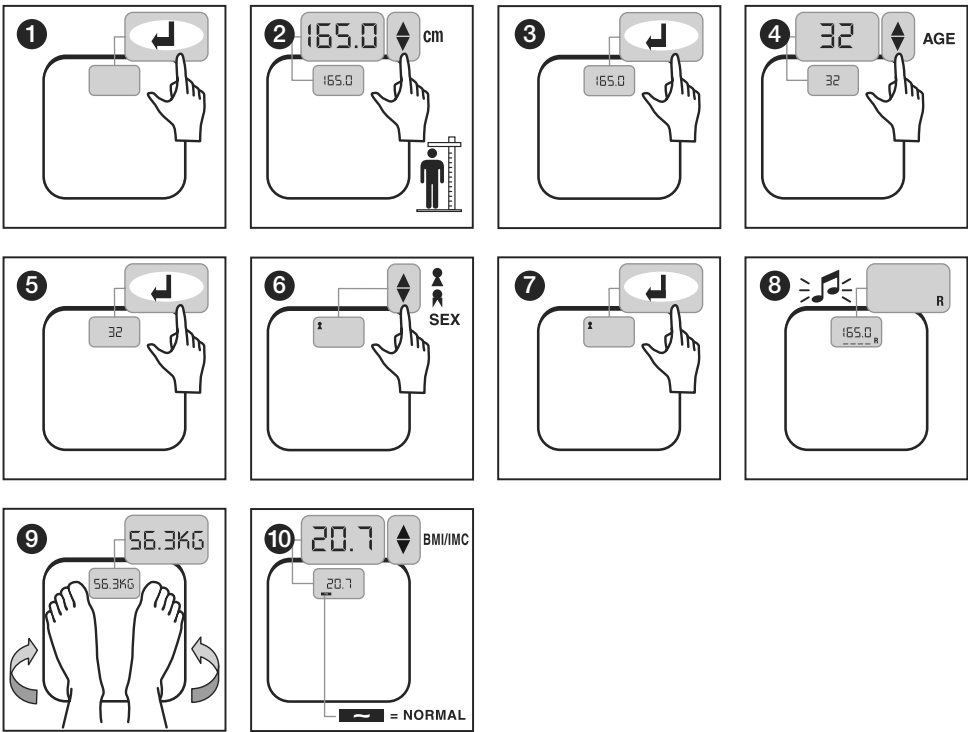
B.2




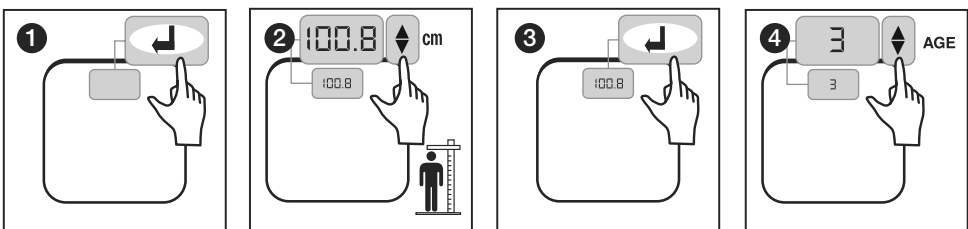
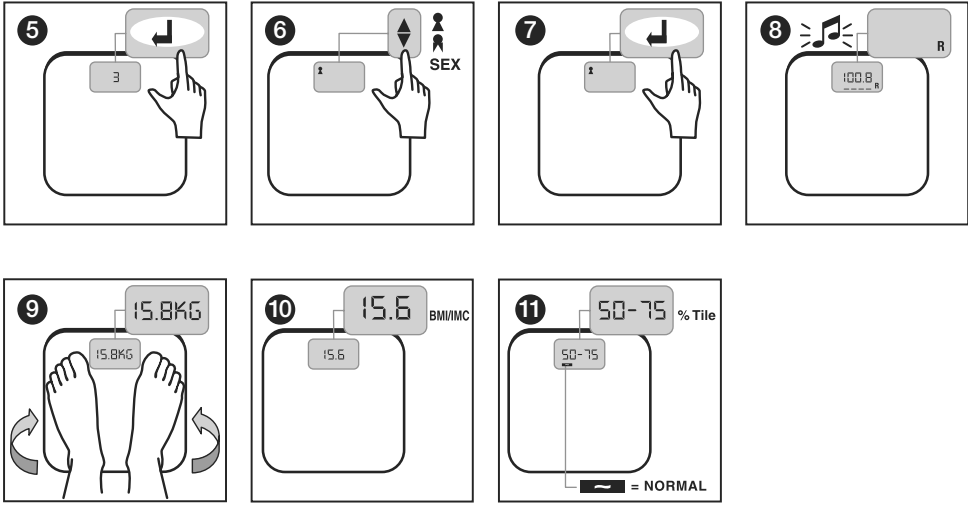
C.1


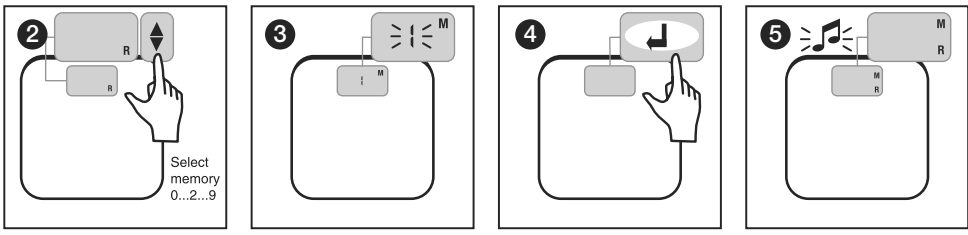
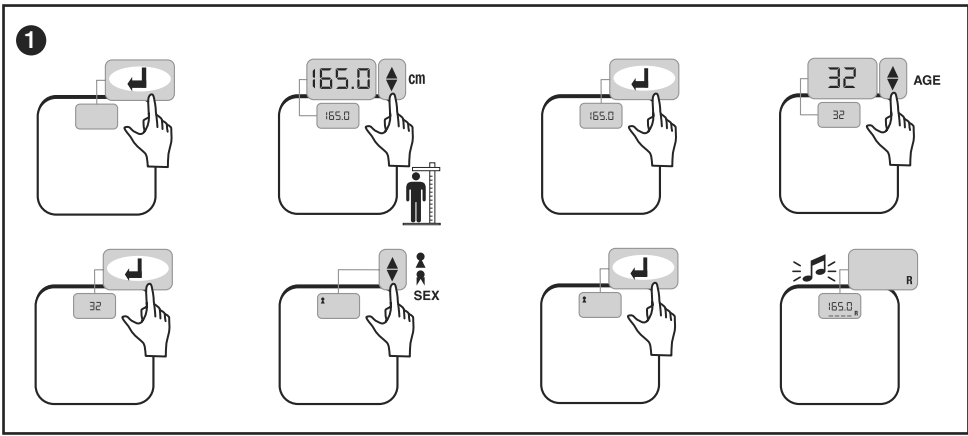
C.2

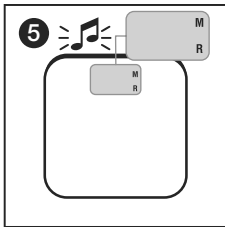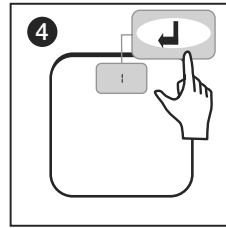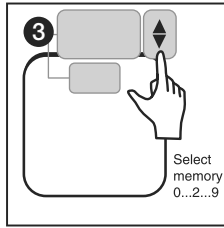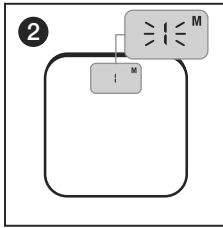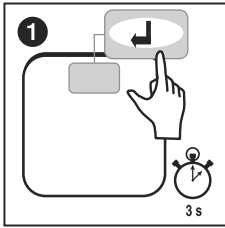cf. B.1.9 → B.1.10
B.2.9 → B.2.11

C.3



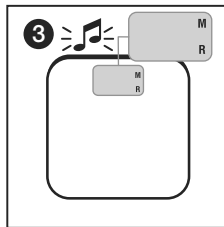
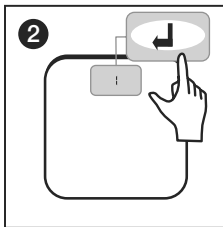
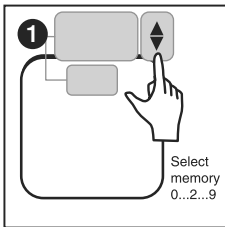
cf. B.1.9 → B.1.10
B.2.9 → B.2.11

D.1


1






1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

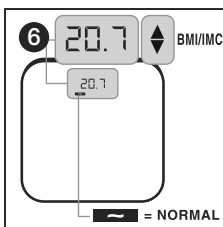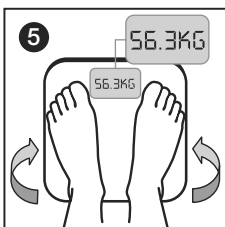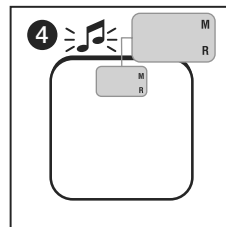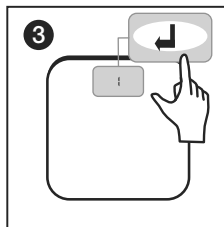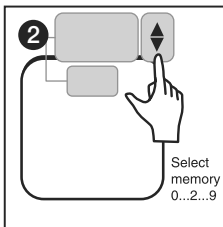7

1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

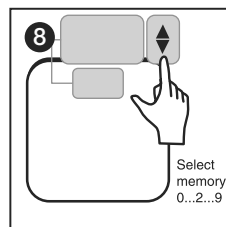
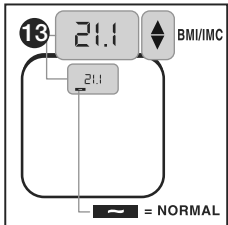
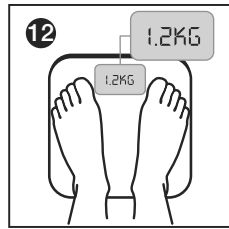
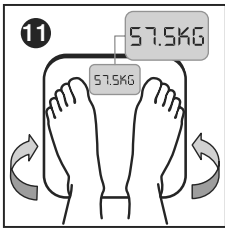
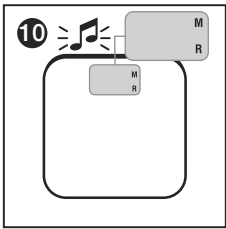
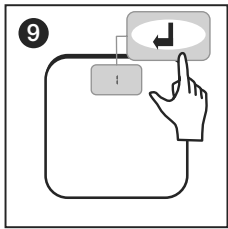
D.2








1

1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

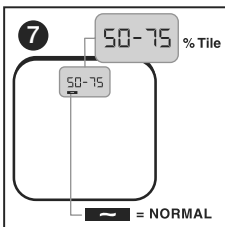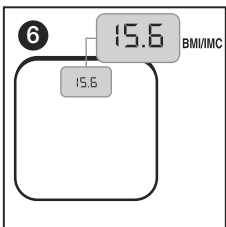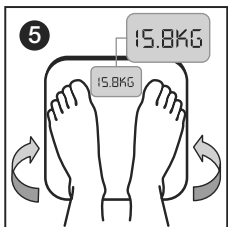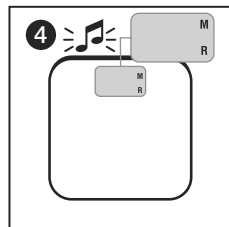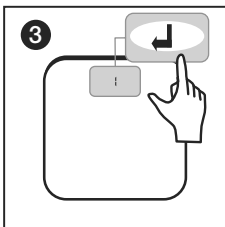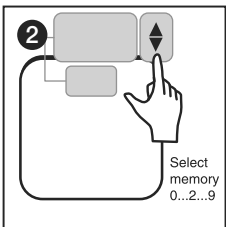8

1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

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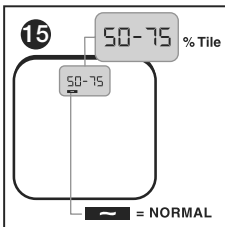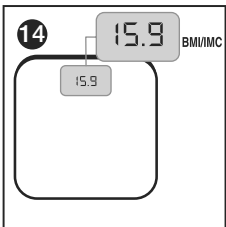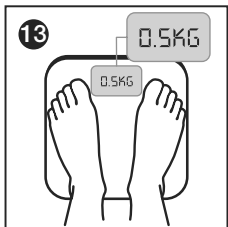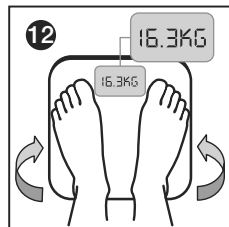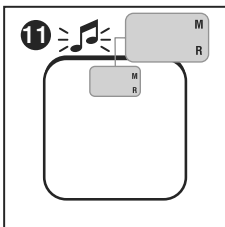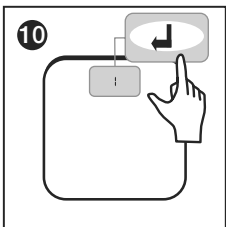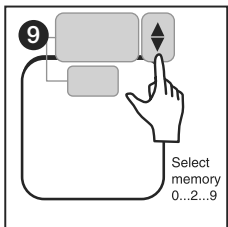ELECTRONIC BATHROOM SCALE FAMILY COACH

WHAT IS THE FAMILY COACH SCALE?

It is a bathroom scale which gives an indication of the body fat of all the members of the family.

The Family Coach scale indicates :

- Your weight,
- Your weight gain or loss compared to your previously saved weight,
- Your Body Mass Index (BMI),
- Your BMI classification according to 4 categories.

Who should use the Family Coach scale?

It is designed for the whole family:

- Children and teenagers from 2 to 18 years
- Adults (from 19 to 99 years*).

** For adults over 65, the BMI classification can sometimes be inaccurate. However, using the Family Coach does not in any way represent a health risk.*

What is BMI (Body Mass Index)?

- For adults:

It is the ratio between a person's weight and height. This index is used to assess a person's body fat.

Calculation: $BMI = \text{weight (kg)} / [\text{height (m)}]^2$

- For children and adolescents:

A child's body fat changes as he/she grows and differs at puberty according to gender.

For this reason, BMI is a function of age and gender. A percentile is then calculated enabling BMI to be classified.

How to use the Family Coach scale?

The Family Coach scale is very easy to use:

- Enter your personal data into the scale (height, age and gender),
- Your BMI (Body Mass Index) will then be calculated according to your weight,
- The scale classifies your BMI according to 4 categories:

- For adults:

> 30	Between 25 and 29.9	Between 18.5 and 24.9	< 18.5
Obese	Overweight	Normal	Underweight

- For children and adolescents:

> 95 P	Between 85 and 95 P	Between 5 and 85 P	< 5 P
Overweight	Risk of overweight	Normal	Underweight

The results bar indicates which category your BMI is in and offers you a clear and immediate result :

Children/Adolescents
Adults

Underweight
Underweight

Normal
Normal

Risk of Overweight
Overweight

Overweight
Obese

ELECTRONIC BATHROOM SCALE FAMILY COACH

Dear Customer,

Thank you for purchasing this product. We hope you will find it satisfactory. Please read these instructions carefully to take full advantage of your product.

READ BEFORE USE

- Do not disassemble the scale as incorrect handling may cause injury.
- The product is not intended for professional usage in hospitals or medical facilities; it is intended for home/consumer use only.
- Clean after use with a lightly dampened cloth without allowing any water to get into the device. Do not use solvents or immerse the product in water.
- Avoid excessive impact or vibration to the unit, such as dropping the product onto the floor.
- The scale must be reset each time it has been moved.
- Pause after each weighing.

PRODUCT SPECIFICATION

Function Keys

3-BUTTON MODEL		4-BUTTON MODEL	
←	<ul style="list-style-type: none"> • Turns the unit on and off. • Enables you to confirm and obtain settings such as height, age and sex. • Displays the 3 criteria (height, age and gender) saved in each memory slot. 	←	<ul style="list-style-type: none"> • Turns the unit on and off. • Enables you to confirm and obtain settings such as height, age and sex. • Displays the 3 criteria (height, age and gender) saved in each memory slot.
		RESET	<ul style="list-style-type: none"> • Clears contents of the user memory in quick recall mode.
▼/▲	<ul style="list-style-type: none"> • Changes the height and age values; toggles between male and female. • Selects memory in Memory mode & Recall mode. • Clears contents of the user memory in Recall mode. 	▼/▲	<ul style="list-style-type: none"> • Changes the height and age values; toggles between male and female. • Selects memory in Memory mode & Recall mode.

Screen Display

	Height	%	Body Mass Index Percentile (for children and adolescents)
	Weight	<	Thin/Underweight
2 - 18	Age for children and adolescents	~	Normal
19-99	Age for adults	>	Overweight/Risk of Overweight
BMI	Body Mass Index (for adults)	>>	Obese/Overweight

ELECTRONIC BATHROOM SCALE FAMILY COACH

PRODUCT USAGE

Initialisation

- Insert the battery(s) in the battery compartment respecting the polarities or pull the red tape to initialise the batteries for the first time. Remove the battery(s) if the product is not in use for an extended period.
- The scale needs to be initialised at the first time of use or after battery replacement. Place the scale on a hard, flat surface. Press briefly and sharply with your foot in the middle of the platform. The display reads OOOO or - - - -, then the scale automatically turns off. Your scale is now ready for use.

Note: it is recommended to take weight a constant time period daily for example: from 6 am to 7 am daily.

Operation

Please refer to the illustrations at the beginning of the instruction manual.

A: WEIGHT CALCULATION ONLY


- Stand on the platform and keep still while the device measures your weight.
- Your weight will be displayed on the screen.
- The scales switch off automatically after use.

B: CALCULATION AND CLASSIFICATION OF BMI

You must first configure your scales by indicating your height, age and sex.

The configuration process is the same for both adults and children:

- Press the button to turn on the scales.
- Press the or buttons to enter your height and the button to confirm this setting.
- Press the or buttons to enter your age and the button to confirm this setting.
- Press the or buttons to enter your sex and the button to confirm this setting.
- Stand on the scales once you hear the confirmation beep.

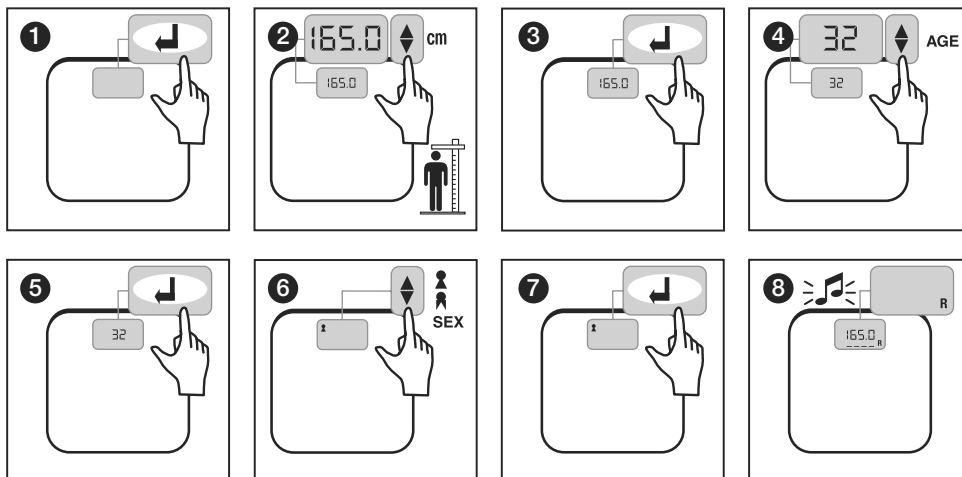
ELECTRONIC BATHROOM SCALE FAMILY COACH

B.1 Indication of results for an adult

The screen first displays your weight then your BMI.


The icon on the bottom of the screen indicates the category your BMI is situated in.

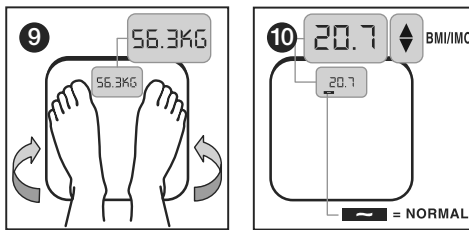B.2 Indication of results for children and teenagers

The screen displays your weight, your BMI and your BMI percentile, in this order.

The icon on the bottom of the screen indicates the category your percentile is situated in.

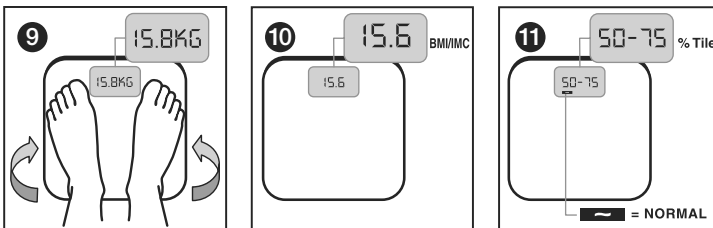
ELECTRONIC BATHROOM SCALE FAMILY COACH

C: MEMORY

C.1 Memorising settings

- First enter your data - height, age and sex (follow the instructions in paragraph B).
- When you hear a beep, select a memory slot by pressing the [▲] or [▼] buttons.
- Press the [↵] button to confirm your choice. A beep will sound.

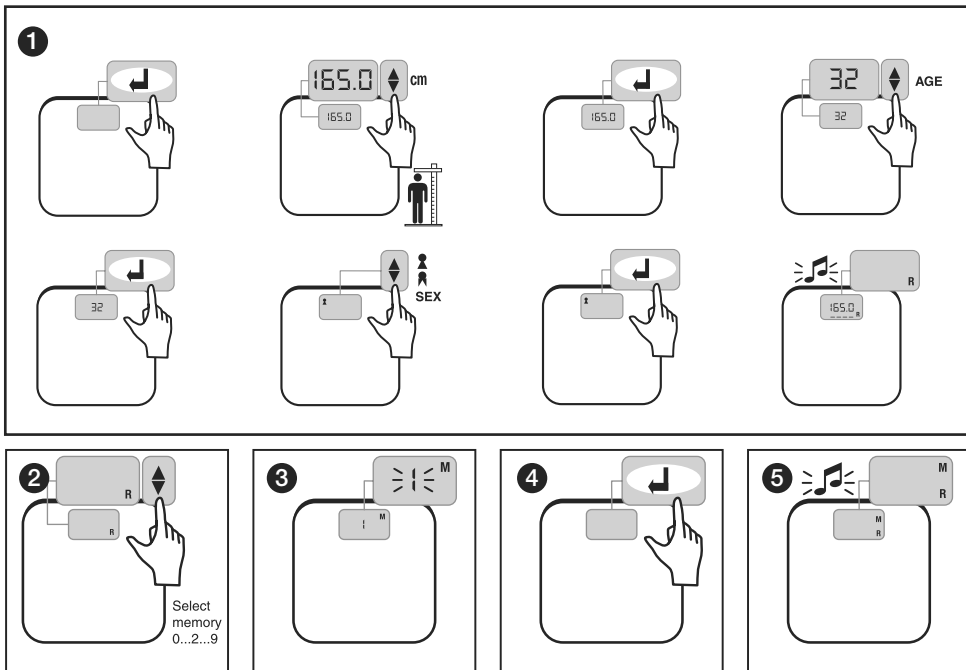C.2 Memory recall

- Press the [↵] button for 3 seconds.
- The figure 1 will flash on the screen. Press the [▲] or [▼] buttons to go to the desired memory slot and the button to confirm your choice.
- Stand on the scales once you hear the confirmation beep.

C.3 Quick memory recall

- With the device switched off, press the [▲] or [▼] buttons to go to the desired memory slot and the [↵] button to confirm your choice.
- Stand on the scales once you hear the confirmation beep.

ELECTRONIC BATHROOM SCALE FAMILY COACH

D: MONITORING CHANGES IN YOUR WEIGHT

The process is the same for adults, children and teenagers. Only the screens displaying the results are different.

When weighing yourself for the first time

- First memorise your data (follow the instructions in paragraph C.1).
- Stand on the scales once you hear the confirmation beep.
- The scales will indicate your weight as well as your BMI (plus the BMI percentile for children and teenagers).
- The icon on the bottom of the screen indicates the category your BMI is situated in.

When weighing yourself subsequently

- Call up your memory slot (follow the instructions in paragraph C.3).
- Stand on the scales once you hear the confirmation beep.
- The scales will indicate your current weight plus the difference between your last weight measured (a "-" sign in front of the figure indicates weight loss).
- The scales will then indicate your BMI (plus the BMI percentile for children and teenagers).
- The icon on the bottom of the screen indicates the category your BMI is situated in.

Warning indicators

Err0	Incorrect initialisation - The scale must be reset.
Err1	Unstable weight: start again.
Err2	Overload: remove weight immediately.
Lo	Battery low: please replace it.
	Memory Full. Press (▼/▲) or RESET to clear the data.

PROTECT THE ENVIRONMENT

Dispose of used batteries in a designated container so that they can be collected and recycled. Do not mix different types of batteries. Do not mix new and used batteries.

After life, dispose of this product at a designated waste recycling point.

GUARANTEE

This product is guaranteed against defects in materials or manufacturing. During the guarantee period, any such defects will be repaired free of charge (proof of purchase must be shown in the case of a complaint under guarantee). The guarantee does not cover damage resulting from improper use or negligence. In the event of a claim, first contact the shop where you purchased the product.

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Terrailon UK Ltd
Terrailon Deutschland GmbH
Terrailon Asia Pacific Ltd
Terrailon Corp