★ ★ ★ ●  
**TRISTAR**

**WG - 2419**

**Opmerkingen:**

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

\*\*\*\*\*


### Notes:

The scale is intended for scaling body weight and household use only.

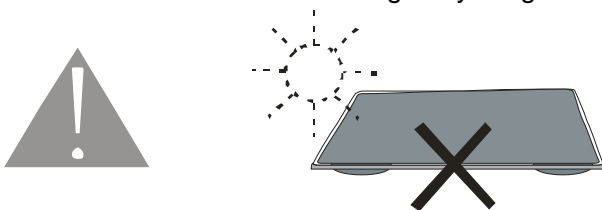
The optimum operating temperature is 20degrees Celsius;  $\pm 10\%$

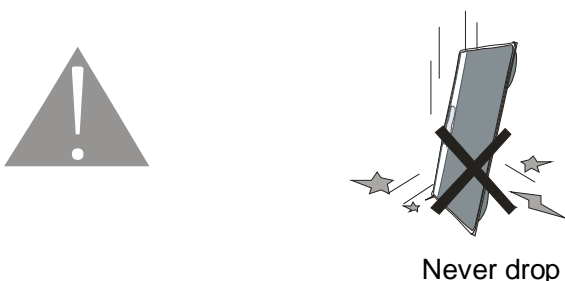
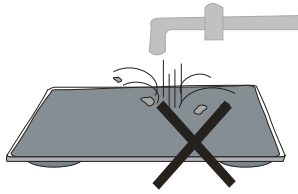
Never wash it with corrosive liquid

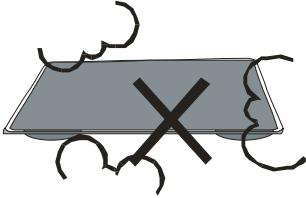
Never use in pollution or high EMI (ElectroMagnetic interference) areas.

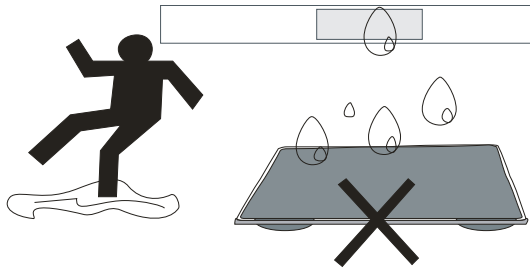Beware of slippery platter/not with wet feet.

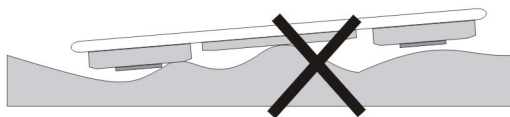The scale should be placed on a firm and level surface (avoid carpet).

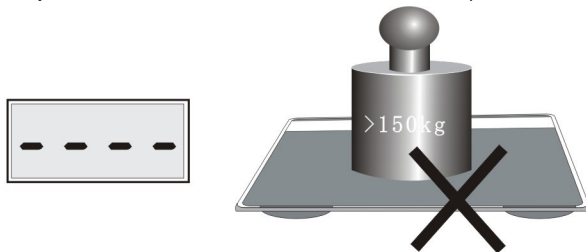Maximum weighting capacity is 150 kg



Clean the scale with soft cloth


**Technical Data:**

Max./Min. capacity:	150kg(330 lb)/3kg(6.6 lb)
Graduation:	100g(0.2 lb)
LCD:	75x35mm(2.9"x1.4")
Tempered glass platform:	285x285x6mm
Batteries:	2xCR2430


**Use of Bathroom Scale**

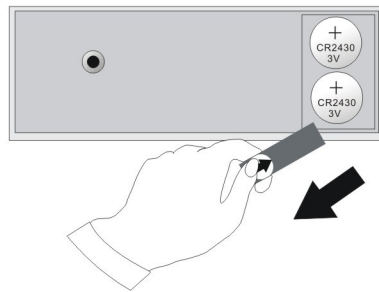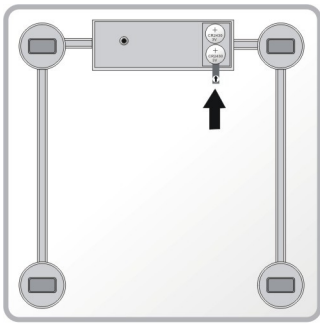
The battery compartments of the scale are situated on the back. Before the first use, open the battery compartments of the scale and pull out the protective slip

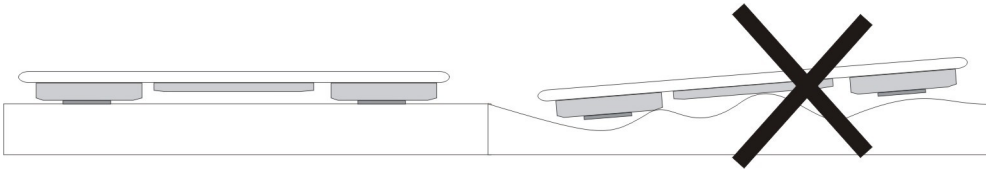Put the scale on a firm and level surface

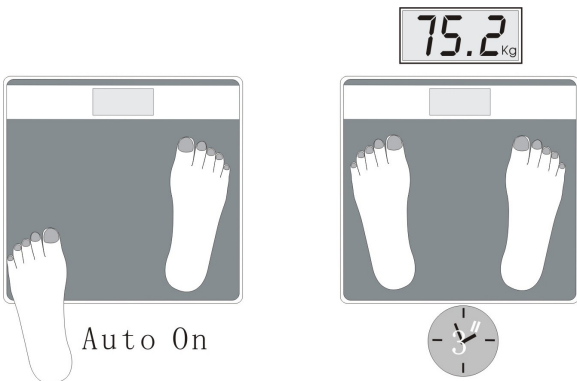
1. As soon as you stand on the scale or press the platform, the scale will switch on automatically.

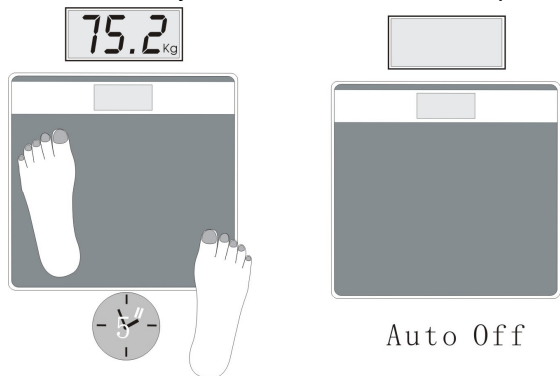
2. Stand steadily on the scale, it will switch on automatically and show the weight. The weight value will go to steady and be locked, when the unit appears following the weight value, and this is the result. For about 5 seconds after the unit appears, it will shut off automatically.

### Changing units

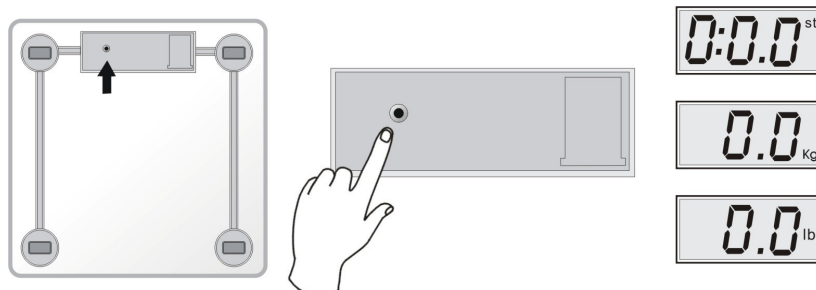You can change the units among kg, lb, st by pressing the button on the back of the scale when it is power on.

### Inserting/ Changing Batteries:

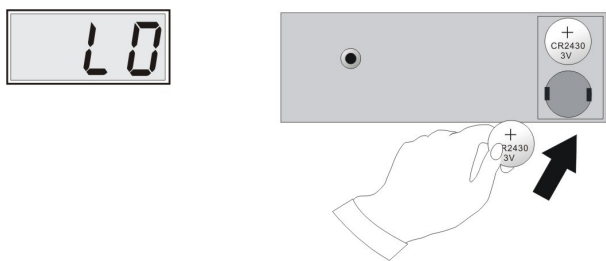The battery compartments of the scale are situated on the back. Before the first use, open the battery compartments of the scale and insert two new CR2430 Lithium batteries respectively in the right polarity. When - LO - displays on the screen, it means batteries are almost worn out, please replace new batteries. Important: When you change your batteries, be careful!! You must eliminate the batteries according to the rules in force. Dump them in an appropriate place in order to assure their recycling and the protection of environment.

### Tips about Body composition

Body composition reflects the internal structure of human body development. Different internal structure, so the mechanism movement differs, furthermore, certain proportion between each composition is required, so the normal physiology mechanism can last. Once some proportion discords, the normal physiology mechanism is destroyed; it will impede the human body's normal growth and health.

Body composition also shows the body status, physique characteristics and stature, in which fat composition reveals the fat or thin level. Therefore, it has great significance on the formulating fatness standard and shape appraisal.

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#### Guidelines for protection of the environment

This appliance should not be put into the domestic garbage at the end of its useful life, but must be disposed of at a central point for recycling of electric and electronic domestic appliances. This symbol on appliance, instruction manual and packaging puts your attention to this important issue. The materials used in this appliance can be recycled. By recycling used domestic appliances you contribute an important push to the protection of our environment. Ask your local authorities for information regarding the point of recollection..

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