
waterpik[®]
treat yourself better™

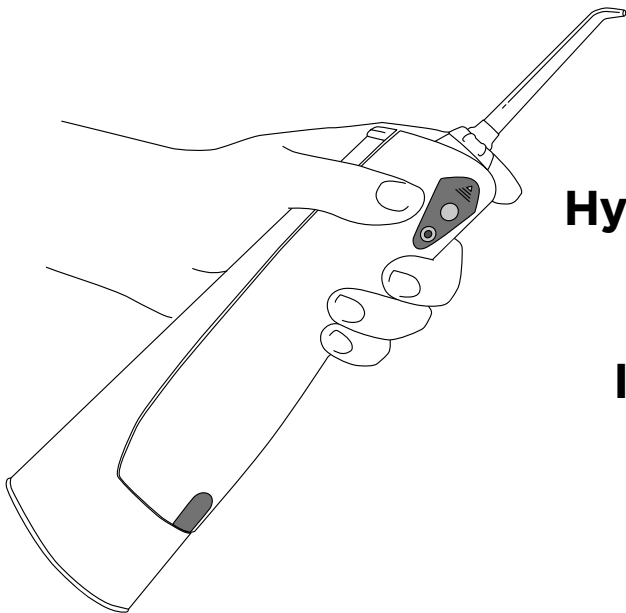**Waterpik[®] Water Flosser
Model WP-360**

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.....
www.waterpik.com

WELCOME

Congratulations! By selecting the Waterpik® Water Flosser you are now on your way to better gum health, and a healthier smile!

With a clinically proven combination of water pressure and pulsations, the Waterpik® Water Flosser removes harmful bacteria deep between teeth and below the gumline – where brushing alone won't reach.


We are so convinced that you will be happy with your new Waterpik® Water Flosser, that we back it up with a 14-day guarantee. If you aren't satisfied that your gums are healthier after using the Waterpik® Water Flosser for 14 consecutive days, feel free to return it to us for a full refund.

Thank you for choosing the Waterpik® Water Flosser.

Traditional Oral Care


Traditional brushing can't easily reach hidden food debris and harmful bacteria.

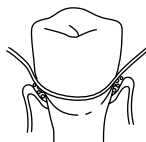String flossing has limited access below the gumline.

Oral Care With The Waterpik® Water Flosser

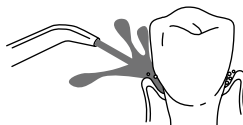The Waterpik® Water Flosser removes trapped food debris and harmful bacteria.

TABLE OF CONTENTS

Important Safeguards	4
Product Description	6
Getting Started	8
Care of Your Waterpik® Water Flosser	13
Troubleshooting Guide	14
Limited One-Year Warranty	15

IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not immerse in water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the charger cord for damage before the first use and during the life of the product.

WARNING:


To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:


- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. or your local distributor, which can be found at www.waterpik.com.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. Without appropriate attachments recommended by Water Pik, Inc., this product is capable of producing pressures that may cause serious damage in these areas. See operating instructions for correct usage.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories that are recommended by Water Pik, Inc.

- Fill reservoir with water or other dental professional recommended solutions only.
- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use mouthwashes/rinses that contain iodine, bleach, or tea tree oil.
- Do not use while wearing any oral jewelry. Remove prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- For household use only. Do not use outdoors.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Children should be supervised to ensure that they do not play with the appliance.


SAVE THESE INSTRUCTIONS.

PRODUCT DESCRIPTION

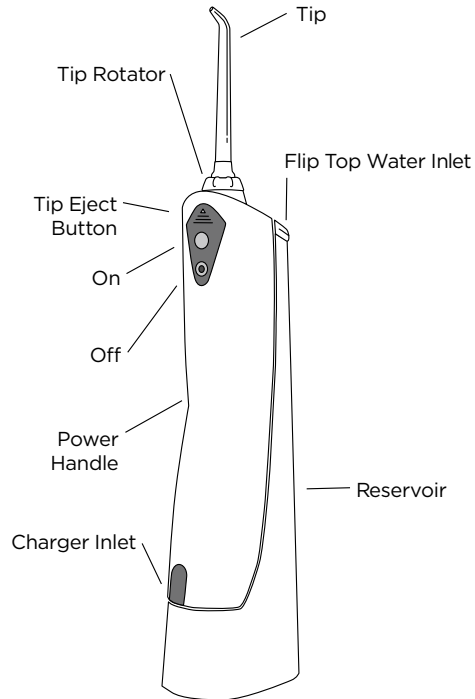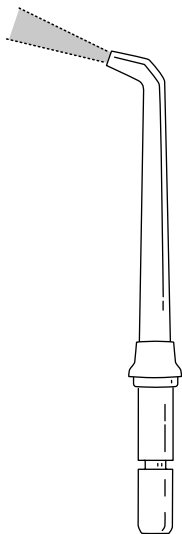
Deep Cleaning

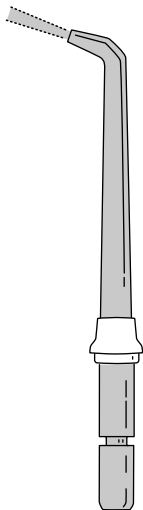Gentle Cleaning

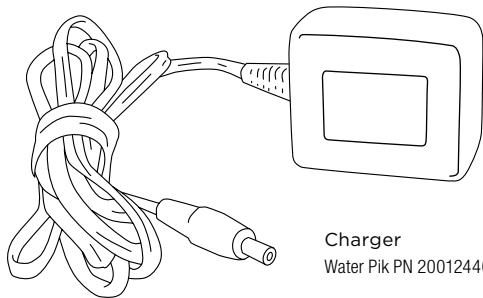Charger
Water Pik PN 20012440


Replacement tips/attachments may be purchased directly from Water Pik, Inc.

Call 1-800-525-2774 (USA) or 1-888-226-3042 (Canada) or visit our web site at www.waterpik.com.

GETTING STARTED


Important! Turn Water Flosser off and charge for 24 hours before first use.

Plug the cord into any standard electrical outlet. Make sure you push the switch into the "OFF" position before charging. The unit will not charge if the switch is in the "ON" position. See "Care of your Waterpik® Water Flosser" for guidance on how to look after your battery. The Waterpik® Water Flosser will stay charged for up to 2 weeks depending on use.

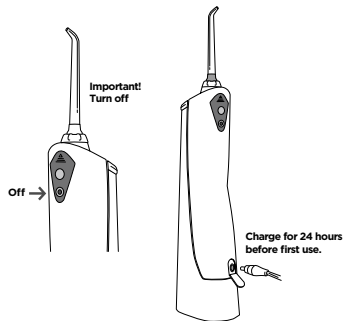
Preparing the Reservoir

Slide Water Flosser power handle into back of reservoir and press firmly to snap into place.


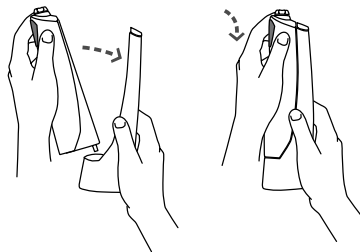
Lift flip top on reservoir and fill the reservoir with lukewarm water.

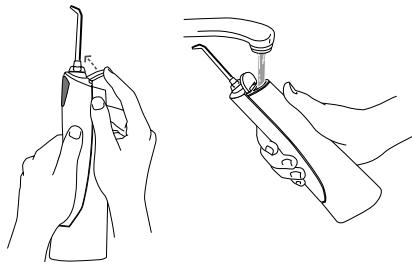Inserting and Removing Tips


Insert tip into the center of the knob at the top of the Water Flosser handle. The colored ring will be flush with the end of the knob if the tip is correctly locked into place.

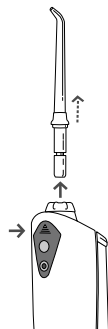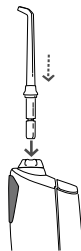To remove the tip from the handle, press the tip eject button and pull tip from the handle.

Adjusting the Pressure Setting


You can adjust the pressure setting on your Water Flosser by simply changing tips. There are two tips available—a high pressure tip for deeper cleaning and a low pressure tip for gentle cleaning.

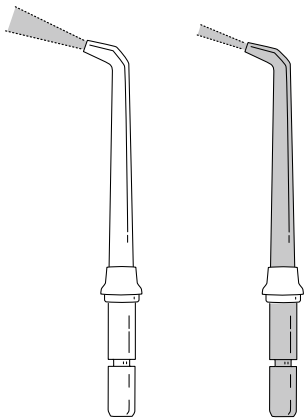
High pressure for deep cleaning

Low pressure for gentle cleaning

Keep unit upright during use for best results.



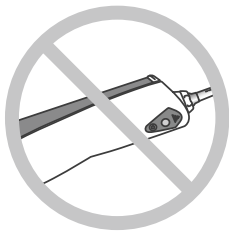
Direct the stream by turning the knob at the top of the handle.

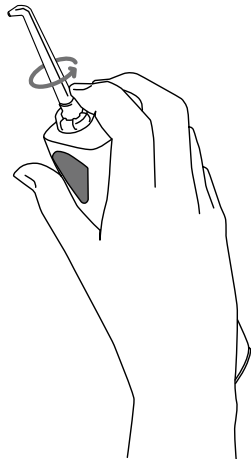
Recommended Technique

Lean low over sink and place the jet tip in mouth. Aim the jet tip toward teeth, and turn the unit on.

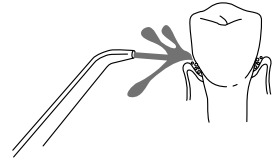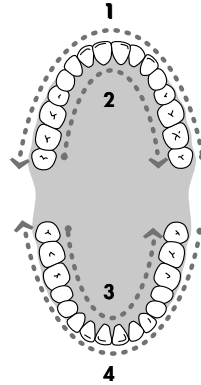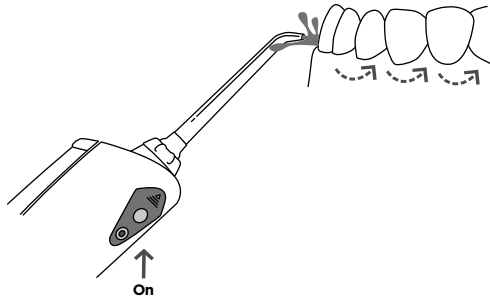
Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working toward the front teeth. Glide jet tip along gumline and pause briefly between teeth. Continue until all areas around and between the teeth have been cleaned.

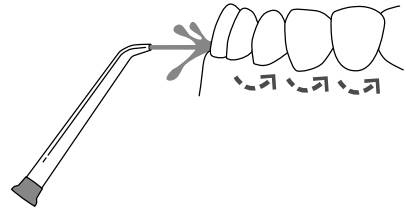
When Finished

Turn the unit off. Empty any liquid left in the reservoir.

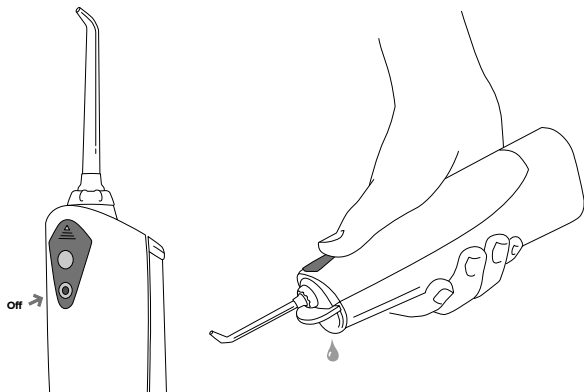If desired, the reservoir may be removed from the power handle by pushing down on the top of power handle and pulling it away from the reservoir.

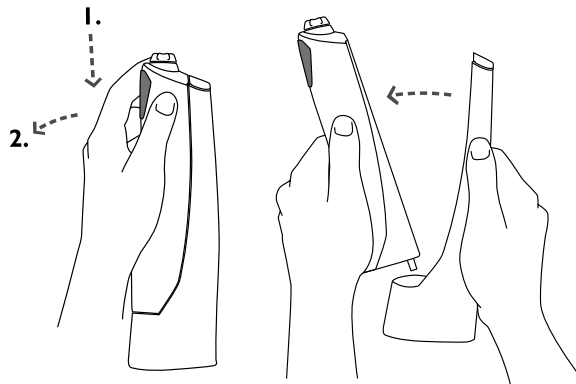Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. Ask your dental professional to recommend the appropriate solution for you. Note: use of some antibacterial solutions could shorten the life of your Waterpik® Water Flosser.

After using any special solution in the reservoir, take a moment after every use to rinse the unit. Partially fill the reservoir with warm tap water. Aim the jet tip into your sink. Turn on the unit and let it run until empty. This will help keep the tip from clogging and help your unit to last longer.

CARE OF YOUR WATERPIK® WATER FLOSSER

Never immerse the unit in water. Wipe the exterior when necessary with a soft cloth and mild non-abrasive cleanser. Reservoir may be washed with soapy water or in top rack of dishwasher.

Before cleaning, unplug from the electrical outlet.

Before exposing the Water Flosser to freezing temperatures, remove the reservoir and run the unit until completely empty.

Battery Care

Before first use charge for 24 hours to obtain full charge.

To charge your Water Flosser, plug the cord into any standard electrical outlet. Make sure you push the switch into the 'OFF' position before charging. The unit will not charge if the switch is in the 'ON' position.

If you use your Water Flosser once a day or less, it is recommended you charge it overnight, once a week. If you use it twice a day or more, charge the unit more frequently. Do not run the batteries down completely as this may shorten the lifespan of the batteries and will require charging the unit for 24 hours to obtain a full charge.

If the Water Flosser is to be stored for an extended period of time, charge it for 24 hours before storing. Never store the Water Flosser for extended periods unless it is fully charged. To maximize the lifespan of the batteries the battery charger may be kept connected to the unit continuously, except for when the product is in use.

Battery Disposal

The unit contains a Nickel Metal Hydride battery. At the end of the battery life, please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable.

Service Maintenance

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all your service needs go to www.waterpik.com.

Some replacement parts and accessories for the Water Flosser are also available from Water Pik, Inc. or your local distributor which can be found at www.waterpik.com. You may also call our toll-free customer service number at 1-800-525-2774 (USA) or 1-888-226-3042 (Canada).

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

TROUBLESHOOTING GUIDE

PROBLEM	CAUSE	CORRECTIVE ACTION
Tips Leakage between the tip and the tip rotator.	1) Tip not fully engaged.	1) Remove tip and reinstall.
Water Flosser Reservoir leaks.	1) Liquid was left in reservoir. 2) Power handle not correctly engaged with reservoir. 3) Water inlet is open.	1) Always empty the reservoir after use. 2) Remove reservoir then reinstall as per instructions. 3) Close water inlet.
Inadequate pressure/No fluid flow.	1) Unit is upside down. 2) Insufficient liquid in reservoir. 3) Reservoir connection clogged. 4) Tip is clogged with debris	1) Tip must be higher than reservoir for correct water flow. 2) Fill with liquid. 3) Fill reservoir with warm tap water and one cap of vinegar. Run until empty. Rinse with full reservoir of tap water. 4) Flush tip under faucet or replace tip.
Unit does not charge.	1) Mains power socket is not live. 2) Switch is in ON position during charging.	1) Turn ON the mains power socket; try an alternative mains power socket. 2) Turn unit switch to OFF position for charging.
Unit doesn't start.	1) Battery not charged.	1) Be sure outlet is functioning; test with another appliance. See "Care of your Waterpik® Water Flosser" for details on how to charge the battery.

Still have questions? Go to www.waterpik.com.

If none of these corrective actions are successful, call the Water Pik, Inc. Customer Service. Toll-Free Number: 1-800-525-2774 (U.S.A.) 1-888-226-3042 (Canada).

LIMITED ONE-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 1 year from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

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Fort Collins, CO 80553-0001 USA

www.waterpik.com

Waterpik® is a trademark of Water, Pik, Inc. registered in Argentina, Australia, Canada, Chile, China, Czech Republic, Hungary, India, Mexico, New Zealand, Pakistan, Poland, Russian Federation, Singapore, South Africa, Switzerland, Taiwan, Ukraine, United Kingdom, and the United States.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Canada, EU, Japan, Mexico, and the United States.

Treat Yourself Better™ is a trademark of Water Pik, Inc.

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644 Eden Park Drive
Cincinnati, Ohio 45202
t 513.961.6225 | f 513.281.3611

JOB INFORMATION

**Product Description: WaterPik
Instruction Manual Model WP-360**

Client: **WaterPik**

UPC No:

Manufacturer:

Last Revision Date: **10/26/09**

Blueprint No: 14x4.5 flat, 7x4.5 folded

STUDIO INFORMATION

Job No: **2272_0003_003**

File No: **2272_0003_003 WP 360 Contless IG.indd**

Studio Contact: **Client Service: Mark Morse**

Creative Lead: Scott Smith

Production Mgr: Doug Bennett

CLIENTS PACKAGING APPROVALS

	INITIALS	DATE
BRAND		
LEGAL		

IMPORTANT NOTICE

laga has checked this artwork for accuracy. Final approval is the client's responsibility. Please double check for errors before reproduction. We recommend that color proofs be submitted to us for approval before production.

TECHNICAL INFORMATION

Application/Version: **Adobe Indesign CS3**

Layers/Pages: **5/1**

Placed Images:

Fonts Used: **Refer to file**

Separator: **Name**

Release Date:

Approved For Release _____ / /

COMPUTER EQUIVALENTS

PANTONE
PROCESS
BLACK

TINTS / BUILT COLORS**brandimage APPROVALS**

gmb	/ /	/ /
Artist/Operator	Date	Production Director Date
Account Executive	/ /	/ /
	Date	Proofreader Date
Design Director	/ /	/ /
	Date	Proofreader Date

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©2009 Water Pik, Inc.