SANITAS

SBF 08

- D** **Glas-Diagnosewaage**
Gebrauchsanweisung 2–7
- GB** **Glass diagnostic scale**
Instruction for Use 8–13
- F** **Pèse-personne impédancemètre
en verre**
Mode d'emploi 14–19
- I** **Bilancia diagnostica in vetro**
Istruzioni per l'uso 20–25
- RUS** **Стеклянные диагностические весы**
Инструкция по применению 26–32

CE

Service-Adresse:
Hans Dinslage GmbH, Riedlinger Str. 28
88524 Uttenweiler, GERMANY
Tel.: +49 (0) 73 74 - 91 57 66
E-mail: service@sanitas-online.de

DEUTSCH

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
ENGLISH

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

1. Device description

1. Electrodes
2. Display
3. „Down“ key
4. „SET“ key
5. „Up“ key

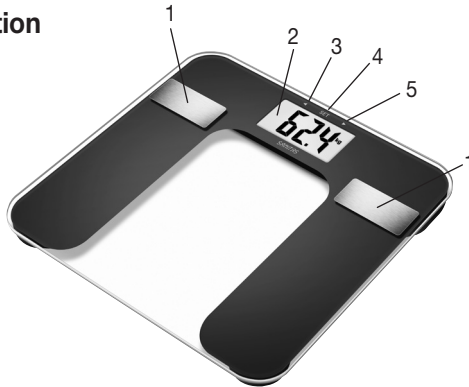2. Notes

Safety notes

- **The scales must not be used by persons with medical implants (e.g. heart pacemakers). Otherwise their function could be impaired.**
- Do not use during pregnancy. Inaccurate measurements may be caused by the amniotic fluid.
- Caution! Do not stand on the scale with wet feet and do not step on the scale when its surface is wet – danger of slipping!
- Batteries are highly dangerous if swallowed. Keep batteries and scale out of reach of small children. If batteries are swallowed, get medical help immediately.
- Keep children away from packing materials (danger of suffocation).
- Batteries must not be charged or reactivated by other means, taken apart, thrown into a fire or short-circuited.

Notes on handling batteries

- Swallowing batteries can be extremely dangerous. Keep the batteries and scale out of the reach of small children. Should a child swallow a battery, seek medical assistance immediately.
- Replace weak batteries before they discharge completely.
- Always replace all batteries at the same time and use batteries of the same type.
- Batteries must not be recharged, taken apart, thrown into an open fire or short circuited.
- Leaking batteries may damage the device. If you do not intend to use the device for longer periods, remove the batteries from the battery compartment.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Batteries can contain toxins that are harmful to health and the environment. Always dispose of batteries in accordance with applicable legal regulations. Do not dispose of batteries with the normal household waste.
- Do not throw batteries into fire. Explosion hazard!

General notes

- The unit is for personal use only and is not intended for medical or commercial applications.
- Please note that technically related measuring tolerances are possible, as the scale is not calibrated for professional medical applications.
- Ages 10 ... 100 years and height settings from 100 ... 220 cm (3-03"–7-03") can be preset. Load capacity: max. 150 kg (330 lb, 24 st). Result in steps of 100g (0.2 lb, 1/4 st). Measurements of body fat, body water and muscle density in steps of 0.1%.

- On delivery the scale is set to the units “cm” and “kg”. A switch located on the back of the scale enables you to change to “inch”, “pound” and “stones” (lb, St).
- Place the scale on a firm level floor; a firm floor covering is required for correct measurement.
- The scale should be cleaned occasionally with a damp cloth. Do not use abrasive detergents and never immerse the scale in water.
- Protect the scale against hard knocks, moisture, dust, chemicals, major temperature fluctuations and heat sources which are too close (stove, heating radiators).
- Repairs may only be carried out by Beurer customer service or authorized dealers. Before submitting any complaint, first check the batteries and replace them if necessary.

Storage and maintenance

The accuracy of the measurements and service life of the device depend on its careful handling:

IMPORTANT

- The unit should be cleaned occasionally. Do not use abrasive detergents and never immerse the unit in water.
- Ensure that no liquid gets on the scale. Never dip the scale into water. Never wash it under running water.
- Do not place any objects on the scale when it is not being used.
- Protect the unit against hard knocks, moisture, dust, chemicals, major temperature fluctuations and heat sources which are too close (stoves, heating radiators).
- Do not press the buttons using excessive force or with pointed objects.
- Do not expose the scale to high temperatures or strong electromagnetic fields (e.g. mobile phones).

3. Information on the diagnostic scale

The measuring principle of the diagnostic scale

This scale operates according to the B.I.A. principle (bioelectric impedance analysis). This enables physical relationships to be measured within seconds by means of an undetectable, completely harmless electric current. The body fat percentage and other physical relationships in the body can be determined by measuring the electrical resistance (impedance), and the inclusion in a calculation of constants and individual parameters (age, height, gender, degree of activity).

Muscle tissue and water have good electrical conductivity, and therefore lower resistance.

Bones and fatty tissue, on the other hand, have low conductivity, as fat cells and bones conduct hardly any current as a result of their very high resistance.

Please be aware that values obtained from the diagnostic scale represent only an approximation of actual analytical medical data. Only a specialist physician can accurately determine body fat, body water, muscle percentage and bone structure using medical procedures (e.g. computed tomography).

General tips

- If possible, always weigh yourself at the same time of day (preferably in the morning), after using the toilet, on an empty stomach and unclothed, in order to achieve comparable results.
- Important for the measurement: body fat can only be measured when barefoot and with the soles of the feet slightly moist. Completely dry soles can result in unsatisfactory measurements due to inadequate conductivity.
- Stand still during the measurement.
- Wait a few hours after unaccustomed strenuous activity.
- Wait approx. 15 minutes after getting up in the morning to allow the water in the body to distribute.

Limitations

When measuring body fat and other values, deviating and implausible results may occur in:

- children under approx. 10 years,
- competitive athletes and body builders,
- pregnant women,
- persons with fever, undergoing dialysis treatment or with symptoms of edema or osteoporosis,
- persons taking cardiovascular medication,
- persons taking vascular dilating or vascular constricting medication,

– persons with considerable anatomic deviations of the legs in relation to total body size (length of the legs considerably shortened or lengthened).

4. Batteries

If provided, pull out the battery insulating strip at the battery compartment or remove the protective film from the battery and insert the battery with the correct polarity. If the scale fails to operate, remove the battery completely and reinsert it.

Your scale is provided with a “Change battery” indicator. If the scale is operating with a battery which is too weak, “Lo” appears on the display and the scale will switch off automatically. In this case, the battery must be exchanged (1 x 3 V lithium battery CR2032).

5. Operation

5.1 Weight measurement

Place the scale on a firm level floor (no carpet); a firm floor covering is required for correct measurement.

Step on the scale. It immediately begins to measure. Stand still with your weight distributed equally on both legs. If the display appears continuously, the measurement has finished.

If you step off the weighing area, the scale switches off after a few seconds.

5.2 Setting user data

To measure your body fat percentage and other physical data, you must enter your personal user parameters.

The scale has 10 memory positions in which you and other members of your family can save and recall personal settings.

Switch on the scale (tap the weighing area). Wait until the display shows “0.0”.

Then press “SET”. The first memory position will appear flashing on the display. Now you can enter the following settings:

| | |
|--------------------|----------------------------------|
| Memory position | 1 to 10 |
| Body size | 100 to 220 cm (3'-03" to 7'-03") |
| Age | 10 to 100 years |
| Sex | male (♂), female (♀) |
| Degree of activity | 1 to 5 |

You can enter the relevant settings with short or long pressure on the up button ▲ or down button ▼. Confirm the settings in each case by pressing “SET”.

The values you've just set are displayed again consecutively. The scale then automatically switches off.

Degree of activity

Selection of the degree of activity must refer to the medium and long term.

– **Degree of activity 1:** No physical activity.

– **Degree of activity 2:** Low physical activity.

A small amount of light physical effort (e.g. short walks, light garden work, gymnastic exercises).

– **Degree of activity 3:** Medium physical activity.

Physical effort for 30 minutes at least 2 to 4 times a week.

– **Degree of activity 4:** High physical activity.

Physical effort for 30 minutes at least 4 to 6 times a week.

– **Degree of activity 5:** Very high physical activity.

Intensive physical effort, intensive training or hard physical work for at least one hour daily.

5.3 Taking measurements

After all parameters have been entered, you can now measure your weight, body fat and the other data.

– Press the “SET” button and select the memory space where you would like to save your basic personal data by pressing the “up” ▲ or “down” ▼ button the appropriate number of times. These are displayed in sequence until “0.0” and the selected gender appear.

- Step onto the scale barefoot and make sure you're standing on both electrodes. First your body weight is determined and displayed.
- Remain standing on the scale; now the body fat and body water analysis is performed. This can take a few seconds.

Important: There may be no contact between feet, legs, calves and thighs. Otherwise the measurement will not be performed correctly.

The following data are automatically displayed consecutively:

1. Weight
2. Body fat in % (FAT)
3. Body water in % (BW)
4. Muscle content (MUS)
5. Body mass index (BMI)

Now all measured values are displayed consecutively and the scale switches off.

6. Evaluation of results

Body fat percentage

The following body fat percentages are for guidance (contact your physician for further information).

Man

| Age | very good | good | average | poor |
|--------|-----------|--------|----------|--------|
| 10-14 | <11% | 11-16% | 16,1-21% | >21,1% |
| 15-19 | <12% | 12-17% | 17,1-22% | >22,1% |
| 20-29 | <13% | 13-18% | 18,1-23% | >23,1% |
| 30-39 | <14% | 14-19% | 19,1-24% | >24,1% |
| 40-49 | <15% | 15-20% | 20,1-25% | >25,1% |
| 50-59 | <16% | 16-21% | 21,1-26% | >26,1% |
| 60-69 | <17% | 17-22% | 22,1-27% | >27,1% |
| 70-100 | <18% | 18-23% | 23,1-28% | >28,1% |

Woman

| Age | very good | good | average | poor |
|--------|-----------|--------|----------|--------|
| 10-14 | <16% | 16-21% | 21,1-26% | >26,1% |
| 15-19 | <17% | 17-22% | 22,1-27% | >27,1% |
| 20-29 | <18% | 18-23% | 23,1-28% | >28,1% |
| 30-39 | <19% | 19-24% | 24,1-29% | >29,1% |
| 40-49 | <20% | 20-25% | 25,1-30% | >30,1% |
| 50-59 | <21% | 21-26% | 26,1-31% | >31,1% |
| 60-69 | <22% | 22-27% | 27,1-32% | >32,1% |
| 70-100 | <23% | 23-28% | 28,1-33% | >33,1% |

A lower value is often found in athletes. Depending on the type of sports, training intensity and physical constitution, values may result which are below the recommended values stated. It should, however, be noted that there could be a danger to health in the case of extremely low values.

Body water

The body water percentage is normally within the following ranges:

Man

| Age | poor | good | very good |
|--------|------|--------|-----------|
| 10-100 | <50% | 50-65% | >65% |

Woman

| Age | poor | good | very good |
|--------|------|--------|-----------|
| 10-100 | <45% | 45-60% | >60% |

Body fat contains relatively little water. Therefore persons with a high body fat percentage have body water percentages below the recommended values. With endurance athletes, however, the recommended values could be exceeded due to low fat percentages and high muscle percentage.

Body water measurement with this scale is not suitable for drawing medical conclusions, for example concerning age-related water retention. If necessary ask your physician. Basically, a high body water percentage should be the aim.

Muscle percentage

The muscle percentage is normally within the following ranges:

Man

| Age | low | normal | high |
|--------|------|---------|-------|
| 10-14 | <44% | 44-57 % | >57 % |
| 15-19 | <43% | 43-56 % | >56 % |
| 20-29 | <42% | 42-54 % | >54 % |
| 30-39 | <41% | 41-52 % | >52 % |
| 40-49 | <40% | 40-50 % | >50 % |
| 50-59 | <39% | 39-48 % | >48 % |
| 60-69 | <38% | 38-47 % | >47 % |
| 70-100 | <37% | 37-46 % | >46 % |

Woman

| Age | low | normal | high |
|--------|------|---------|-------|
| 10-14 | <36% | 36-43 % | >43 % |
| 15-19 | <35% | 35-41 % | >41 % |
| 20-29 | <34% | 34-39 % | >39 % |
| 30-39 | <33% | 33-38 % | >38 % |
| 40-49 | <31% | 31-36 % | >36 % |
| 50-59 | <29% | 29-34 % | >34 % |
| 60-69 | <28% | 28-33 % | >33 % |
| 70-100 | <27% | 27-32 % | >32 % |

Body weight/Body mass index

The body mass index (BMI) is a number that is often called upon to evaluate body weight. The number is calculated from body weight and height, the formula is – body mass index = body weight : height². The measurement unit for BMI is [kg/m²].

Man <20 years

| Age | Under-weight | Normal weight | Over-weight |
|-----|--------------|---------------|-------------|
| 10 | 12,7–14,8 | 14,9–18,4 | 18,5–26,1 |
| 11 | 13,0–15,2 | 15,3–19,1 | 19,2–28,0 |
| 12 | 13,3–15,7 | 15,8–19,8 | 19,9–30,0 |
| 13 | 13,7–16,3 | 16,4–20,7 | 20,8–31,7 |
| 14 | 14,2–16,9 | 17,0–21,7 | 21,8–33,1 |
| 15 | 14,6–17,5 | 17,6–22,6 | 22,7–34,1 |
| 16 | 15,0–18,1 | 18,2–23,4 | 23,5–34,8 |
| 17 | 15,3–18,7 | 18,8–24,2 | 24,3–35,2 |
| 18 | 15,6–19,1 | 19,2–24,8 | 24,9–35,4 |
| 19 | 15,8–19,5 | 19,6–25,3 | 25,4–35,5 |

Woman <20 years

| Age | Under-weight | Normal weight | Over-weight |
|-----|--------------|---------------|-------------|
| 10 | 12,7–14,7 | 14,8–18,9 | 19,0–28,4 |
| 11 | 13,0–15,2 | 15,3–19,8 | 19,9–30,2 |
| 12 | 13,3–15,9 | 16,0–20,7 | 20,8–31,9 |
| 13 | 13,7–16,5 | 16,6–21,7 | 21,8–33,4 |
| 14 | 14,2–17,1 | 17,2–22,6 | 22,7–34,7 |
| 15 | 14,6–17,7 | 17,8–23,4 | 23,5–35,5 |
| 16 | 15,0–18,1 | 18,2–24,0 | 24,1–36,1 |
| 17 | 15,3–18,3 | 18,4–24,4 | 24,5–36,3 |
| 18 | 15,6–18,5 | 18,6–24,7 | 24,8–36,3 |
| 19 | 15,8–18,6 | 18,7–24,9 | 25,0–36,2 |

Man/Woman ≥20 years

| Age | Under-weight | Normal weight | Over-weight |
|-----|--------------|---------------|-------------|
| ≥20 | 15,9–18,4 | 18,5–24,9 | 25,0–40,0 |

Source: Bulletin of the World Health Organization 2007;85:660-7; Adapted from WHO 1995, WHO 2000 and WHO 2004.

Note that with very muscular physiques (body builders), the BMI interpretation will show overweight. The reason for this is that no account is taken of well above average muscle mass in the BMI formula.

Results in relation to time

Note that it is only the long-term trend which counts. Short period fluctuations in weight over a few days are mostly the result of a loss of fluid.

The interpretation of the results will depend on changes in your: overall weight and body fat, body water and muscle percentages, as well as on the period during which these changes take place. Rapid changes within days must be distinguished from medium term changes (over weeks) and long term changes (months).

A basic rule is that short term changes in weight almost exclusively represent changes in water content, whereas medium and long term changes may also involve the fat and muscle percentages.

- If your weight reduces over the short term, but your body fat percentage increases or remains the same, you have merely lost water – e.g. after a training session, sauna session or a diet restricted only to rapid weight loss.
- If your weight increases over the medium term and the body fat percentage falls or stays the same, then you could have built up valuable muscle mass.

If your weight and body fat percentage fall simultaneously then your diet is working – you are losing fat mass. Ideally you should support your diet with physical activity, fitness or power training. By this means, you can increase your muscle percentage over the medium term.

7. Incorrect measurement

If the scale detects an error during weighing, „FAT Err“ or “Err” appears in the display.


If you step onto the scale before “0.0” appears in the display, the scale will not operate properly.

| Possible causes of error: | Remedy: |
|--|---|
| – The maximum load-bearing capacity of 150 kg (330 lbs) was exceeded. | – Only weigh the maximum permissible weight. |
| – The electrical resistance between the electrodes and the soles of your feet is too high (e.g. with heavily callused skin). | – Repeat weighing barefoot. Slightly moisten the soles of your feet if necessary. Remove the calluses on the soles of your feet if necessary. |
| – Your body fat lies outside the measurable range (less than 5 % or greater than 80 %). | – Repeat weighing barefoot. – Slightly moisten the soles of your feet if necessary. |
| – Your percentage of water lies outside the measurable range (less than 43 % or greater than 75 %) | – Repeat weighing barefoot. – Slightly moisten the soles of your feet if necessary. |

8. Disposal

Used, fully discharged batteries must be disposed of in a specially labeled collection container, at toxic waste collection points or through an electrical retailer. You are under legal obligation to dispose of batteries correctly.

Note: You will find these markings on batteries containing harmful substances: Pb = battery containing lead, Cd = battery containing cadmium, Hg = battery containing mercury.

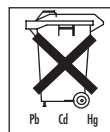Please dispose of the device in accordance with the directive – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please refer to the local authorities responsible for waste disposal.

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